



BurDental 博登科技



博登科技有限公司
B&D Technologies Co.,Ltd

Dental Burs 牙科车针、磨头

Vol. 8 : 2024

About Us

B&D Technologies Company Limited are specialized in Rotary things of dental, Jewelry, Beauty or Hobbies & DIY (Do it Yourself) .

Our products include Diamond Burs and discs, Tungsten Carbide Burs or cutters, Flexible Rubber Burs, Wheels or Points, Polishing Brushes, Mounted Stones, CAD/CAM Milling Burs or any Diamond Stones or Polishers.

Our Rotatory Instruments are widely applied to Cut, Grind, Finish, Polish any materials of Zirconia, Lithium Disilicates, Porcelain (Ceramics) , Alloy or Metal, Acrylic (Resin) , PEEK, Composites etc . They can match any Handpiece of FG (Friction Grip) , RA/CA(Right Angle), HP conveniently and precisely .

We have obtained the certificates of ISO 13485 and CE since the year of 2015 . We provide innovative and professional products and service.

Whoever you are dentist or technician, we always make you relaxing, safe and efficient. We are your best friend in dentistry, Jewellery or Hobbies.



ISO 13485



CE 0197



Established in 2007



Laboratory Part 技工类

Page 33-66 第33-66页



Diamond bur 金刚石磨头

Page 34-35 第34-35页

Fully-sintered Diamond bur 全砂烧结磨头

Page 36-37 第36-37页

Sintered Diamond Disc Diamond Disc

烧结切片、金刚石切片

Page 38-40 第38-40页

CAD / CAM Milling Bur 铣刀

Page 41-43 第41-43页

Tungsten Carbide bur 钨钢车针

Page 44-47 第44-47页

Diamond Stones & Polisher 陶瓷(橡胶) 结合剂金刚石磨头 (全瓷、全锆磨头)

Page 48-52 第48-52页

Silicon Polisher 硅胶抛光磨头

Page 53-54 第53-54页

Mounted stone, Sandpaper 车石、砂纸棒

Page 55 第55页

Polishing Brushes, Felt burs Mandrels

抛光刷、毛毡磨头、夹针

Page 56-58 第56-58页

Hard Z Cutter etc . 硬锆磨头等

Page 59 第59页

Kit 套装

Page 60-62 第60-62页

Samples Collection 样品册

Page 63 第63页

FlexSoft Partial Denture 隐形义齿胶

Page 64 第64页

Diamond Bur HP (Electroplated) 金刚石磨头 (电镀)



Sharp and Durable,
Grind or polish Ceramics
锋利、耐用
磨削平稳、不偏摆
不留黑印、适合车瓷、



Grit


- Black / Super Coarse
- Green / Coarse
- Blue / Medium
- Red / Fine

Fine
Red
细砂
红圈



编号 No. This grit can grind Zirconia 可用于车氧化锆


Medium
Blue
标准砂
蓝圈

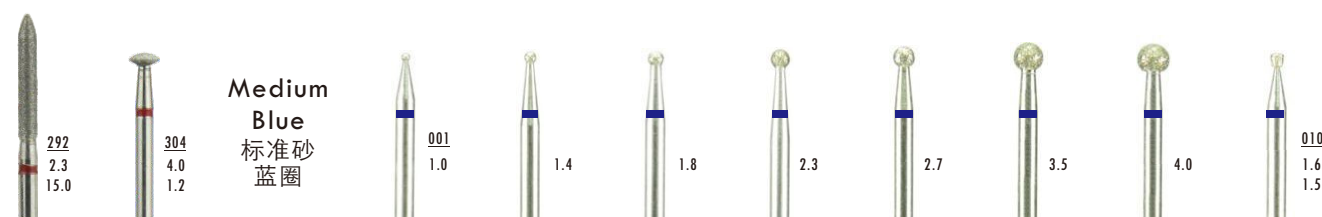

编号 No.


Round / 圆球

Inverted Cone

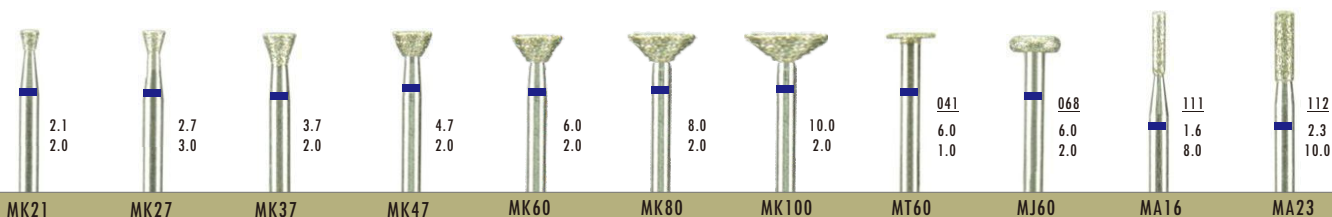


Wheel / 平轮

Cylinder / 平柱

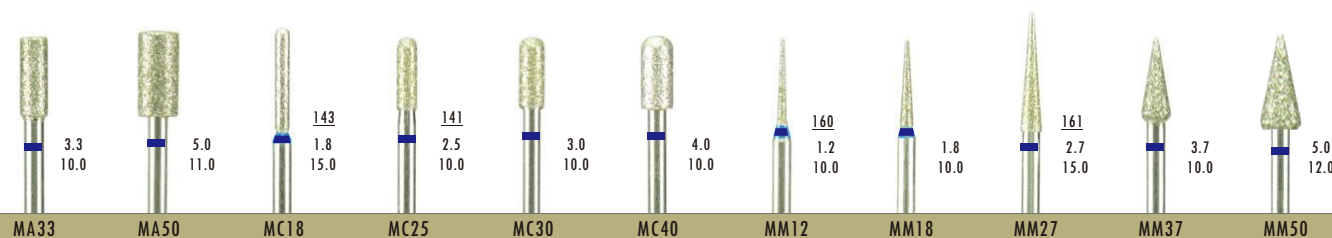

Round Cylinder / 圆头柱


Needle / 尖针

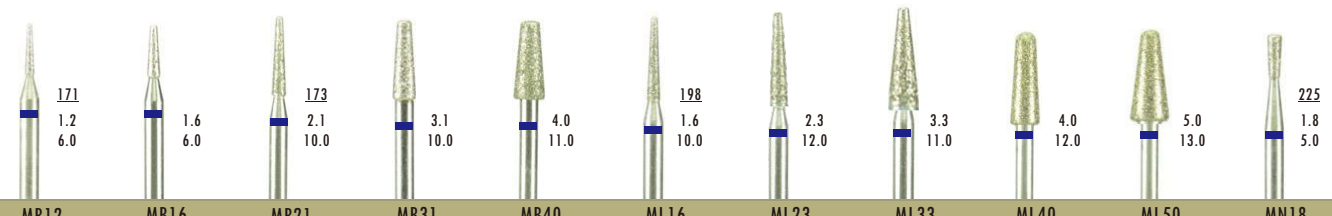

Taper / 平锥

Round Taper / 圆头锥

Inverted Cone

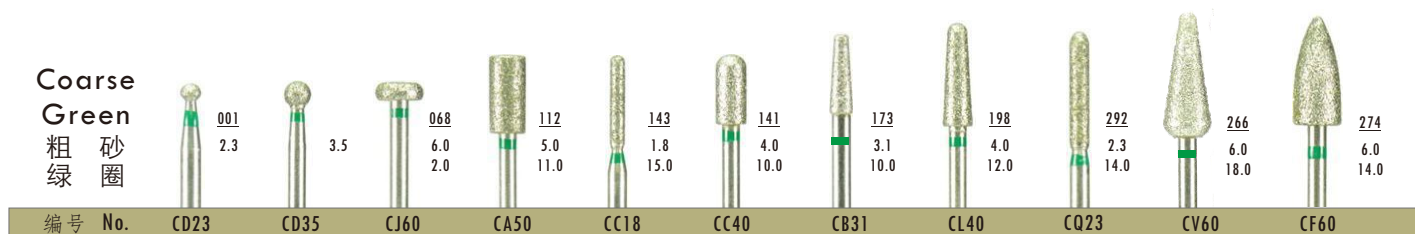


Transparent Showing plates 展示台 120x120x8mm



No. 1216

7 x 7 holes
1.60mm, for FG burs

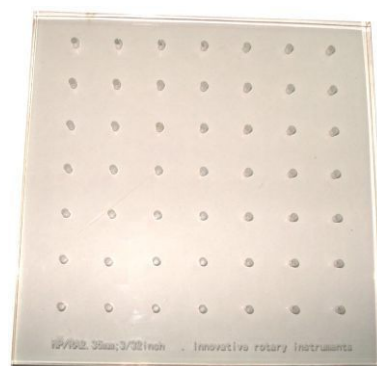

No. 1223

7 x 7 holes
2.35mm, for RA or HP burs.


Fully-Sintered Diamond Bur . HP

全砂烧结磨头



Extremely Durable,
Grind or polish smoothly
Especially for Metal

从里到外都是砂
平稳、耐用，适合车金

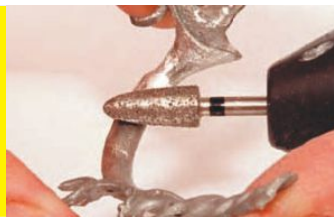


Grit

- Black/Super Coarse
- Green/Coarse
- Blue/Medium
- Red/Fine
- Yellow/Extra Fine

Super
Coarse
Black
特粗砂
黑圈



编号 No. D060S D070S K080S K100S J070S J100S A050S A060S C050S C060S

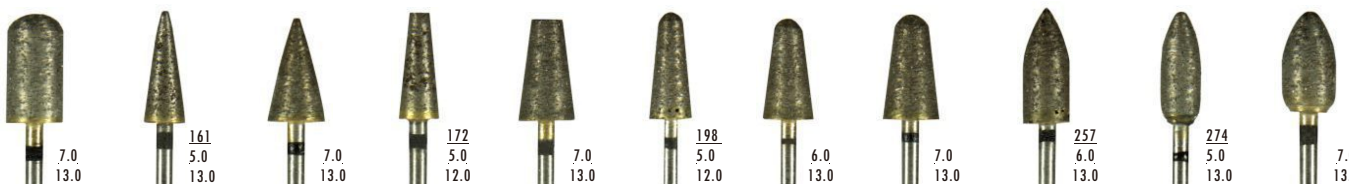

C070S M050S M070S B050S B070S L050S L060S L070S G060S F050S F070S

Coarse
Green
粗砂
绿圈



R050S R070S D050C D060C J060C J070C A037C A050C A060C C037C






C050C C060C B037C B050C B060C L037C L050C L060C N050C G050C F037C

Medium
Blue
标准砂
蓝圈



F050C R037C R050C D031 D037 D050 D060 K031 K037 K050

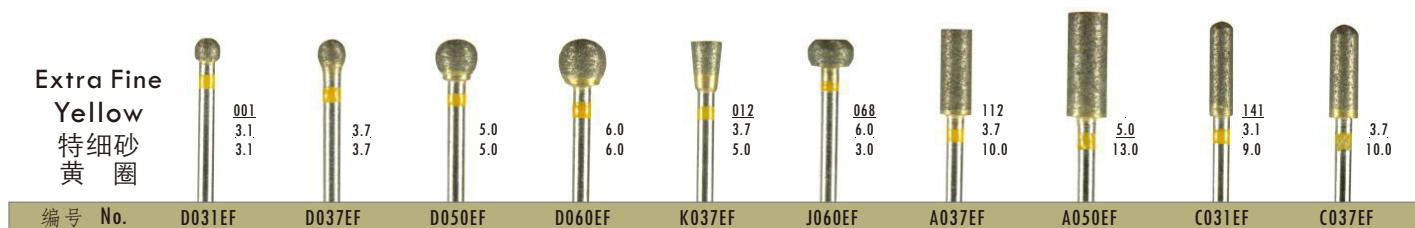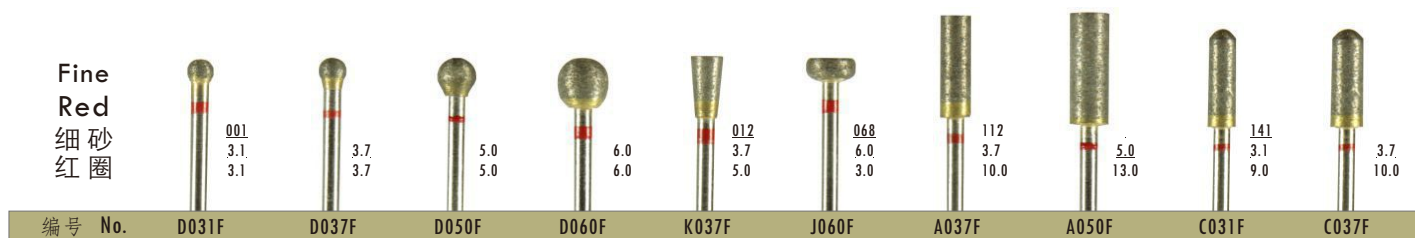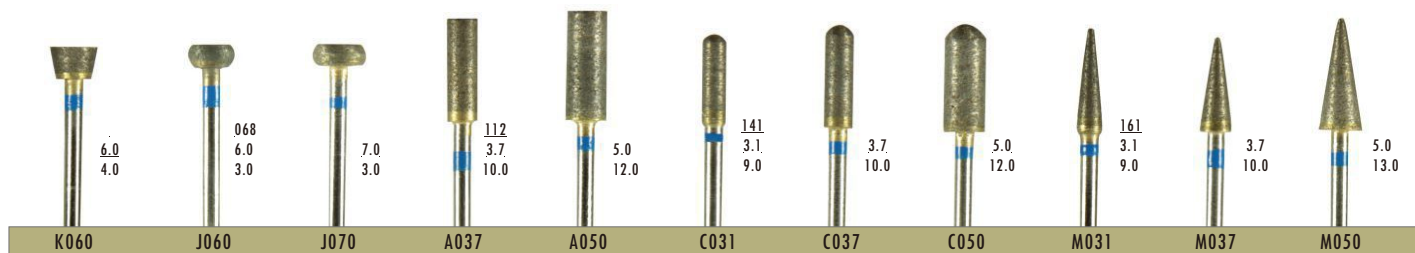


Sintered Diamond Disc

烧结切割片

Grinding or Trimming of
sintered Zirconia or any other
Extra hard workpiece

用于切割烧结后的氧化锆或者
各种超硬材质工件。

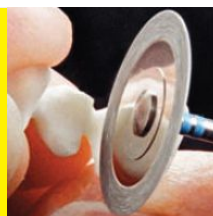

The galvanically sintered edge of
the diamond discs offers excellent
cutting performance.

提供三种不同粒度的砂粒，供您试
用。



70 µm

Item No.



S394/190/0.25mm

S394/220/0.25mm

S394/190/0.30mm

S394/220/0.30mm

Size (mm)

19 x 0.25mm

22 x 0.25mm

19 x 0.30mm

22 x 0.30mm



60 µm

Item No.



S395/190/0.25mm

S395/220/0.25mm

S395/190/0.30mm

S395/220/0.30mm


Size (mm)

19 x 0.25mm

22 x 0.25mm

19 x 0.30mm

22 x 0.30mm

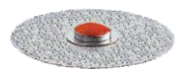

40 µm

Item No.



S396/190/0.25mm

S396/220/0.25mm

S396/190/0.30mm

S396/220/0.30mm

Size (mm)

19 x 0.25mm

22 x 0.25mm

19 x 0.30mm


22 x 0.30mm

Diamond Disc

金刚石切片



How to order these Discs:

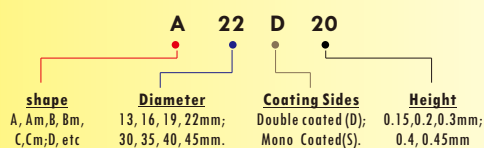

Extra Thin, 0.15mm
Double sided
Extra Fine
超薄, 双面砂



编号 No.	A19D15	A22D15	Am13D15	Am16D15	Am19D15	Am22D15
--------	--------	--------	---------	---------	---------	---------

Thin, 0.20 x 19mm
Double sided
Fine Grit
较薄, 双面砂



编号 No.	A19D20	B19D20	C19D20	N19D20	Am19D20	Bm19D20	Em19D20	Hm19D20
--------	--------	--------	--------	--------	---------	---------	---------	---------

Thin, 0.20 x 22mm
Double sided
Fine Grit
较薄, 双面砂



编号 No.	Am22D20	Bm22D20	Cm22D20	Em22D20	Im22D20	Jm22D20	Hm22D20
--------	---------	---------	---------	---------	---------	---------	---------




A22D20	B22D20	C22D20	E22D20	F22D20	H22D20	J22D20	N22D20
--------	--------	--------	--------	--------	--------	--------	--------

Hard, 0.3 - 0.4mm
Double sided
Medium Grit
厚硬, 双面砂



编号 No.	Am19D30	Am22D30	Am22D40	A19D40	A22D40
	19.0 x 0.3 mm	22.0 x 0.3 mm	22.0 x 0.4 mm	19.0 x 0.4 mm	22.0 x 0.4 mm

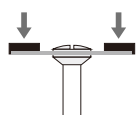


Single-sided
仅单面镀砂



编号 No.	A22S20	Am22S20	Am22S30
	22.0 x 0.2 mm, Fine	22.0 x 0.2 mm	22.0 x 0.3 mm, Medium


Fast Cutting
Plaster Disc ,
0.3mm,
Double sided
石膏模切片
平整、无偏摆

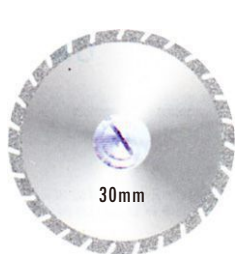

30mm



35mm

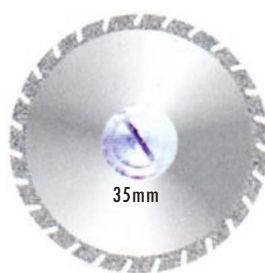

35mm



35mm


编号 No.

C30D30

A35D30

C35D30

J35D30

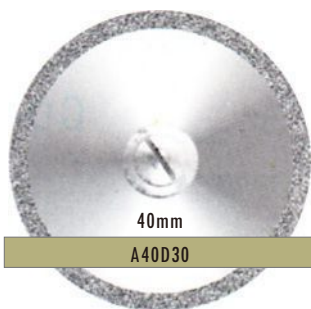

40mm

A40D30



40mm

C40D30




40mm

J40D30




45mm

Cb45D30

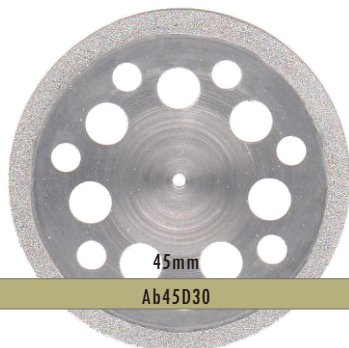

45mm

Ab45D30

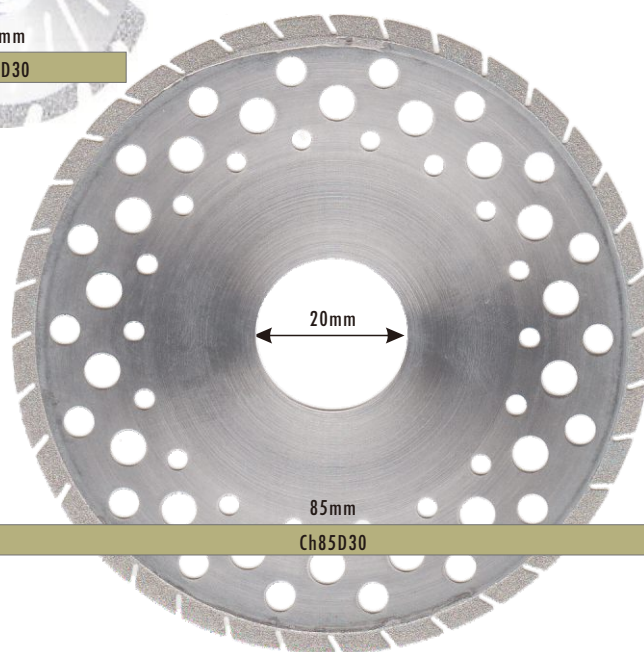

20mm

85mm

Ch85D30

Tungsten Carbide Burs . HP 钨钢车针



Shape

- S1, ISO 001, Round
- S2, ISO 002, Inverted cone
- S21, S31, ISO 107, Cylinder
- S33, ISO 168, Taper
- S31R, ISO 137, Round, Cylinder

										
S1	001 0.5	001 0.6	001 0.7	001 0.8	001 0.9	001 1.0	001 1.2	001 1.4	001 1.6	001 1.8
球型/Round										
编号 No.	1005	1006	1007	1008	1009	1010	1012	1014	1016	1018
参考号 Ref.	HP1/4	HP1/2	HP1s	HP 1	HP 1b	HP 2	HP 3	HP 4	HP 5	HP 6

										
001 2.1	001 2.3	001 2.5	001 2.7	001 3.1	S2	010 0.8	010 1.0	010 1.2	010 1.4	010 1.6
					倒锥/Inv. Cone					
1021	1023	1025	1027	1031		2008	2010	2012	2014	2016
HP 7	HP 8	HP 9	HP 10	HP 12		HP 34	HP 35	HP 36	HP 37	HP 38


										
010 1.8	010 2.1	010 2.3	S21	107 0.8	107 0.9	107 1.0	107 1.2	107 1.4	107 1.6	107 1.8
			圆柱/Cylinder							
2018	2021	2023		21008	21009	21010	21012	21014	21016	21018
HP 39	HP 40	HP 41		HP 55	HP 56	HP 57	HP 58	HP 59	HP 60	HP 61

										
107 2.1	107 2.3	S31	107 0.8	107 0.9	107 1.0	107 1.2	107 1.4	107 1.6	107 1.8	107 2.1
		圆柱/Cylinder								
21021	21023		31008	31009	31010	31012	31014	31016	31018	31021
HP 62	HP 63		HP 555	HP 556	HP 557	HP 558	HP 559	HP 560	HP 561	HP 562

										
107 2.3	S33	168 0.8	168 0.9	168 1.0	168 1.2	168 1.6	168 2.1	168 2.3	S33L	171 0.9
	平锥/Taper								平锥加长/Taper L	
31023		33008	33009	33010	33012	33016	33021	33023		33L009
HP 563		HP 698	HP 699	HP 700	HP 701	HP 702	HP 703	HP 704		HP 699 L

										
171 1.0	171 1.2	171 1.6	171 2.1	S31R						
				圆头柱 Round Cyl.						
33L010	33L012	33L016	33L021		31R010	31R012	31R014			
HP 700 L	HP 701 L	HP 702 L	HP 703 L		HP 1557	HP 1558	HP 1559			

Tungsten Carbide Cutter . HP 钨钢磨头

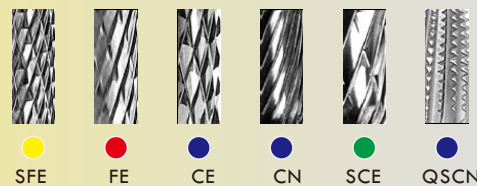

Super Fine
Cross Cut
(SFE)
特细双纹
黄圈
ISO 500.104.110.



编 号 No.	C023SFE	F040SFE	F060SFE	G023SFE	L016SFE	L023SFE	L040SFE	L060SFE	M060SFE	P023SFE
参考号 Ref.	0230.023	7230.040	5430.060	1030.023	1230.016	0130.023	5730.040	5530.060	5930.060	7730.023

Fine
Cross Cut
(FE)
较细双纹
红圈
ISO 500.104.140.



Q023SFE	Q023SFE	C023FE	C060FE	F060FE	G040FE	G060FE	L016FE	L023FE
0830.023	0430.023	0220.023	5820.060	5420.060	6120.040	5320.060	1220.016	0120.023

Medium
Cross Cut
(CE)
标准双纹
蓝圈
ISO 500.104.190.



L040FE	L060FE	M060FE	P023FE	P060FE	Q023FE	A023CE	C023CE	C060CE
5720.040	5520.060	5920.060	7720.023	5120.060	0420.023	0310.023	0210.023	5810.060



D023CE	D030CE	D040CE	D050CE	F040CE	F060CE	G023CE	G040CE	G060CE	L023CE
7110.023	7110.030	7110.040	7110.050	7210.040	5410.060	1010.023	6110.040	5310.060	0110.023

								Medium Plain Cut (CN) 标准单纹 蓝圈 ISO 500.104.175.	
201 040 14.0	201 060 15.0	167 023 14.50	167 040 14.5	161 060 15.0	237 023 5.0	237 060 11.0	292 023 16.0		001 030 3.0
L040CE	L060CE	M023CE	M040CE	M060CE	P023CE	P060CE	Q023CE		D030CN
5710.040	5510.060	—	—	5910.060	7710.023	5110.060	0410.023		7170.030

							Coarse Cross Cut (SCE) 粗糙双纹 绿圈 ISO 500.104.220			
001 040 4.0	001 050 5.0	201 040 14.0	201 060 15.0	161 060 14.0	257 060 12.0	274 060 14.0		201 060 14.0	274 060 14.0	257 060 12.0
D040CN	D050CN	L040CN	L060CN	M060CN	G060CN	F060CN		C060SCE	F060SCE	G060SCE
7170.040	7170.050	5770.040	5570.060	5970.060	5370.060	5470.060		5850.060	5450.060	5350.060

			Super Coarse Cross Cut (XCE) 特粗双纹 黑圈 ISO 500.104.223		Straight Edge with Spiral Cut (QSCN) 直刃螺旋纹 蓝圈 ISO 500.104.176				
201 040 14.0	201 060 14.0	237 060 11.0		263 070 14.0		274 060 15.0	201 023 14.5	200 045 13.0	263 060 12.0
L040SCE	L060SCE	P060SCE		L070XCE		F060QSCN	L023QSCN	L045QSCN	LB060QSCN
5750.040	5550.060	5150.060		6055.070		7275.060	0175.023	5675.045	5275.060

Straight Edge with Cross Cut (GTi) 直刃双纹 黑圈 ISO 500.104.194						Spiral Cut Fine (QF) 较细螺旋纹 红圈 ISO 500.104.134					
		141 023 8.0	194 023 8.0	201 040 14.0	289 023 8.0			274 040 11.5	274 060 15.0	201 023 14.5	201 040 13.0
		CS023GTi	LS023GTi	L040GTi	QS023GTi			F040QF	F060QF	L023QF	L040QF
		1726.023	0726.023	5726.040	0826.023			7224.040	5424.060	0124.023	5724.040

Spiral Cut Extra Fine (QSF) 特细螺旋纹 黄圈 ISO 500.104.137							 	
		141 023 8.0	274 060 14.0	198 023 8.0	201 040 14.0	289 023 8.0		
		CS023QSF	F060QSF	LS023QSF	L040QSF	QS023QSF		
		1725.023	5425.060	0725.023	5725.040	0825.023		
							118.060	
							467.023	

Tungsten Carbide Burs for Milling Machines


钨钢车针(精密研磨仪用)



Precisely match Various System

Extra - long working life & high cutting capacity ;
Used for trimming, polishing implant teeth shoulders ;
Milling precious or Non-precious metals or Titanium .

Medium
Cross Cut
(CM, LM)
标准双纹
ISO 500.103.190

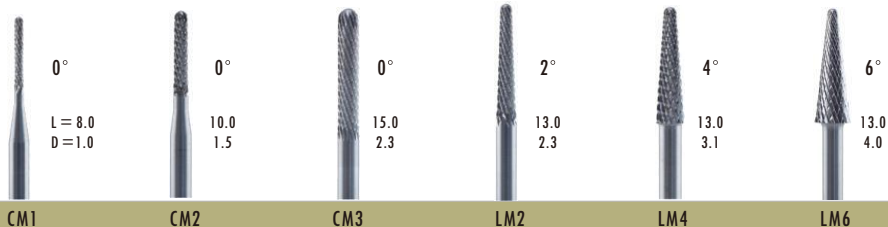


Fine
Cross Cut
(CF, LF)
较细单纹
ISO 500.103.135



Coarse
Cross Cut
(CC, LC)
粗糙双纹
ISO 500.103.220



Max Mill
(CW, LW)
蜡型针
ISO 500.103.217

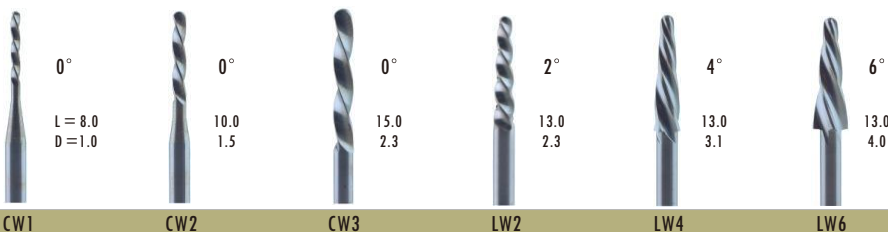

GZH

PRE-SINTERED ZIRCONIA POLISHING BURS

预烧结氧化锆
抛光磨头

- Polishing pre-sintered zirconia simplifies finishing after final sintering
抛光预烧结氧化锆可简化最终烧结后的加工步骤
- No contamination of zirconia, as no colour pigments are used
由于未使用颜色颜料，氧化锆无污染
- Developed for the lower strength of pre-sintered zirconia
为预烧结氧化锆的较低强度而开发的



Set GZHS-01

Pre-sintered Zirconia Finishing

预烧结氧化锆
打磨

NO.	GZH243-5C	GZH275-6C	GZH040-15C
D x H mm	5.5 x 13	6.5 x 13	15 x 0.4

HP Shank, Coarse
HP 长柄, 粗糙

Pre-sintered Zirconia Polishing 2-in-1 Burs

预烧结氧化锆
抛光
二合一磨头

NO.	GZH507-10MF	GZH243-10MF	GZH243-6MF	GZH507-7MF	GZH204-10MF
D x H mm	10 x 22.5	10 x 20	6 x 16	7 x 19	10 x 24

HP Shank, Medium + Fine
HP 长柄, 中等 + 精细

1 STEP
MEDIUM



2 STEP
FINE



10 Pieces

DuCoBur® 陶瓷结合剂金刚石磨头 (全瓷、全锆磨头)

Diamond Stones for Zirconia or AllCeramics
(Ceramic Pregnant Diamond Grinding Burs)



Widely Used for Zirconia, Lithium Disilicate, Ceramics etc .
COMFORTABLE at the minimal pressure;
HEATLESS no need water cooling ;
MORE DURABLE under our suggested PRM

不产热、无需水冷 更耐磨、经济实惠 更锋利、手感好 不绷瓷、更轻松 型号齐、用途广、啥都磨得动

Medium Blue
蓝色、中等粒度
1

Medium Blue 蓝色、中等粒度 1									
Item No.	DC009	DC010	DC010B	DC011	DC014	DC015	DC016	DC017	DC003A
Size (mm)	4 x 4	4 x 7	4 x 11	3 x 7	2.5 x 4.5	3.5 x 10	2.5 x 7	3.5 x 11	6.5 x 2
									
DC001	DC002	DC002B		DC012	DC012B		DC013		DC004A
4 x 13	5 x 13	6.5 x 13		13 x 2	13 x 3.5		13 x 2		16 x 2/0.7
									
DC005A	DC004	DC005		DC006		DC007			
16 x 2	22 x 2/0.7	22 x 2		22 x 4.5		15 x 3.5			

Medium Black
黑色、中等粒度
1

Medium Black 黑色、中等粒度 1							
Item No.	DX001	DX002	DX002B	DX012	DX012B	DX013	DX004A
Size (mm)	4 x 13	5 x 13	6.5 x 13	13 x 2	13 x 3.5	13 x 2	16 x 2/0.7
							
DX005A	DX004	DX005	DX006	DX007			
16 x 2	22 x 2/0.7	22 x 2	22 x 4.5	15 x 3.5			



Medium
Yellow

黄色、中等粒度



1



Medium
Yellow

黄色、中等粒度


1

							
Item No.	DY001	DY002	DY002B	DY012	DY012B	DY013	DY004A
Size (mm)	4 x 13	5 x 13	6.5 x 13	13 x 2	13 x 3.5	13 x 2	16 x 2/0.7



DY005A

16 x 2



DY004

22 x 2/0.7



DY005

22 x 2



DY006

22 x 4.5



DY007

15 x 3.5



Medium
Green

绿色、中等粒度



1



Medium
Green

绿色、中等粒度


1

									
Item No.	DG009	DG010	DG010B	DG011	DG014	DG015	DG016	DG017	DG003A
Size (mm)	4 x 4	4 x 7	4 x 11	3 x 7	2.5 x 4.5	3.5 x 10	2.5 x 7	3.5 x 11	6.5 x 2



DG001

4 x 13



DG002

5 x 13



DG002B

6.5 x 13



DG012

13 x 2



DG012B

13 x 3.5



DG013

13 x 2



DG004A

16 x 2/0.7



DG005A

16 x 2



DG004

22 x 2/0.7



DG005

22 x 2



DG006

22 x 4.5




DG007

15 x 3.5

Extra Fine
Yellow
Gold shank


Suitable to adjusting Lithium Disilicate, All Ceramics, Zirconia of Occlusal, Lingual, Edge, Intracoronar, Proximal etc. Dentists will feel very convenient if they can prepare a few pieces beside the chairs.

专门针对诊所的全瓷、玻璃陶瓷打磨抛光、氧化锆牙调咬合、发育沟、窝沟、边缘、内冠、邻面、舌侧等。黄色、特细粒度。粒度细、不绷瓷。

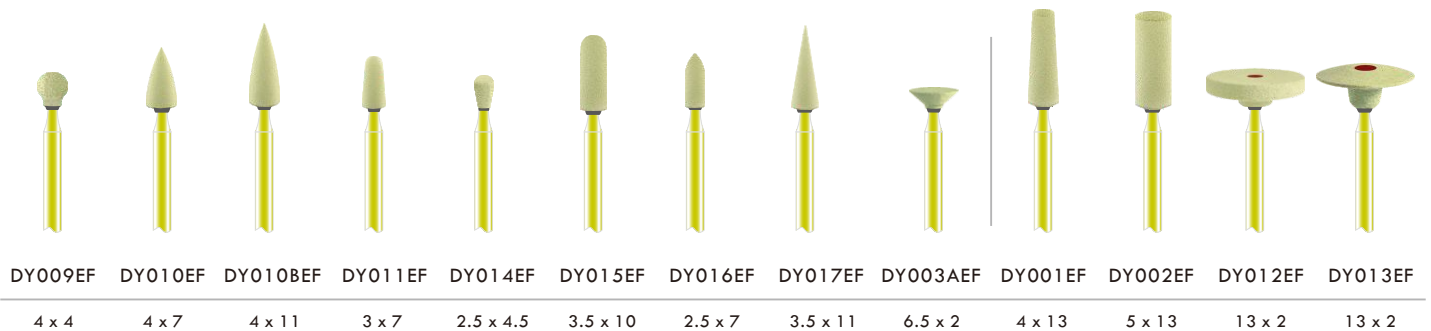

DuCoBur® 豆酷巴

Diamond Polishers for Zirconia or AllCeramics
(Rubber Pregnant Diamond Polishing Burs)

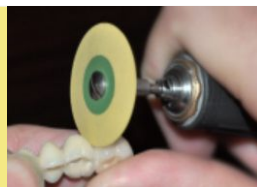

Widely Used for Zirconia, Lithium Disilicate, Ceramics etc.

COMFORTABLE at the minimal pressure;
HEATLESS no need water cooling ;
MORE DURABLE under our suggested PRM

不产热、无需水冷 更耐磨、经济实惠 手感好、很服帖 抛就亮、更轻松 型号齐、用途广、啥都抛得亮

Super Coarse
Black
黑色、特粗粒度



Item No.

D1S

D2AS

D2S

D3S

D3BS

Size (mm)

26 x 2

20 x 2

26 x 2

22 x 3

26 x 2

Coarse
Green
绿色、粗糙粒度



Item No.

D1C

D2AC

D2C

D3C

D3BC

Size (mm)

26 x 2

20 x 2

26 x 2

22 x 3

26 x 2

Medium
Blue

蓝色、中等粒度



Item No.

D1M

D2AM

D2M

D3M

D3BM

Size (mm)

26 x 2

20 x 2

26 x 2

22 x 3

26 x 2



D201M

5 x 16



D203M

6 x 16



D205M

12 x 2



D205BM

12 x 4



D207M

12 x 2.5



D13M

14 x 2



D17M

17 x 2.2



D26M

26 x 2.2

Fine
Red

红色、较细粒度



Item No.

D1F

D2AF

D2F

D3F

D3BF

Size (mm)

26 x 2

20 x 2

26 x 2

22 x 3

26 x 2



D201F

5 x 16



D203F

6 x 16



D205F

12 x 2



D205BF

12 x 4



D207F

12 x 2.5



D13F

14 x 2



D17F

17 x 2.2



D26F

26 x 2.2

Extra Fine
Yellow
黄色、特细粒度



Item No.

D1EF

D2AEF

D2EF

D3EF

D3BEF

Size (mm)

26 x 2

20 x 2

26 x 2

22 x 3

26 x 2



D201EF

5 x 16



D203EF

6 x 16



D205EF

12 x 2



D205BEF

12 x 4



D207EF

12 x 2.5



D13EF

14 x 2



D17EF

17 x 2.2



D26EF

26 x 2.2

Ultra Fine
Grey
灰白色、特细粒度



Item No.

D1UF

D2AUF

D2UF

D3UF

D3BUF

Size (mm)

26 x 2

20 x 2

26 x 2

22 x 3

26 x 2



D201UF

5 x 16



D203UF

6 x 16



D205UF

12 x 2




D205BUF

12 x 4



D207UF

12 x 2.5

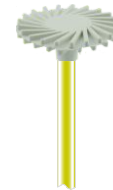

D13UF

14 x 2



D17UF

17 x 2.2



D26UF

26 x 2.2

**Silicon Rubber Wheels ,
Points & Burs
For Polish Ceramics
胶轮、胶粒、硅胶磨头
(烤瓷牙抛光)**




White, Coarse , to adjust and remove Porcelain
白色, 粗糙, 用于烤瓷工件的修整。
Pink, Medium , to Smooth or polish Porcelain,
粉色, 中等, 用于烤瓷毛面平整和预抛光。
Blue, Fine , to high-lustre polishing Porcelain,
蓝色, 较细, 用于烤瓷工件的细抛光。


Coarse
White
白色、粗糙



Item No.

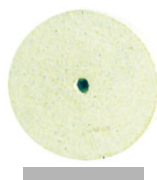

C 104 C

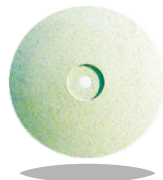

C 107 C



C 101 C



C203C



C201C



C205C



C207C

Size (mm)

22 x 3

22 x 3

6 x 22

6 x 13

5 x 13

12 x 2

13 x 2

Medium
Pink
粉色、中等粒度



Item No.



C 104 M



C 107 M



C 101 M



C203M



C201M



C205M



C207M

Size (mm)

22 x 3

22 x 3

6 x 22

6 x 13

5 x 13


12 x 2

13 x 2

Fine
Blue
蓝色、较细



Item No.

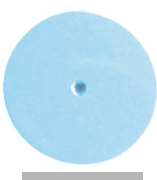

C 104 F



C 107 F



C 101 F



C203F



C201F



C205F



C207F

Size (mm)

22 x 3

22 x 3

6 x 22

6 x 13

5 x 13

12 x 2

13 x 2

**Rubber Wheels , Points & Burs
For Polish Metal or Alloy
胶轮、胶粒、橡胶磨头
(金属、合金抛光)**



Black, Coarse , to adjust Cobalt-Chromium
Alloy and hard metal
黑色, 粗糙, 用于磨削钴铬合金和硬金属。
Brown, Medium , to remove metal and
semi-precious metal
棕色, 中等, 用于烤瓷毛面平整和预抛光。
Green, Fine , to fast smooth scratch or polish metal,
绿色, 较细, 用于快速去除划痕和抛光金属。

Coarse
Black
黑色、较粗



Item No.



A102 C



A103C



A104 C



A203C



A201C



A206C



A207C

Size (mm)

22 x 3

22 x 3

6 x 22


4 x 12

5 x 13

12 x 2

13 x 2

Medium
Brown
棕色、中等粒度
100

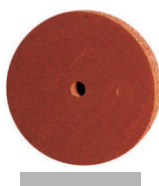

Item No.	A102 M	A103 M	A104 M	A203 M	A201 M	A206 M	A207 M
----------	--------	--------	--------	--------	--------	--------	--------

Size (mm)	22 x 3	22 x 3	6 x 22	4 x 12	5 x 13	12 x 2	13 x 2
-----------	--------	--------	--------	--------	--------	--------	--------

Fine
Green
绿色、较细
100



Item No.	A102 F	A103 F	A104 F	A203 F	A201 F	A206 F	A207 F
----------	--------	--------	--------	--------	--------	--------	--------

Size (mm)	22 x 3	22 x 3	6 x 22	4 x 12	5 x 13	12 x 2	13 x 2
-----------	--------	--------	--------	--------	--------	--------	--------

Silicon Polishing Burs
For Composites, Resins
硅胶磨头
(树脂胶托抛光)



Green, Coarse, to reduce, finish Acrylic, Resin.
绿色, 粗糙, 用于修整和成型树脂胶托。
Gray, Medium, to Smooth Acrylic, Resin.
灰黑色, 中等, 用于平整和预抛光树脂胶托。
Yellow, Fine, to final polish Acrylic, Resin.
黄色, 较细, 用于树脂胶托的细抛光。

Coarse
Green
绿色、粗糙
100



Item No.	Rh 101 C	Rh 103 C	Rh 106 C	R103 C	R106 C	R209 C
----------	----------	----------	----------	--------	--------	--------

Size (mm)	13 x 20	10 x 24	10 x 20	10 x 24	10 x 20	6 x 13
-----------	---------	---------	---------	---------	---------	--------

Medium
Gray
灰黑、中等粒度
100



Item No.	Rh 101 M	Rh 103 M	Rh 106 M	R103 M	R106 M	R209 M
----------	----------	----------	----------	--------	--------	--------

Size (mm)	13 x 20	10 x 24	10 x 20	10 x 24	10 x 20	6 x 13
-----------	---------	---------	---------	---------	---------	--------

Fine
Yellow
黄色、较细
100



Item No.	Rh 101 F	Rh 103 F	Rh 106 F	R103 F	R106 F	R209 F	R202 F
----------	----------	----------	----------	--------	--------	--------	--------

Size (mm)	13 x 20	10 x 24	10 x 20	10 x 24	10 x 20	6 x 13	4 x 12
-----------	---------	---------	---------	---------	---------	--------	--------

Green Stones
Contour Ceramic



绿车石
调烤瓷



Pink Stones
for Metal
粉车石
调金属,合金



White Stone
for Metal
Composite,
白车石
调金属,树脂

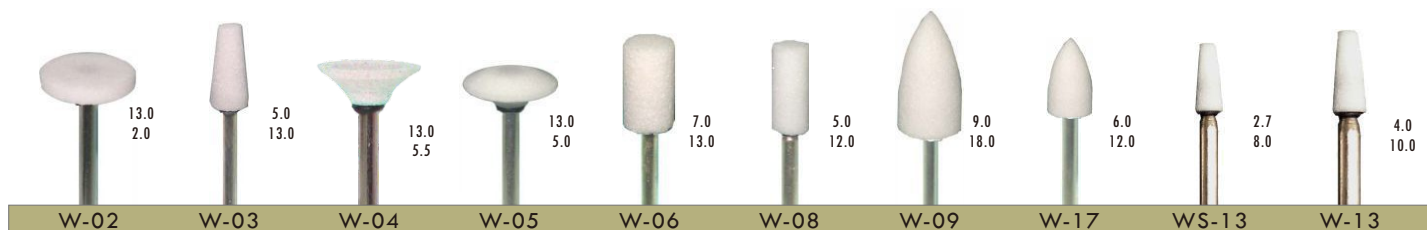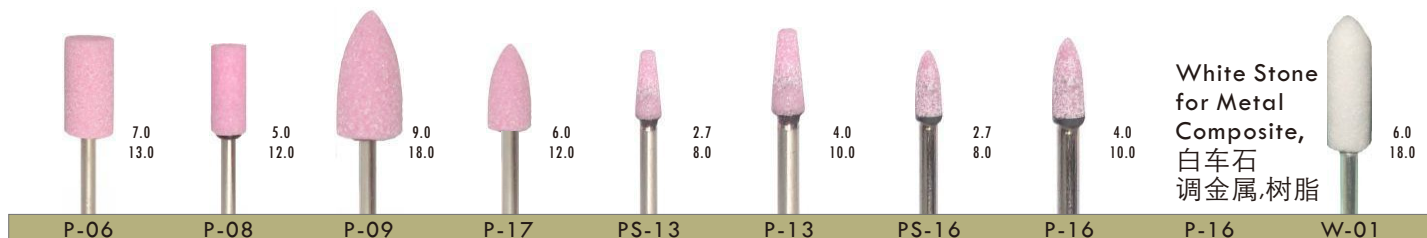

Mounted Sandpaper 砂纸棒



Polishing Brushes 抛光刷

				
百洁布 暗红,中粗	百洁布 红,硬	百洁布 棕,中粗	百洁布 紫,软	百洁布 普通
01MT25	01Mp25	01My25	01Ms25	01M22
Red.Soft	Red.coarse	Brown. Medium	Purple.fine	Blue

						
碳化硅 磨料丝	氧化铝 磨料丝	曲钢丝	曲铜丝	黑猪鬃	灰马毛	灰羊毛
01SiC22	01Al22	01A22	01B22	01D22	01E22	01F22
Silicon Carbide	Al. Oxide	Crimped steel wire	Crimped brass wire	Chungking Bristle	Grey Pony hair	Grey goat hair

					
白羊毛	灰羊毛 +超纤皮	黑猪鬃 19mm	灰马毛	灰羊毛	白羊毛
01G22	01Fr22	01D19	01E19	01F19	01G19
White goat hair	Goat hair / chamois	Chungking Bristle	Grey Pony hair	Grey goat hair	White goat hair

							
黄棉布	白棉布	橘红 羊皮	双面绒 羊皮	咖啡 超纤	白厚 超纤	羊绒	白棉线
01P22	01Q22	01R22	01Ry22	01R22.2	01Ry22.2	01S22	01T22
Yellow cloth	White cloth	Chamois leather	Sheep Leather	Brown Micro-Fiber	White Microfiber	Cashmere	Cotton yarn



						
黑猪鬃	灰马毛	白羊毛	曲钢丝	曲铜丝	黑猪鬃	灰羊毛
02D22	02E22	02G22	04A12	04B12	04D12	04F12
Chungking Bristle	Pony hair	White goat hair	Steel	Brass	Bristle	Goat Hair

						
碳化硅 磨料丝	氧化铝 磨料丝	曲钢丝	曲铜丝	黑猪鬃	灰马毛	灰羊毛
05SiC05	05Al05	05A05	05B05	05D05	05E05	05F05
SiC	Al. Oxide	Steel	Brass	Bristle	Pony Hair	Goat Hair

Polishing brush & Wheel in Bristle, Cotton Cloth 抛光刷



These brushes are made of bristle hair or White, yellow pure cotton cloth, assembled on polishing latch to finish or polish Plaster Moulds, Acrylic
台式抛光机配用的抛光刷/布轮,用于抛光塑料,树脂,模型等。

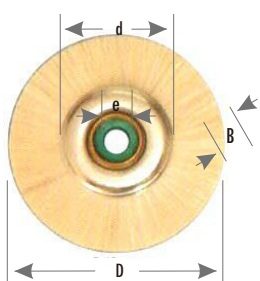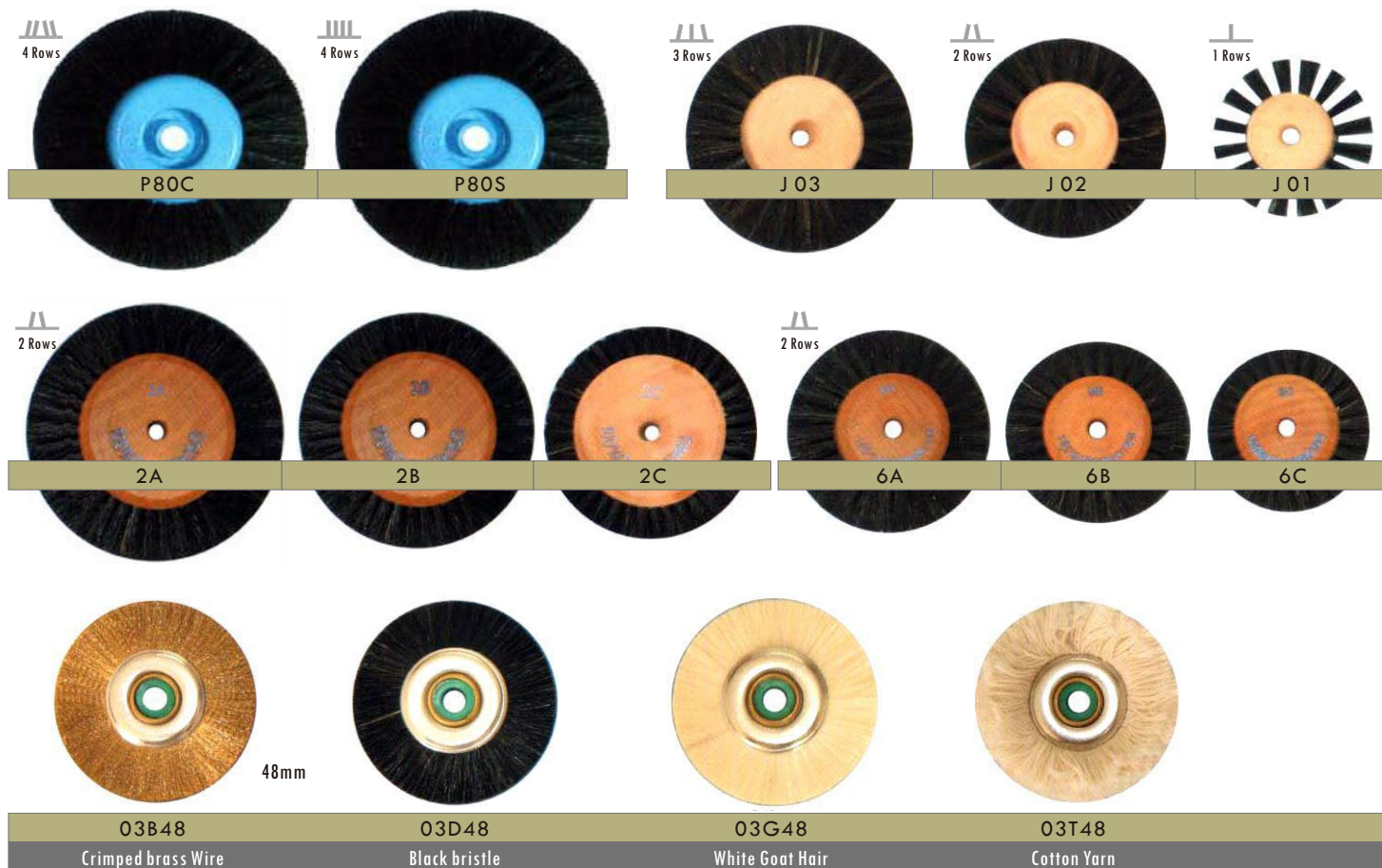

Item	2A	2B	2C	6A	6B	6C	J03	J02	J01	P80	03B48
D	79	73	68	64	60	54	70	58	51	80	48
d	48	48	48	35	35	35	35	30	25	40	25
e	6.5	6.5	6.5	6.5	6.5	6.5	8.5/5	8/5	7/5	6.5	6.35
B	16	16	16	16	16	16	24	16.5	15	36	9
Rows	2	2	2	2	2	2	3	2	1	4	1



Converging




Straight


Felt Bobs & Burs 毛毡磨头




03N25	03Ns25	3P	2P	1P	1R
25 x 5 mm, 硬质	25 x 5 mm, 软质	15 x 25 mm	20 x 40 mm	22 x 50 mm	22 x 50 mm

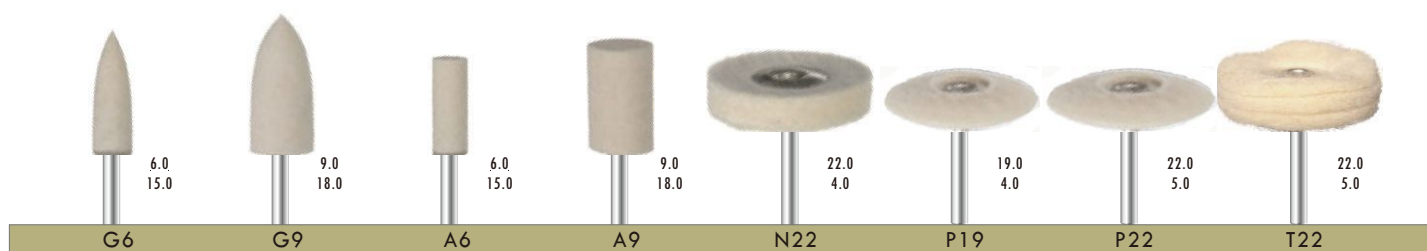

Mandrels Stainless Reinforced 不锈钢夹针

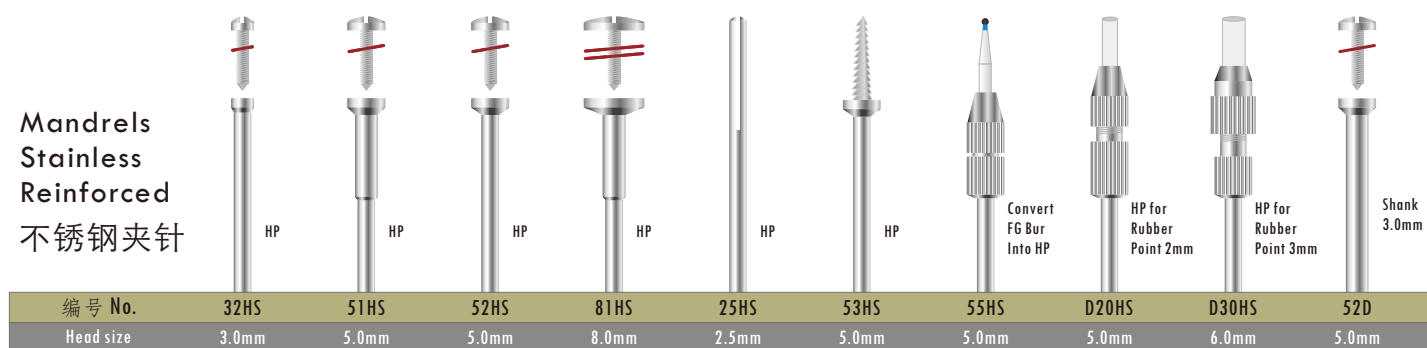

Mandrels Steel 普通夹针

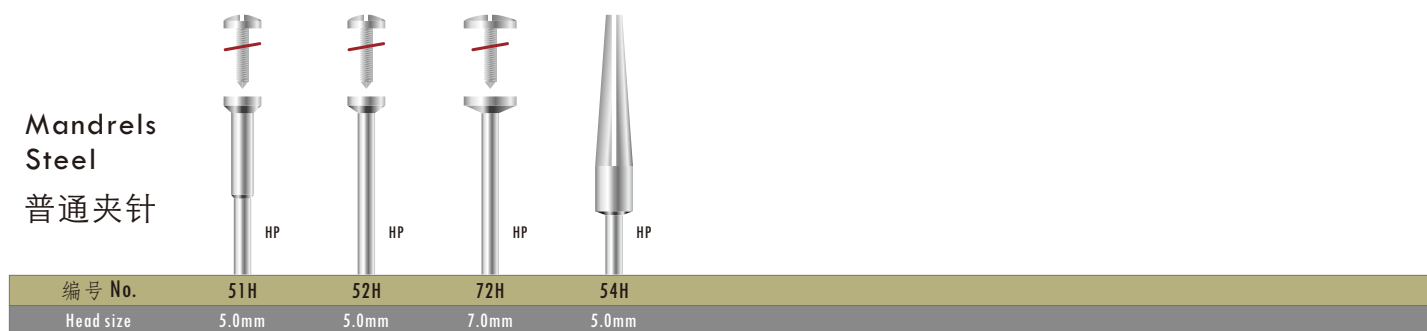

Hard Z Cutter 硬锆磨头


Made of Nano-grade Zirconia, used to grinding, finishing or polishing Ceramics and Acrylic, Resin, Zirconia etc. Also applied to Soft materials such as wood, Ceramic or stones etc

用于打磨,成型,抛光烤瓷,树脂,氧化锆(烧结前)。特殊的纳米级氧化锆胚体精制而成。



Sharpening Stones 磨刀石

Adjust to minimum speed;
Touch the stones softly.
Using water will be better

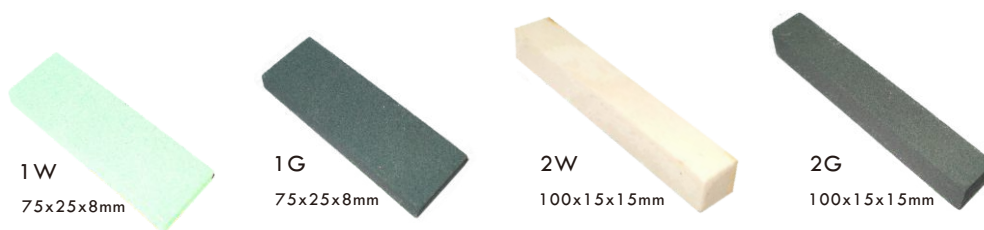

Cleaning Brushes By Hand

清洁刷(手用)



Stainless Steel Disc 不锈钢小锯片

for Soft Plastic, Acrylic, Wood etc.




Stone Disc 刚玉车片



3Z09.1


Zirconia Luster Polishing Kit
氧化锆抛光套装

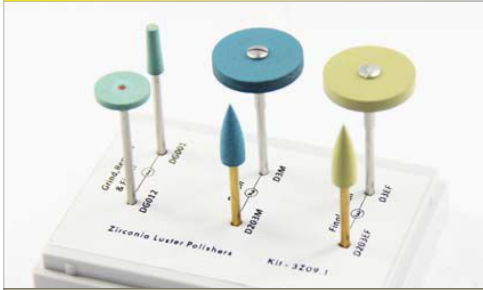

- + 2 pieces of Diamond Stones DG001, DG012
- + 2 pieces of Diamond Polishers, Blue, D203M, D3M
- + 2 pieces of Diamond Polishers, Yellow, D203EF, D3EF

3Z00.2


Grinding and polishing Kit of Zirconia
氧化锆打磨抛光套装

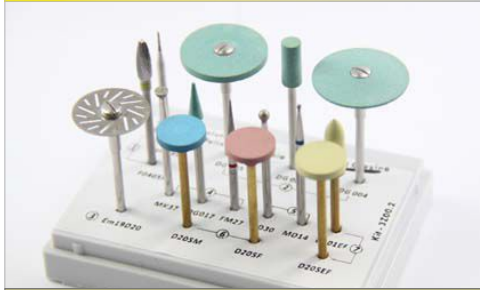

- + 2 pieces of Tungsten carbide Burs HP 699, F040SFE
- + 1 set of Diamond Disc Em19D20 Mounted
- + 4 pieces of HP Diamond Burs: MK37, FM27, FD30, MD14
- + 4 pieces of Diamond Stones: DG017, DG002, DG004, DG005
- + 4 pieces of Diamond Polishers: D205M, D205F, D205EF, D201EF

3Z00.3


Grinding and Polishing Kit for Zirconia
氧化锆打磨抛光套装

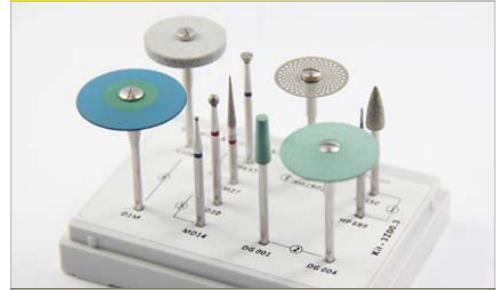

- + 1 piece HP Tungsten Carbide Burs HP 699
- + 5 pieces Diamond Burs: MG50, MD14, FD30, FM27, MK37;
- + 1 set of Diamond Disc Hm19D20 Mounted
- + 1 piece of Diamond Polisher D1M
- + 1 set of Silicon polisher C104C Mounted

3L04.2


Grinding and Polishing Kit for Lithium Disilicate
玻璃陶瓷打磨抛光套装

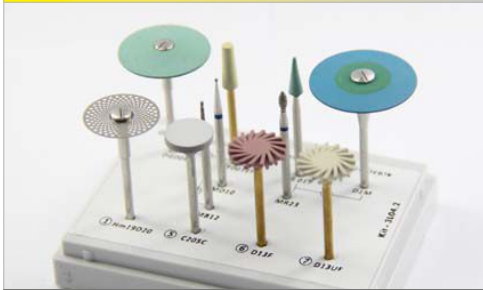

- + 1 set of Diamond Disc HM19D20, Mounted,
- + 3 pieces of Diamond Stones DG017, DG004, DY001EF
- + 3 pieces of HP Diamond Burs MD10, MB12, MR23
- + 3 pieces of Diamond Polisher D1M, D13F, D13UF
- + 1 piece of Silicon Polisher C205C

3L04.3


Grinding and Polishing Kit for Lithium Disilicate
玻璃陶瓷打磨抛光套装

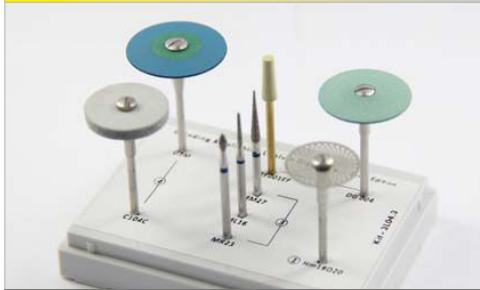

- + 1 set of Diamond Disc HM19D20, Mounted,
- + 2 pieces of Diamond Stones DG004, DY001EF,
- + 3 pieces of HP Diamond Burs MR23, ML16, MM27,
- + 1 piece of Diamond Polisher D1M
- + 1 piece of Silicon Polisher C205C

3MC02.3

Adjusting and Polishing Kit for Ceramics
烤瓷工件打磨抛光套装

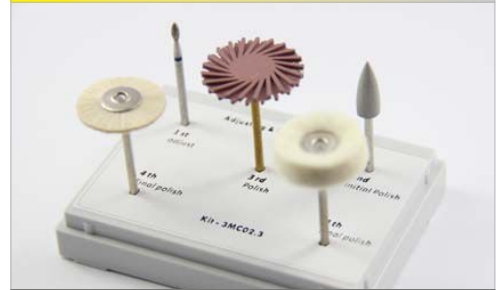

- + 1 piece of HP Diamond Burs MR23
- + 1 piece of HP Silicon Polisher C203C
- + 1 piece of Diamond Polisher D26f
- + 2 pieces of Polishing Brushes 01G22, N22

3MR09.1

Adjust & Polish Invisible Braces
隐形牙套调整抛光套装



- + 1 piece of Tungsten Carbide Burs AL02354
- + 1 set of Diamond Disc C19D20, Mounted
- + 2 pieces of Polishing Brushes 01Ms25, 01S22

3A00.1

Popular HP Kit for grinding or polish Metal or Alloy
常用打磨抛光金属、合金的套装



* Universally grind, remove, polish and finish of Alloy, metal, Jewellery, stones etc.
It widely applied to Laboratory, DIY or Hobbies.
+ 5 pieces of HP sintered diamond burs,
+ 4 pieces of HP Tungsten Carbide Cutters
+ 10 pieces of HP Tungsten Carbide burs

+ 4 pieces of White stones
+ 2 pieces of Pink stones
+ 10 pieces of HP Silicon Polishers
+ 3 pieces of Rubber wheel mounted
+ 1 pieces of Woolen Felt wheel mounted,
+ 5 pieces of HP Polishing brushes

3MC00.1

Popular HP Kit for grinding or polish Ceramics / Porcelain
常用打磨抛光陶瓷的套装



* Universally grind, remove, polish and finish of Ceramics or Porcelains etc.
It widely applied to Laboratory, DIY or Hobbies.
+ 10 pieces of HP Diamond burs,
+ 4 pieces of Diamond Discs mounted
+ 6 pieces of HP Green Stones

+ 3 pieces of Silicon Rubber wheel mounted
+ 8 pieces of Silicon Rubber Polishers
+ 4 pieces of Polishing brushes

3LR00.1

Popular HP Kit for grinding or polish Acrylic, Resin
常用打磨抛光树脂胶托的套装



* Universally grind, remove, polish and finish of Acrylic, Resin etc.
It widely applied to Laboratory, DIY or Hobbies.
+ 2 pieces of HP Tungsten Carbide Cutters
+ 32 pieces of HP Silicon rubber Polishers

+ 10 pieces of Polishing brushes
+ 3 pieces of Woolen Felt polishing brushes

3Z00.1

Popular HP Kit for grinding or polish Zirconia, Lithium Disilicate
常用打磨抛光氧化锆、玻璃陶瓷的套装



* Universally grind, polish and finish of Zirconia and Lithium Disilicate etc.
It widely applied to Laboratory, DIY or Hobbies.
+ 4 pieces of HP Diamond Burs
+ 2 pieces of Diamond Discs Mounted

+ 6 pieces of Diamond Stones
+ 4 pieces of Diamond polishing wheels
+ 2 pieces of Diamond polishing burs
+ 5 pieces of Diamond Twist polishers
+ 1 pieces of Polishing brush.

ZS100

Hot HP Mounted stones mixed
车石套装, 常用。




* 100 Pcs, Hot Mounted stones in Green, Pink, White, applied to Alloy, metal, ceramics, Jewellery, stones, hard wood etc.
+ 30 pieces of HP Green Stones
+ 30 pieces of HP Pink Stones

+ 30 pieces of White Stones
+ 10 pieces of Brown Stones

ZD 101

Popular HP diamond burs KIT
常用金刚石磨头套装

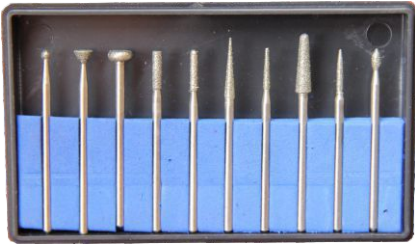


MD23, MK47, MJ60, MA23, MC18, MM27, MB21, ML40, MG16, MR23

10 Hot Spec. of HP Diamond, universally grind, remove, polish and finish of Alloy, metal, Jewellery, stones etc. It widely applied to Laboratory, DIY or Hobbies.

用途: 选取10支最常用的金刚石磨头, 形状全, 大小适中, 可以用来处理陶瓷工件。

ZD 101F

Popular HP diamond burs KIT
常用金刚石磨头套装

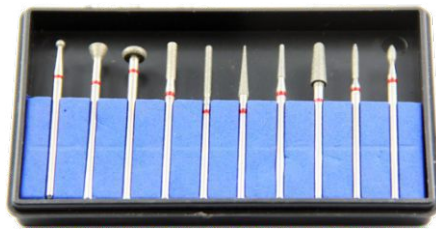


FD23, FK47, FJ60, FA23, FC18, FM27, FB21, FL40, FG16, FR23

10 Hot Spec. of HP Diamond of Fine grit, universally grind, remove, polish and finish of Zirconia, Lithium Disilicate, Porcelain etc. It widely applied to Laboratory, DIY or Hobbies.

用途: 选取10支最常用的金刚石磨头, 形状全, 大小适中, 可以用来处理氧化锆、玻璃陶瓷, 烤瓷等工件。

ZD 102

Small & Tiny Cut, HP diamond burs KIT
金刚石磨头, 小号套装

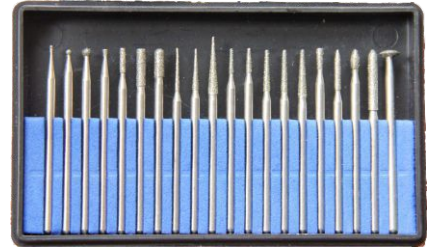

MD10, MD14, MD23, MK16, MA16, MA23, MC25, MM12, MM18, MM27, MB12, MB16, MB21, ML16, ML23, MN18, MG16, MR23, MQ23, MU50

20 Spec. small or sharp burs. To create tiny hole, slim curves or tiny, deep grind smartly. Used to Alloy, metal, ceramics, Jewellery, stones, hard wood etc. It widely applied to Laboratory, DIY or Hobbies.

用途: 20支金刚石磨头, 用来开小洞, 细边, 细槽, 或者小的磨削。很轻巧, 灵活。适用于各种金属, 合金, 陶瓷, 珠宝, 宝石, 木头等。

ZD 103

Fast Trimming, HP diamond burs KIT
金刚石磨头套装, 快速磨削




CQ23, CL40, SD35, SA50, SC40, SP70, SG50, SV70, SV60, SF60

20 Spec. HP Diamond for fast trimming. Grinding of Plastic, Jewellery, stones, Hard wood, etc. Save time when removing and work efficiently

用途: 10支常用的金刚石磨头, 头部相对比较大, 可以快速磨削塑料, 木头, 陶瓷和石头等等, 适用于牙科, 珠宝加工和DIY爱好者。

ZD 104

Hot Fully-Sintered diamond KITs
全砂烧结磨头套装, 常用

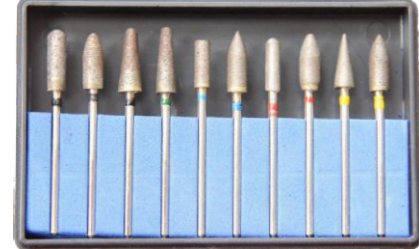


C 050S, F050S, L050S, B050C, A037, G050, C037F, F050F, M050F, G050F

10 spec. Hot Fully-sintered Diamond HP burs, universal grind, remove, polish and finish of Alloy, metal, Jewellery, stones etc. It widely applied to Lab, DIY or Hobbies.

用途: 10个常用全砂烧结磨头, 5种不同粗细。广泛适用于研磨, 抛光金属, 合金, 珠宝, 宝石。经久耐用, 方便实惠。

ZD 105

Fast Trimming Fully-Sintered diamond KIT
全砂烧结磨头套装, 大磨削

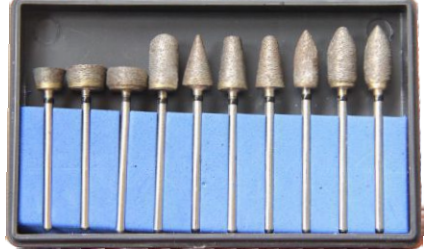


K080S, K100S, J100S, C070S, M070S, B070S, L070S, G060S, F070S, R070S

If workpieces are extremely hard and waste too much time when grinding, this fast trimming kit may help to work efficiently. Applied to Alloy, metal, Jewellery, stones, hard wood, etc.

用途: 10个大头的全砂烧结磨头, 用于粗磨, 去块, 效率高, 经济实惠。

ZT 101

Popular HP Tungsten carbide burs KIT
常用低速钨钢针

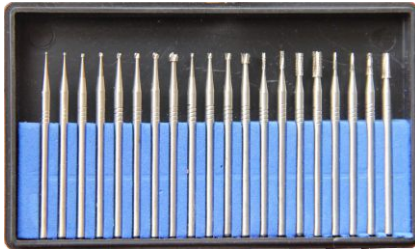

HP 1/4, HP 1/2, HP 1, HP 2, HP 3, HP 4, HP 5, HP 8, HP 34, HP 36, HP 38, HP 41, HP 556, HP 557, HP 560, HP 563, HP 699, HP 701, HP 702, HP 704

20 Spec. HP Tungsten Carbide burs are used to create tiny hole, slim curves or tiny, deep grind very smartly. It is applied to Alloy, metal, ceramics, Jewellery, stones, hard wood etc.

用途: 20个型号钨钢车针, 用于雕牙花, 钻小孔, 洞, 开细槽, 细线等, 金属, 合金, 宝石, 坚木都可以用。

ZT 102

Hot Tungsten carbide cutters KIT
钨钢磨头套装, 常用



D050SP, L040SC, L060SC, F060SC, B023SC, L023SC, L040FC, C023FC, L040EFC, F060CC

10 spec. Tungsten carbide Cutter, universal grind, remove, polish and finish of Acrylic, Resin, stones, soft metal or wood.

用途: 10个常用钨钢磨头, 形状全, 大小适中。主要用来研磨, 抛光塑料, 树脂, 软金属, 石头等, 方便实惠。

ZF 101

Acrylic Polishing Kit
树脂基托抛光套装



F060SC, L040SC, U12 x 2, U13 x 2, U16 x 2

Simple and straight way to adjust Acrylic without the problem associated with pumice.

用途: 2个型号的钨钢磨头快速磨削, 三个不同粒度的硅胶磨头依次进行抛光。

ZA 101

Non-precise metal, Alloy Adjust and polish
非贵金属, 铬钴调整与抛光



W-01, W-09, W-08, W-16, W-03, W-06
A201C, A203C, A206C, A201F, A203F

To adjust and polish Cobalt-chromium, nickel chromium Alloy. White stones for adjusting or grinding, Rubber burs for polishing.

用途: 试用于铬钴、镍铬合金等非贵金属的研磨和抛光 6支白玉石用来研磨和调整, 5支硅胶磨头用来抛光。

ZP 100

Sandpaper Polishing burs Kit
砂纸棒套装




240 #, 400 #, 600 #, 800 #, 1000 #, 1500 #, 2000 #, 2500 #, 3000 #, 5000 #

10 different Grit Sandpaper burs are used for polishing Alloy, metal, Jewellery stones, hard woods etc. Working step by step, mirror surface in coming soon, water cooling or dry are acceptable.

用途: 10个不同粒度的砂纸棒, 从240目到5000目不等, 请根据需要, 选取合适的粒度, 耐心细致地使用。

Z 2001

Popular Felt Polishing bur Kit
毛毡抛光磨头套装

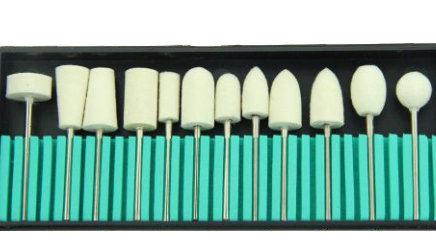

12 pieces of High quality Felt polishing burs are widely used to mould, Alloy, metal, Jewellery, stones, Hard woods etc.

12 pieces of High quality Felt polishing burs are widely used to mould, Alloy, metal, Jewellery, stones, Hard woods etc.

用途: 12支常用的毛毡磨头, 广泛适用于各种模具, 金属, 陶瓷, 珠宝, 宝石和木头等。

Samples Collection

样品册

- Flip-Cart Plastic Pages ;
- Smart, Easily fill all burs ;
- Economic and Recyclable ;
- Nice-looking for showing to Customers ;
- Convenient to take to anywhere .



No. 171B5

Empty Collection Wallet
样品本(不带车针)



4 Colors Available ,
100 holes per page,
5 pages per set,
Smart to HP Burs



No. 1704D

Sintered Diamond Burs Collection
全砂烧结磨头样品本



115 specifications ,
Fully Sintered Diamond,
HP shanks



No. 1705D

HP Diamond Burs, Mounted Stones
金刚石磨头, 车石样品本



More than 100 specifications ,
Electroplating Diamond burs
HP Shanks



No. 1707T

Tungsten carbide Burs Collection
高品质钨钢磨头样品本



More than 100 specifications ,
High Quality TC burs,
HP Shanks



No. 1708P

Polishing Burs, Mandrels Collection
抛光磨头, 夹针样品本



More than 100 specifications ,
Silicon, Rubber, Diamond burs
HP Shanks



隐形义齿胶

FlexSoft Partial Dentures



The alternative to conventional Partial Dentures ;
Flexible . Soft. Easy-to-finish.Easy-to-polish
弹性好，易成型，易抛亮，美观大方。



编号 No.	规格 Spec.
D2	Crystal ; 18 gs
C1	Ligth Pink ; 9 gs
C2	Ligth Pink ; 18 gs
B1	Standard Pink ; 9 gs
B2	Standard Pink ; 18 gs
A1	Deep Pink ; 9 gs
A2	Deep Pink ; 18 gs
For bulk pak of 1 kilo.	
C3	Ligth Pink ; 1000 gs
B3	Standard Pink ; 1000 gs
A3	Deep Pink ; 1000 gs



Instructions for Use

Suitable for: Partial and complete dentures, occlusal splints and nightguards.

Not Suitable for: Small unilateral restorations.

Warnings: Use eye protection and heat-resistant gloves when injecting.

Use proper ventilation, vacuum system and masks when grinding prosthodontic resins.

Precautions: Cartridge must be injected within 1 hour of opening the foil pouch.

I. PREPARING THE MASTER MODEL FOR DUPLICATION

All undercuts that will affect the path of insertion should be blocked-out with wax prior to duplication. Do not relieve necklines of teeth or areas desirable for retention. For free-end saddles, it is not necessary to remove the distal undercut of the last natural tooth, unless the area is too deep. For free-end saddle cases where anterior teeth are also being replaced, undercuts in the anterior should also be blocked out.

II. SET-UP

Duplicate blocked-out master model. Articulate the duplicate model and the master model (you will need the master model articulated at finishing steps to check occlusion). Set up the case on the duplicate model using only wax (no base plate).

Allow 1 mm of thickness under the teeth to flow to the labial side and into the retention diatorics in the teeth. Use a small round carbide bur to create mechanical retention in the teeth, from mesial to distal and connected to the center diatoric hole.

III. WAXING

Important: Wax the case with anatomical features to reduce grinding during the finishing stage.

Recommended thickness of wax-up:

1. Palatal thickness 1.2 - 1.5 mm
2. Buccal & Labial Flange 1.2 - 2.0 mm
3. Clasps 1.0 - 1.5 mm
4. Lingual Flange 1.2 - 2.0 mm
5. Lingual Major Connector (Mandibular) 1.7 - 2.0 mm

IV. INVESTING AND SPRUING

Invest the case in the center of the flask, aiming the posterior portion toward the sprue opening. After the stone (or plaster) has set, use two sprues from the flask opening to the most posterior portion of the wax on both sides (see diagram). Sprues should be approximately 7 mm. width. Do not create sharp angles with the sprues, as the heated resin must flow in a straight or slightly curved path. Apply a thin layer of petroleum jelly on the entire stone surface. Then invest top half of flask.

V. REMOVING WAX

After the stone has set, place flask in boiling water for 8 minutes. Open the flask and wash out the wax, making sure all wax residues are eliminated. Allow the flask to cool. Apply thin layers of separator on the entire surface of the stone and allow it to dry completely. Do not apply separator on teeth. When dry, close flask.

VI. INJECTING

Wear safety goggles and heat-protective gloves for this procedure.

Spray inside of the cartridge sleeve with mold release spray and place it into furnace. Turn on the furnace and allow it to reach approximately 550°F (290°C), preheating cartridge sleeve for approximately 15 minutes.

Remove cartridge from foil-lined bag (do not remove cartridge from bag until this step). Spray the cartridge and bronze disk with mold release spray. Insert the cartridge into the cartridge sleeve with the cap end first. Insert the bronze disc into the cartridge sleeve against the end of the cartridge. Use the thick bronze disc for medium tubes and the thin bronze disc for the large tubes. Heat cartridge for 12 minutes.

Position the flask with the sprue end up, directly under the shaft of the press. After 12 minutes, remove the cartridge sleeve from the furnace. Keep the cartridge sleeve in a horizontal position until resting on top of the flask, and then turn the cartridge sleeve into a vertical position. Immediately use the press to inject the heated cartridge into the flask using a steady fast motion. Continue applying pressure until tube bursts and continue with maximum downward pressure for 3 minutes.

Remove the cartridge sleeve from the flask with a slight twisting motion. After removing the cartridge sleeve, let the flask cool for at least 30 minutes. Do not put flask into a cold water bath. Use a plaster nipper to de-flask. Do not use a hammer. Remove the compressed cartridge from the cartridge sleeve with a knock-out punch. Important: Cartridge sleeve must be pre-heated for 10 minutes BEFORE the next injection.

VII. FINISHING AND POLISHING

Remove appliance from the model. Cut off the sprues with a cut-off wheel. Trim the edges of appliance with a coarse grinding wheel. Avoid overheating. Use rubber wheels to smooth scratches. Place the appliance on the articulated master model to check the fit, retention, thickness and occlusion. Clean the tissue surface with a small bristle brush.

To polish, use a coarse dental abrasive, such as pumice, with a wet rag wheel. Apply final shine with a hi-shine compound on low speed with a final buff wheel. Be careful not to overheat appliance.




铝桶

Empty Aluminum Cartridge

编号 No.	规格 D x H, mm
25085	25 x 85 mm
25100	25 x 100 mm
25CAP	24.2 x 16 mm, Cap
22065	22 x 65 mm
22085	22 x 85 mm
22100	22 x 100 mm
22120	22 x 120 mm
22CAP	21.2 x 16 mm, Cap



BurDental 博登科技

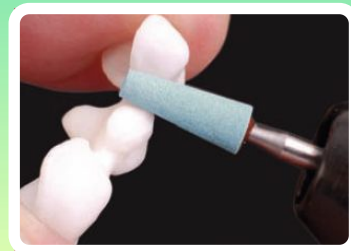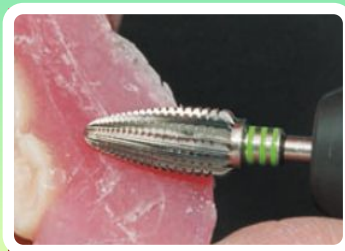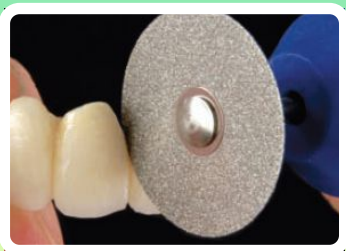

博登科技有限公司

B&D Technologies Co.,Ltd

地址: 广东省广州市番禺区洛浦街恒达产业园A-522

Address : A-522, Heng Da Industrial Park, Luo Pu Street, Pan Yu District, 510000,
GuangZhou, GuangDong Province, China

网站(Web): www.burdental.com ; 邮箱(Email): info@burdental.com

手机(Mobile, WhatsApp) : 0086 1866 6599 173