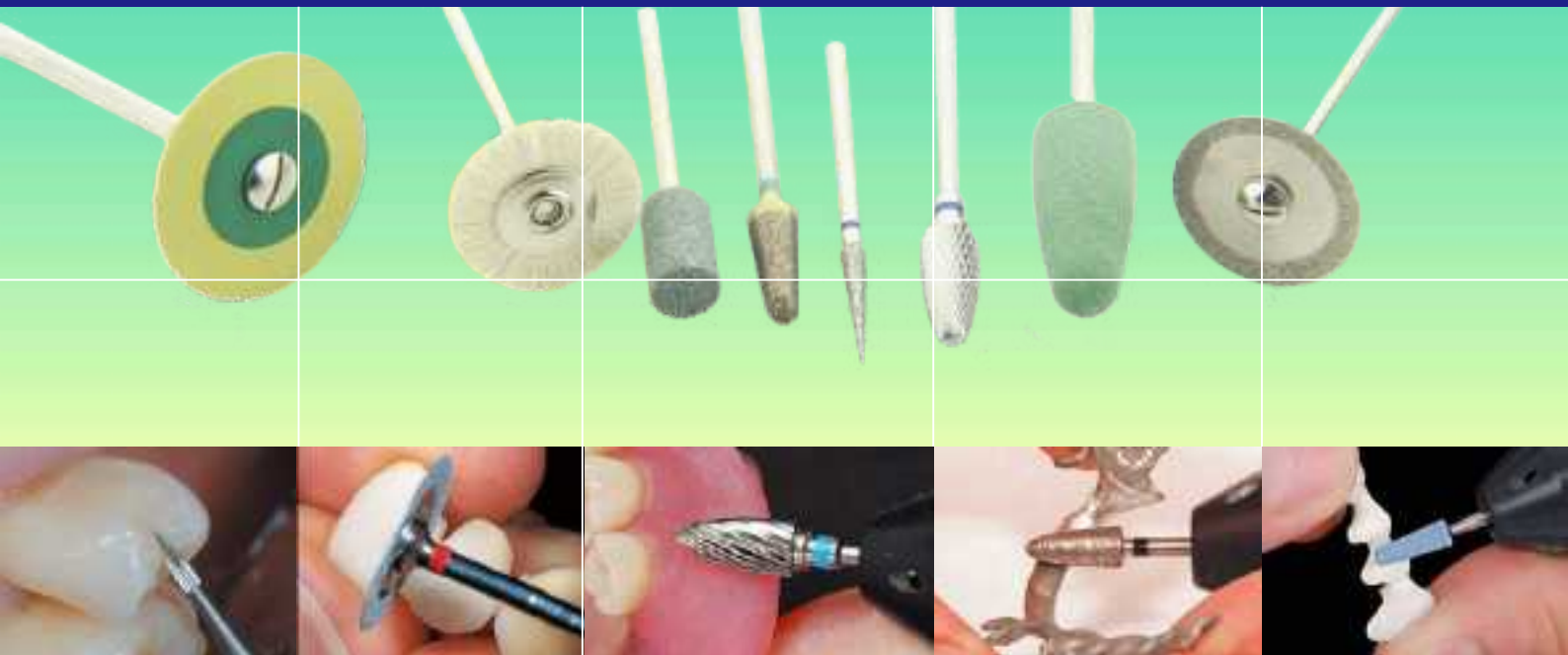




BurDental



Rotary things

Vol. 7

About Us

B&D Technologies are specialized in innovative dental rotary instruments.

We provide all range of rotary instruments from Diamond, Tungsten carbide, Rubber to stone tools.

Our products include Diamond Burs, Diamond Discs, Tungsten Carbide burs, Flexible Polishers, Polishing Brushes, Mounted Stones, Mandrels and other accessories.

B&D diamond burs are made from the finest grade surgical steel and are multilayer coated with premium quality diamonds. Our burs are **ISO13485** and **CE** certificated.

You receive diamond burs that cut quickly and smoothly. All B&D diamond dental burs are guaranteed against defects in materials and workmanship.

We have experiences of more than 15 years to research and develop abrasive technology for dental, Jewelry and others. Providing excellent cutting, grinding, polishing and finishing tools are our mission.

B&D provides excellent service of communication, logistics, training and technical supports. We always deliver good quality tools to customer's hands in very short time and provide professional after-sale support and services. You will feel extremely relaxed and comfortable to choose our tools.

Under the backgrounds of globalization, we develop local and oversea markets actively and grow up very fast. We trust B&D have been your best friend in dentistry.



ISO 13485



CE 0197



Established in 2007

Content 目录

Clinic Part 临床类

Page 3 - 24 第3至24页

Diamond burs 金刚石车针

Page 4 - 8 第4至8页

Diamond Stones Diamond Polishers

陶瓷金刚车石、橡胶金刚石磨头

Page 13 第13页

Cleaning Burs after Debonding Silicon Polishers

(脱胶后) 洁牙磨头
硅胶抛光磨头

Page 14-15 第14-15页

White Stones

Green Stones

白车石、绿车石

Page 15 第15页

Polishing Brushes 抛光刷

Page 16 第16页

Kit 套装

Page 18 - 23 第18至23页

Tungsten Carbide burs

钨钢车针

Page 9 - 12 第9至12页



Clinic

Mandrels, Blocks, Sterilisable

夹针、消毒盒

Page 17 第17页

Samples Collection

样品册

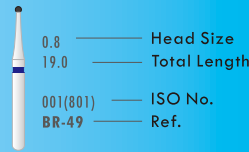
Page 24 第24页

Diamond Bur

金刚砂车针

Advantages





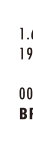
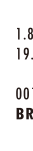
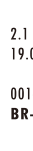
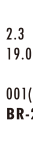
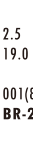
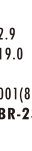





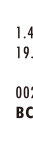

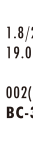
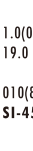





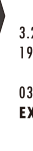
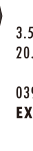

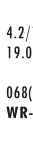









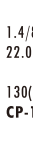









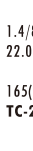











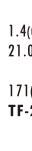
1. Strictly comply with ISO Standard ;
2. Sharper and More Efficient ;
3. Durable and Valuable ;
4. Full Range of shapes and grits.
5. Fastly develop new items.
6. Individual Design and OEM Service








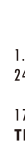














001-008M Item No.

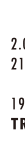Grit







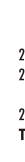
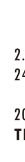

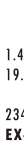
- Black / Super Coarse
- Green / Coarse
- Blue / Medium
- Red / Fine
- Yellow / Extra Fine

 0.8 19.0 001(801) BR-49	 1.0 19.0 001(801) BR-45	 1.2 19.0 001(801) BR-46	 1.4 19.0 001(801) BR-41	 1.6 19.0 001(801) BR-40	 1.8 19.0 001(801) BR-31	 2.1 19.0 001(801) BR-28	 2.3 19.0 001(801) BR-26	 2.5 19.0 001(801) BR-25	 2.9 19.0 001(801) BR-23
001-008M 001-010M 001-012M 001-014M 001-016M 001-018M 001-021M 001-023M 001-025M 001-029M									
Ball , Round /BR									
 6mm * 1.2 22.0 697(801L) BR-L46	 * 1.4 22.0 697(801L) BR-L41	 * 1.6 22.0 697(801L) BR-L40	 * 1.8 22.0 697(801L) BR-L31	 1.2/2.3 19.0 002(802) BC-43	 1.4/2.3 19.0 002(802) BC-42	 1.6/3.5 19.0 002(802) BC-32	 1.8/2.3 19.0 002(802) BC-31	 1.0(0.8)/1.0 19.0 010(805) SI-45	 1.2(0.9)/1.2 19.0 010(805) SI-46
697-012M 697-014M 697-016M 697-018M				002-012M 002-014M 002-016M 002-018M			010-010M 010-012M		
Round with Long Neck /BR				Round with Collor /BC			Inverted Conical / SI		
 1.4(1.0)/1.6 19.0 010(805) SI-47	 1.6(1.1)/1.8 19.0 010(805) SI-48	 1.0/2.3 19.0 019(806) DI-41	 1.4/3.2 19.0 019(806) DI-42	 3.2/5.0 19.0 039(811L) EX-11	 3.5/7.0 20.0 039(811L) EX-12	 4.5/0.9 19.0 041(818) WR-44	 4.2/1.8 19.0 068(909) WR-13	 1.0/4.0 19.0 109(836) SF-41	 1.2/4.0 19.0 109(836) SF-31
010-014M 010-016M 019-010M 019-014M				039-032M 039-035M		041-045M 068-042M		109-010M 109-012M	
Double Conical				Wheel		Flat End Cylinder			
 1.4/4.0 19.0 109(836) SF-42	 1.4/6.0 21.0 110(837) SF-21	 1.6/6.0 21.0 110(837) SF-15	 1.2/8.0 22.0 111(837L) SF-11	 1.4/8.0 22.0 111(837L) SF-12	 1.6/8.0 22.0 111(837L) SF-13	 1.2/8.0 22.0 130(885) CP-11	 1.4/8.0 22.0 130(885) CP-12	 1.0/4.0 19.0 139(838) SR-41	 1.2/4.0 19.0 139(838) SF-31
109-014M 110-014M		110-016M 111-012M		111-014M 111-016M		130-012M 130-014M		139-010M 139-012M	
Cylinder, end cone, pointed						Round End, Cylinder			
 1.4/4.0 19.0 139(838) SR-42	 1.4/6.0 21.0 140(880) SR-21	 1.2/8.0 22.0 141(881) SR-11	 1.4/8.0 22.0 141(881) SR-12	 1.6/8.0 22.0 141(881) SR-13	 1.2/0.2 21.0 150(839) EX-18	 1.2/6.0 21.0 164(852) TC-26	 1.4/8.0 22.0 165(858) TC-21	 1.0/10.0 23.0 166(859) TC-27	 1.2/10.0 23.0 166(859) TC-28
139-014M 140-014M		141-012M 141-014M		141-016M 150-012M		164-012M 165-014M		166-010M 166-012M	
End Cut						Needle			
 1.4/10.0 23.0 166(859) TC-10	 1.6/10.0 23.0 166(859) TC-11	 1.8/10.0 23.0 166(859) TC-16	 1.0(0.6)/3.5 19.0 169(845) TF-41	 1.2(0.7)/4.0 19.0 170(845) TF-42	 1.4(0.8)/4.0 19.0 170(845) TF-43	 1.6(0.9)/4.0 19.0 170(845) TF-31	 1.8(1.0)/5.0 20.0 170(845) TF-23	 2.1(1.2)/5.0 20.0 170(845) TF-22	 1.4(0.7)/6.5 21.0 171(846) TF-20
166-014M 166-016M		166-018M 169-010M		170-012M 170-014M		170-016M 170-018M		170-021M 171-014M	
Flat End Cone									











									
1.6(0.8)/6.5 21.0	1.0(0.6)/8.0 22.0	1.2(0.7)/8.0 22.0	1.4(0.7)/8.0 22.0	1.6(0.8)/8.0 22.0	1.8(0.9)/8.0 22.0	2.3(1.3)/9.0 22.0	1.4(0.7)/10.0 24.0	1.6(0.8)/10.0 22.5	1.8(0.9)/10.0 22.5
171(846) TF-21	172(847) TF-07	172(847) TF-08	172(847) TF-09	172(847) TF-10	172(847) TF-56	172(847) TF-14	173(848) TF-11	173(848) TF-12	173(848) TF-13
171-016M	172-010M	172-012M	172-014M	172-016M	172-018M	172-023M	173-014M	173-016M	173-018M

									
1.0(0.6)/4.0 19.0	1.2(0.7)/4.0 19.0	1.4(0.8)/4.0 19.0	1.6(0.9)/4.0 19.0	2.0(1.0)/5.0 20.0	1.4(0.7)/6.5 21.0	1.6(0.8)/6.5 21.0	1.8(1.0)/6.5 21.0	2.0(1.1)/6.5 21.0	1.4(0.8)/8.0 22.0
196(855) CR-41	196(855) CR-31	196(855) CR-42	196(855) CR-12	196(855) CR-11	197(849) TR-20	197(849) TR-21	197(849) TR-62	197(849) TR-63	198(850) TR-25
196-010M	196-012M	196-014M	196-016M	196-020M	197-014M	197-016M	197-018M	197-020M	198-014M








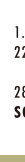


Round End Cone

									
1.6(0.9)/8.0 22.0	1.8(1.0)/8.0 22.0	1.4(0.7)/10.0 22.0	1.6(0.8)/10.0 22.0	1.8(0.9)/9.0 22.0	2.2(1.3)/9.0 22.0	2.2(1.3)/11.0 24.0	2.4(1.4)/11.0 24.0	1.0(0.8)/2.8 19.0	1.4(1.1)/4.0 19.0
198(850) TR-24	198(850) TR-23	199(852) TR-11	199(852) TR-12	199(852) TR-13	199(852) TR-14	200(852L) TR-15	200(852L) TR-19	233(808) EX-41	234(808L) EX-31
198-016M	198-018M	199-014M	199-016M	199-018M	199-022M	200-022M	200-024M	233-010M	234-014M

Pear









									
1.8(1.2)/5.0 19.0	2.1(1.6)/5.0 19.0	3.2(1.9)/5.0 19.0	1.2/6.0 21.0	1.4/6.0 21.0	1.4/8.0 22.0	1.6/8.0 22.0	1.2/10.0 24.0	1.4/10.0 24.0	1.6/10.0 24.0
237(808R) EX-20	237(808R) EX-21	237(808R) EX-26	248(861) FO-41	248(861) FO-20	249(862) FO-21	249(862) FO-22	250(863) FO-11	250(863) FO-12	250(863) FO-13
237-018M	237-021M	237-032M	248-012M	248-014M	249-014M	249-016M	250-012M	250-014M	250-016M

Flame

									
1.8/3.6 19.0	2.3/5.0 19.0	2.8/5.0 19.0	3.2/5.0 19.0	2.0/4.0 19.0	3.0/5.0 19.0	1.2/6.0 21.0	1.4/8.0 22.0	1.6(1.0)/8.0 22.0	1.8(1.1)/8.0 22.0
257(368,830) FO-32	257(368,830) FO-34	257(368,830) FO-25	257(368,830) FO-27	277(379,833) TR-22	277(379,833) TR-28	288(867,877) SO-20	289(868,878) SO-21	298(878K) 878K16	298(878K) 878K18
257-018M	257-023M	257-028M	257-032M	277-020M	277-030M	288-012M	289-014M	298-016M	298-018M





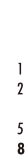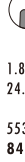



Egg








Torpedo

							
1.6(1.0)/10.0 24.0	1.4/8.0 22.0	1.6/8.0 22.0	cutting depth 0.3mm 1.6(1.0)/6.0 21.0	cutting depth 0.4mm 1.8(1.0)/6.0 21.0	cutting depth 0.5mm 2.1(1.1)/6.0 21.0	1.6(0.8)/8.0 22.0	1.8(0.9)/10.0 24.0
299(879K) 879K16	115(842SP) BC-11	178(848SP) BC-12	552(834) VP-20	552(834) VP-21	552(834) VP-22	546(847KR) 847KR16	553(848KR) 848KR18
299-016M	115-014M	178-016M	552-016M	552-018M	552-021M	546-016M	553-018M

Depth Marking


Taper, Round Edge

							
Super Coarse Black	1.4 19.0	1.6 19.0	1.8 19.0	2.1 19.0	2.3 19.0	2.5 19.0	1.6(1.1)/1.8 19.0
001(801) BR-41S	001(801) BR-40S	001(801) BR-31S	001(801) BR-28S	001(801) BR-26S	001(801) BR-25S	001(801) BR-23S	010(805) SI-48S
001-014S	001-016S	001-018S	001-021S	001-023S	001-025S	001-029S	010-016S

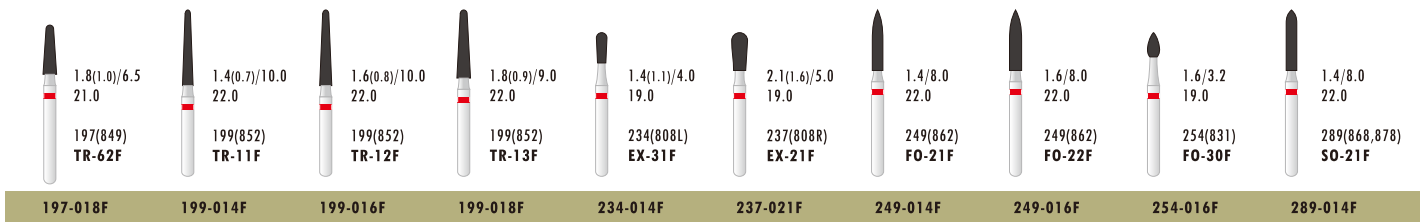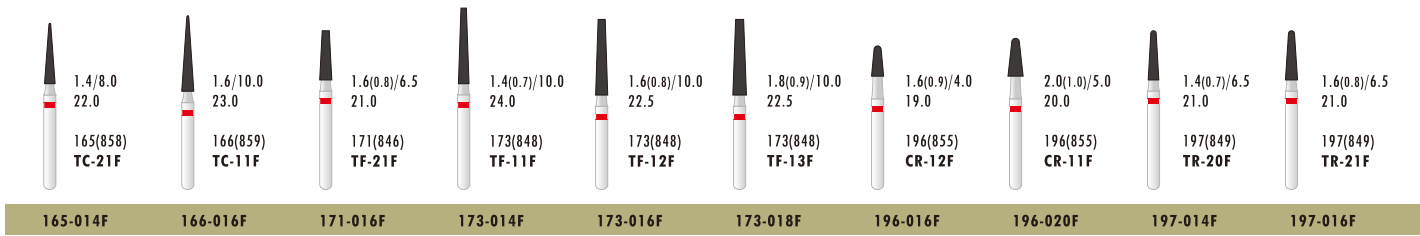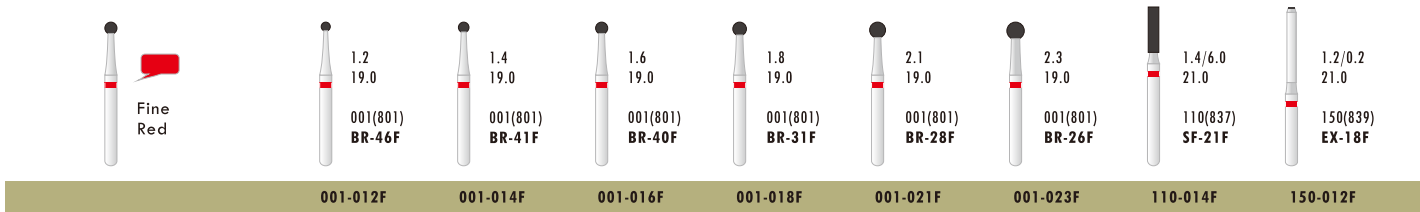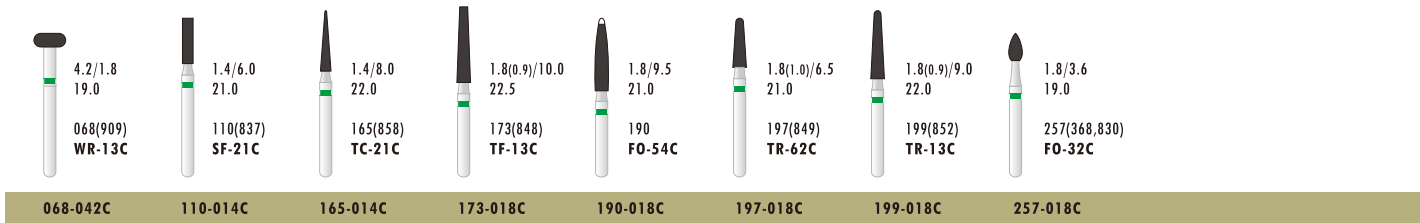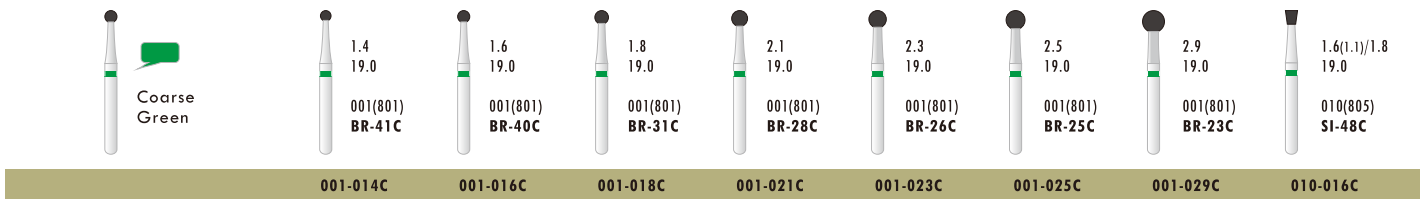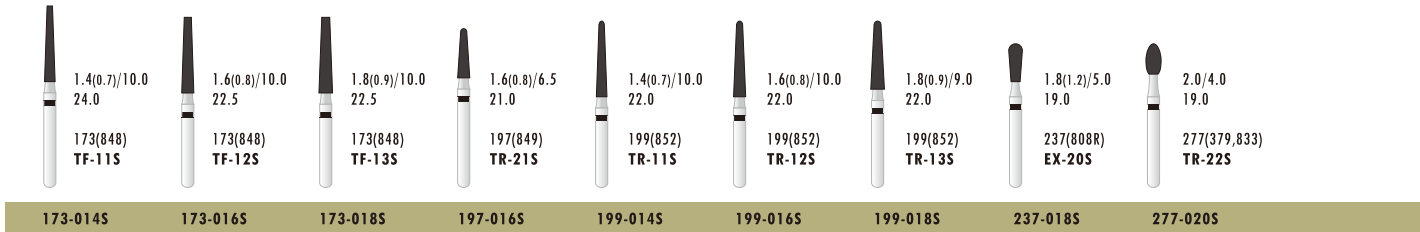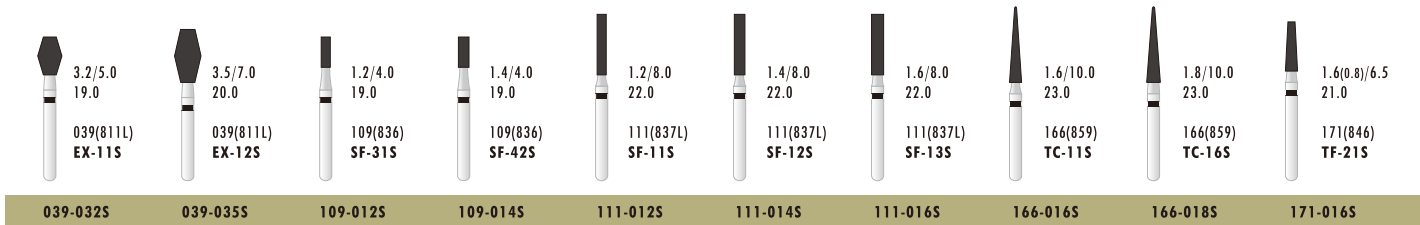




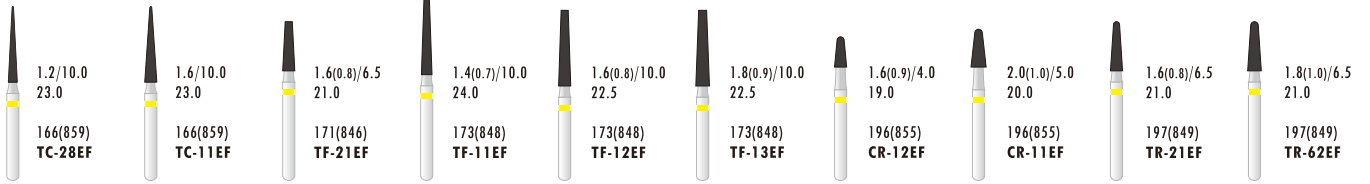
001-012EF 001-014EF 001-016EF 001-018EF 001-021EF 001-023EF 165-014EF 166-010EF




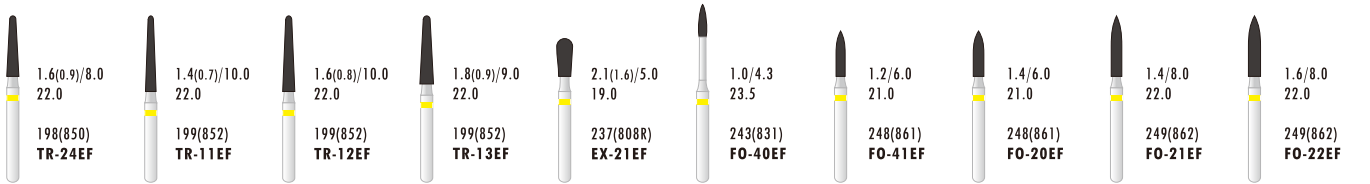
166-012EF 166-016EF 171-016EF 173-014EF 173-016EF 173-018EF 196-016EF 196-020EF 197-016EF 197-018EF



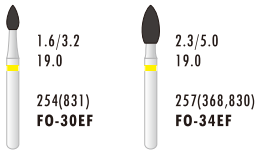
198-016EF 199-014EF 199-016EF 199-018EF 237-021EF 243-010EF 248-012EF 248-014EF 249-014EF 249-016EF





254-016EF 257-023EF

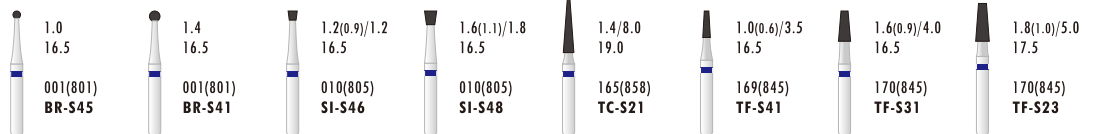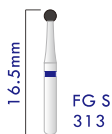





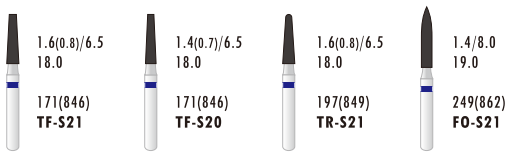
FG Short 001-S010M 001-S014M 010-S012M 010-S016M 165-S014M 169-S010M 170-S016M 170-S018M



171-S016M 171-S014M 197-S016M 249-S014M

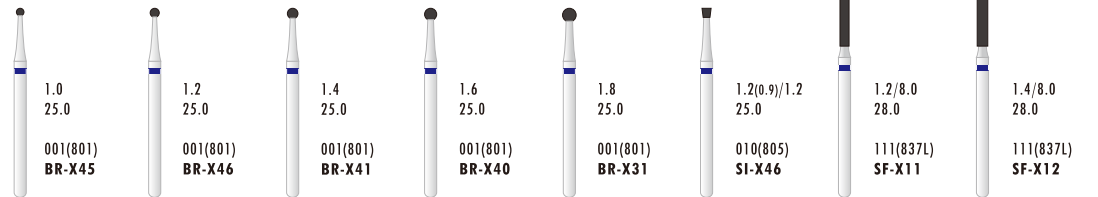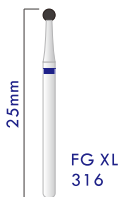

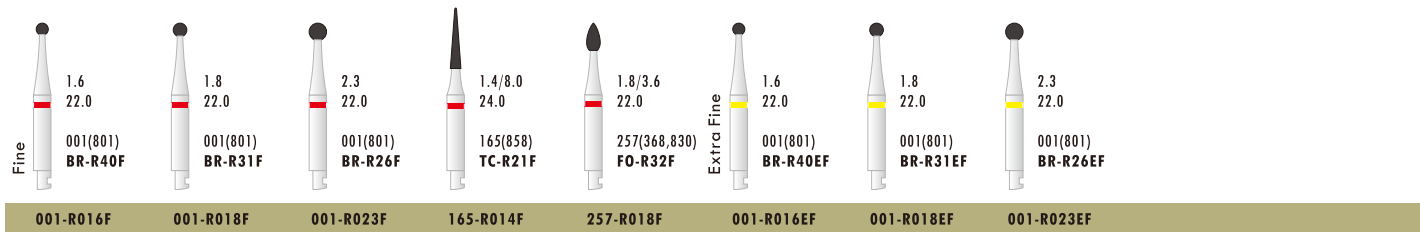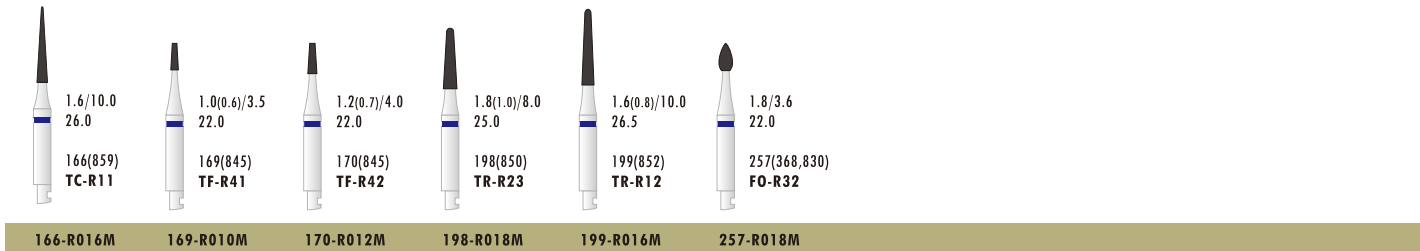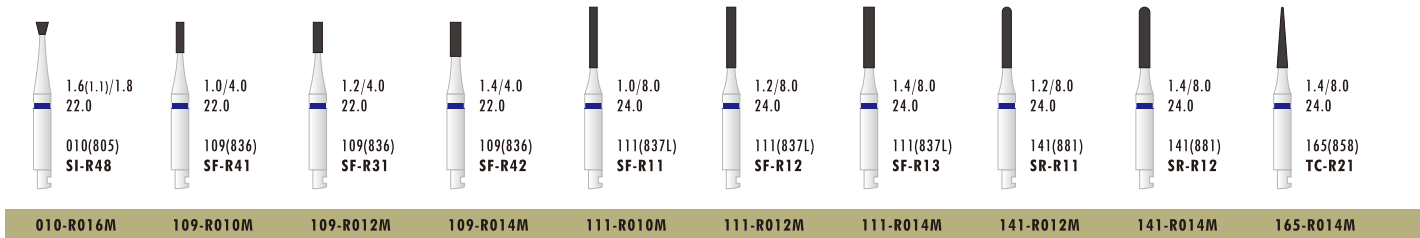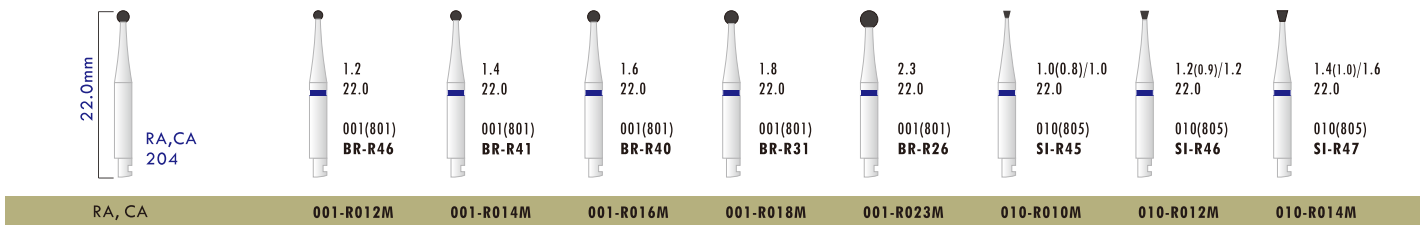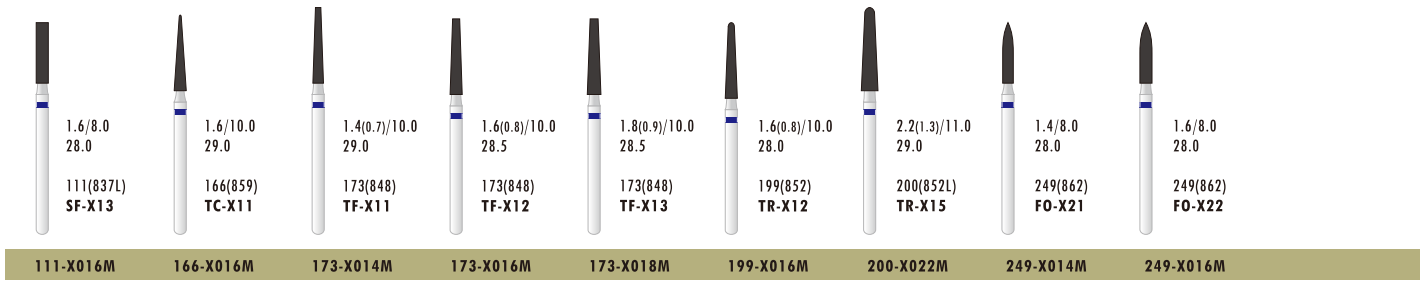FG Surgical 001-X010M 001-X012M 001-X014M 001-X016M 001-X018M 010-X012M 111-X012M 111-X014M





Tungsten Carbide Bur

钨钢车针




Advantages

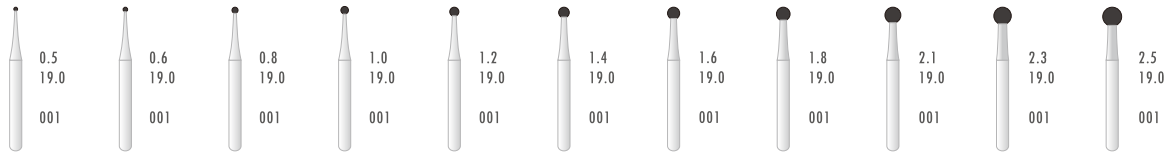
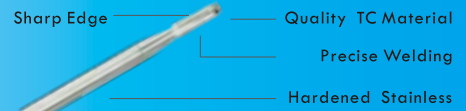
1. Strictly comply with ISO Standard
2. Sharper and More Efficient
3. Durable and Valuable
4. FG, FG Short, FG Long, RA, RA Long etc




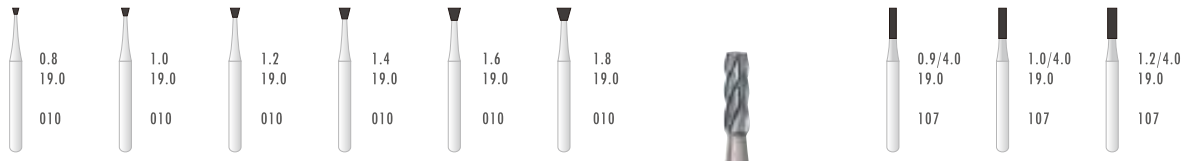
Ball, Plain Cut FG 1/4 FG 1/2 FG 1 FG 2 FG 3 FG 4 FG 5 FG 6 FG 7 FG 8 FG 9



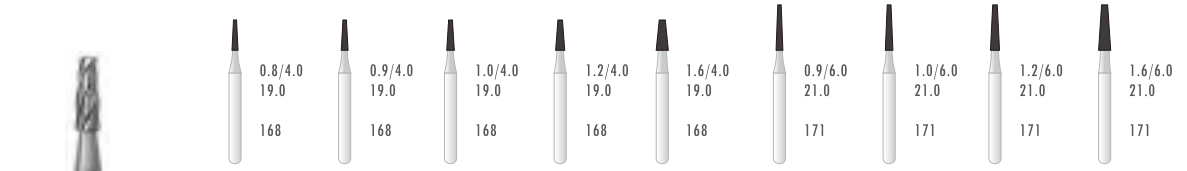
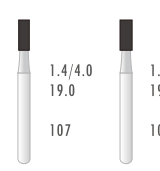
Inverted Cone FG 34 FG 35 FG 36 FG 37 FG 38 FG 39 Cylinder, Plain Cut FG 56 FG 57 FG 58



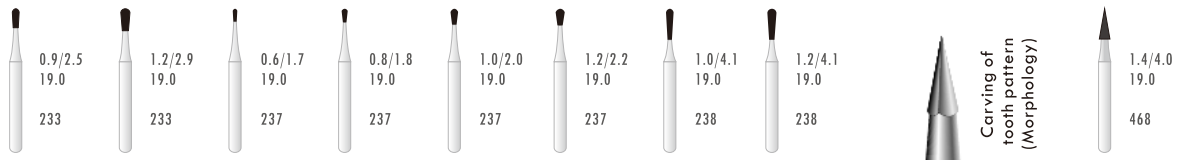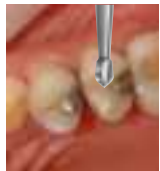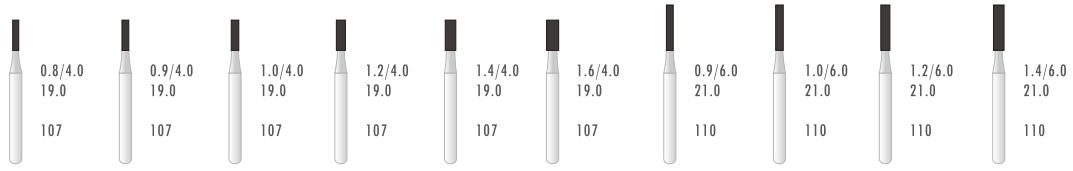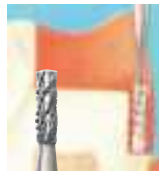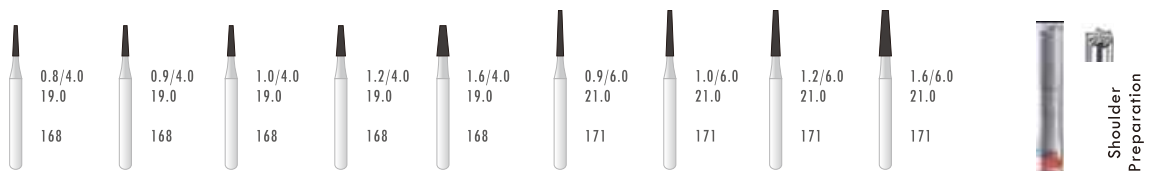
FG 59 FG 60 Flat End Taper, Plain Cut FG 168 FG 169 FG 170 FG 171 FG 172 FG 169L FG 170L FG 171L FG 172L














Pear Plain Cut FG 245 FG 246 FG 329 FG 330 FG 331 FG 332 FG 331L FG 332L X-Mas Tree, 3 slided FG 468























Cylinder, Cross Cut FG 555 FG 556 FG 557 FG 558 FG 559 FG 560 FG 556L FG 557L FG 558L FG 559L















Flat End Taper, Cross Cut FG 698 FG 699 FG 700 FG 701 FG 702 FG 699L FG 700L FG 701L FG 702L End Cutting


 1.0/4.0 19.0 150	 1.2/4.0 19.0 150		 0.8/4.0 19.0 137	 0.9/4.0 19.0 137	 1.0/4.0 19.0 137	 1.2/4.0 19.0 137	 1.4/4.0 19.0 137	 1.6/4.0 19.0 137		 0.9/4.0 19.0 194
FG 957	FG 958	Round Cylinder, Plain Cut	FG 1155	FG 1156	FG 1157	FG 1158	FG 1159	FG 1160	Round Taper, Plain Cut	FG 1169

 1.0/4.0 19.0 194	 1.2/4.0 19.0 194	 1.6/4.0 19.0 194		 0.9/4.0 19.0 137	 1.0/4.0 19.0 137	 1.2/4.0 19.0 137	 1.4/4.0 19.0 137		 1.2/4.0 19.0 194	 1.6/4.0 19.0 194
FG 1170	FG 1171	FG 1172	Round Cylinder, Cross Cut	FG 1556	FG 1557	FG 1558	FG 1559	Round Taper, Cross Cut	FG 1701	FG 1702


	 1.0/2.0 19.0 237	 1.2/2.2 19.0 237		 1.0/4.0 19.0 137	 1.2/4.0 19.0 137		 1.0/4.0 19.0 107	 1.2/4.0 19.0 107
Pear, Cross Cut	FG 1931	FG 1932	Round Cylinder, Cross Cut	FG 1957	FG 1958	Round Taper, Cross Cut	FG 2057	FG 2058

	 0.9/3.0 19.0 6-side Perio	 1.2/5.5 21.0 6-side Flame
Scaling to remove Tartar	ET3	ET5


 1.6.5mm FG S 313	 0.8 001	 1.0 001	 1.2 001	 1.4 001	 1.6 001	 0.9/2.5 233	 0.8/1.8 237	 1.0/2.0 237	 0.9/4.0 107	 1.0/4.0 107	 1.0/4.0 137
FG Short	FG S1	FG S2	FG S3	FG S4	FG S5	FG S245	FG S330	FG S331	FG S556	FG S557	FG S1157

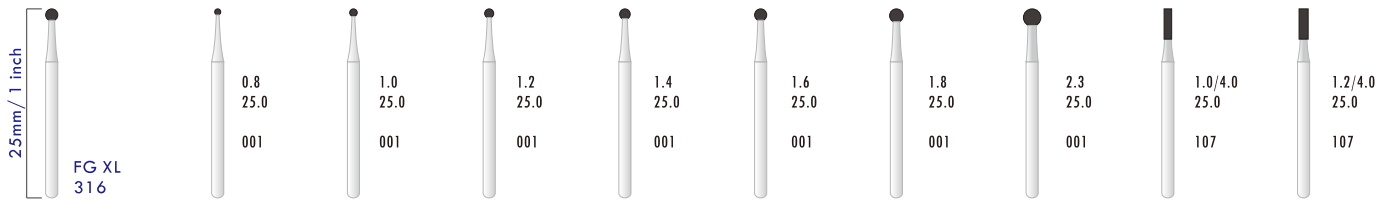

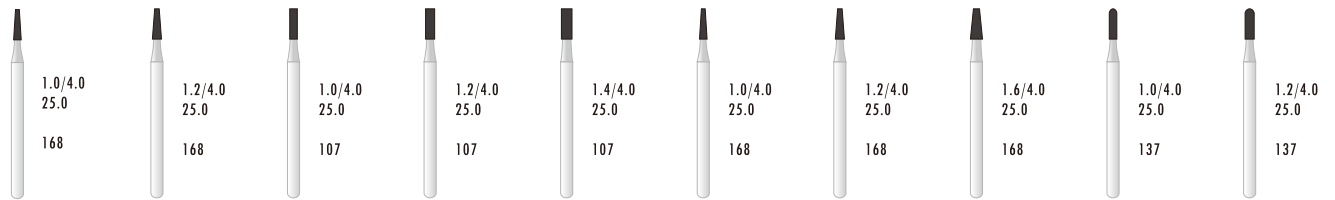
FG Surgical FG X1 FG X2 FG X3 FG X4 FG X5 FG X6 FG X8 FG X57 FG X58



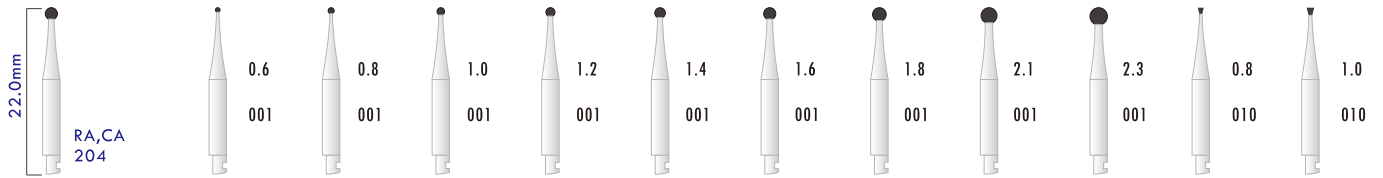
FG X170 FG X171 FG X557 FG X558 FG X559 FG X700 FG X701 FG X702 FG X1557 FG X1558



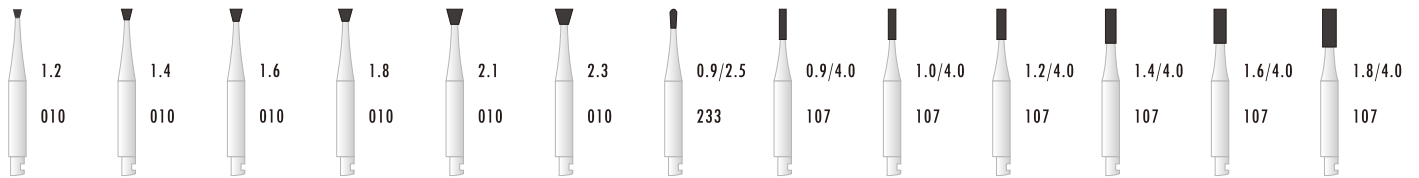
RA, CA RA 1/2 RA 1 RA 2 RA 3 RA 4 RA 5 RA 6 RA 7 RA 8 RA 34 RA 35



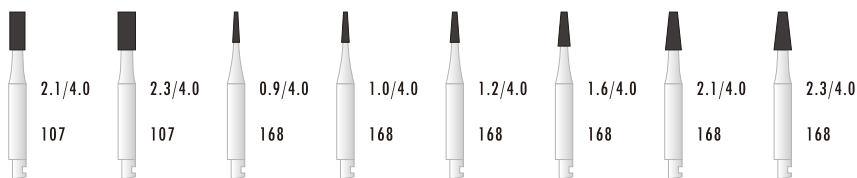RA 36 RA 37 RA 38 RA 39 RA 40 RA 41 RA 245 RA 556 RA 557 RA 558 RA 559 RA 560 RA 561



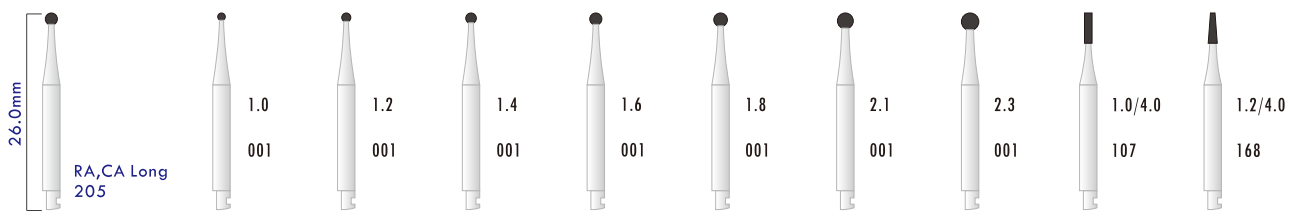
RA 562 RA 563 RA 699 RA 700 RA 701 RA 702 RA 703 RA 704





RA, CA RA X2 RA X3 RA X4 RA X5 RA X6 RA X7 RA X8 RA X557 RA X701

Zekrya













Round End Taper	151	X151	XL151	Bone Cutter	161	Bone Cutter	162	H253E	H254E
	1.6(0.9)/11 23.0 199	1.6(0.9)/11 25.0 199	1.6(0.9)/11 28.0 199		1.6(1.0)/10.0 25.0 408		1.6(1.0)/10.0 25.0 408	1.2(0.7)/6.0 28.0 415	1.2(0.7)/6.0 34.0 415













Flat End Taper , Cylinder	EC25	EC28	EC32	ES25	ES28	Round End Taper	152	X152	XL152
	1.6(1.2)/8.0 25.0 172	1.6(1.2)/8.0 28.0 172	1.6(1.2)/8.0 32.0 172	1.6/8.0 25.0 111	1.6/8.0 28.0 111		1.6(0.9)/9.0 21.0 210	1.6(0.9)/9.0 23.0 210	1.6(0.9)/9.0 25.0 210







Finishing Carbide Burs

Round	FG 7004	FG 7005	FG 7006	Flat End Taper	FG 7204	FG 7205	Pointed Taper	FG 7611	FG 7612
	1.4 19.0 001	1.6 19.0 001	1.8 19.0 001		1.4/9.0 23.0 184	1.6/9.0 23.0 184		1.0/4.2 19.0 159	1.4/6.0 21.0 164

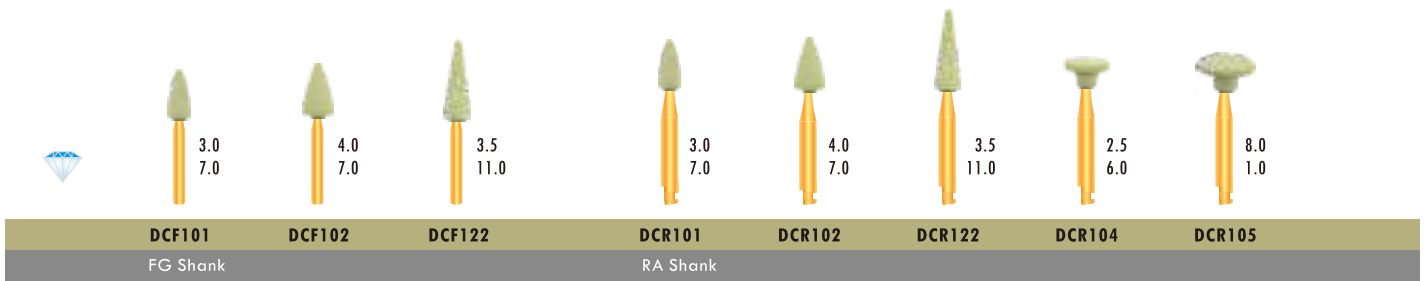










Egg	FG 7404	FG 7406	Flat End Taper	FG 7114	FG 7214	FG 7714	Round End Taper	FG 7675
	1.4/3.3 19.0 277	1.8/3.5 19.0 277		0.9/3.5 19.0 168	1.0/5.5 21.0 170	1.4/9.0 23.0 172		1.6(1.0)/8.0 22.0 198





8 - 12 Blades	Flame	FG 7901	FG 7902	FG 7903	FG 7903L
		0.9/3.6 19.0 496	1.0/3.6 19.0 496	1.2/3.6 19.0 496	1.2/9.0 23.0 250

Diamond Stone 陶瓷金刚车石 Finish Zirconia, Lithium Disilicate, Ceramics . 研磨氧化锆、玻璃陶瓷、烤瓷。



Diamond Rubber Polishers 橡胶金刚石抛光磨头

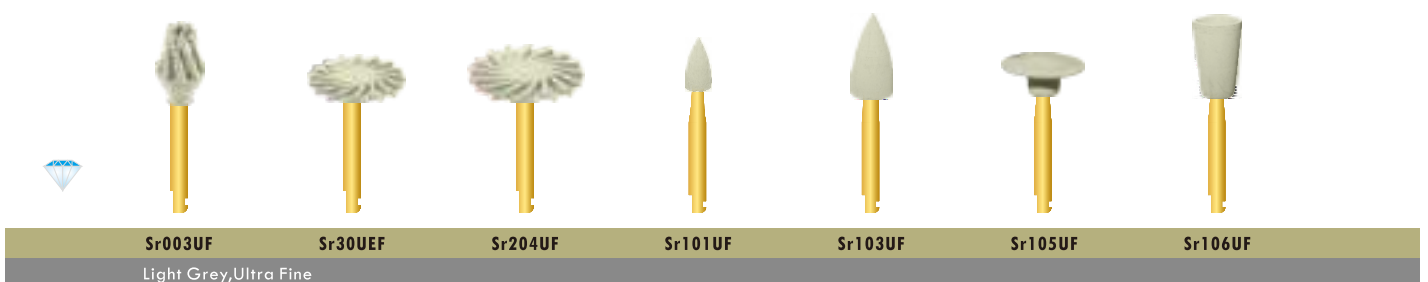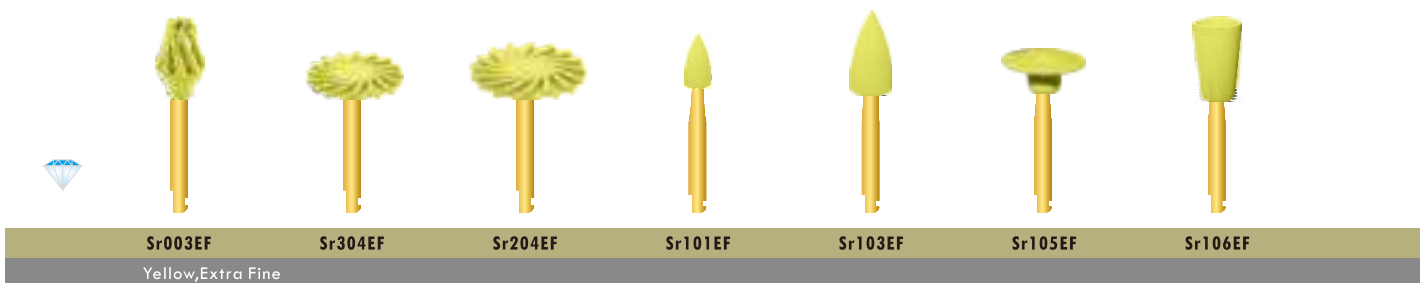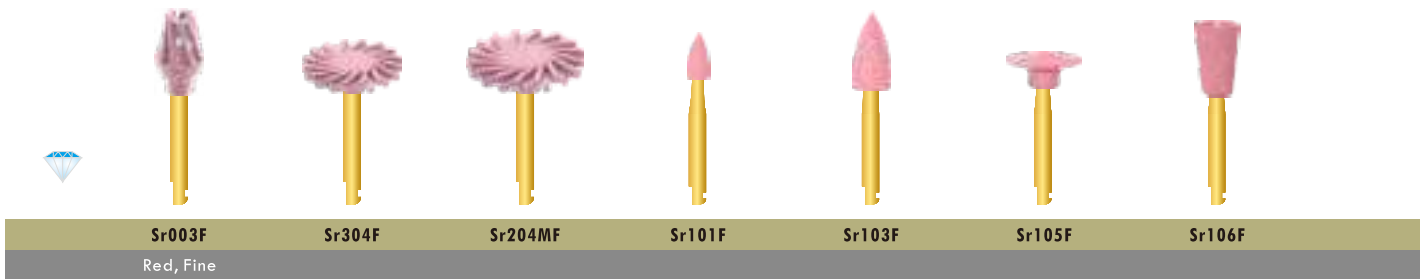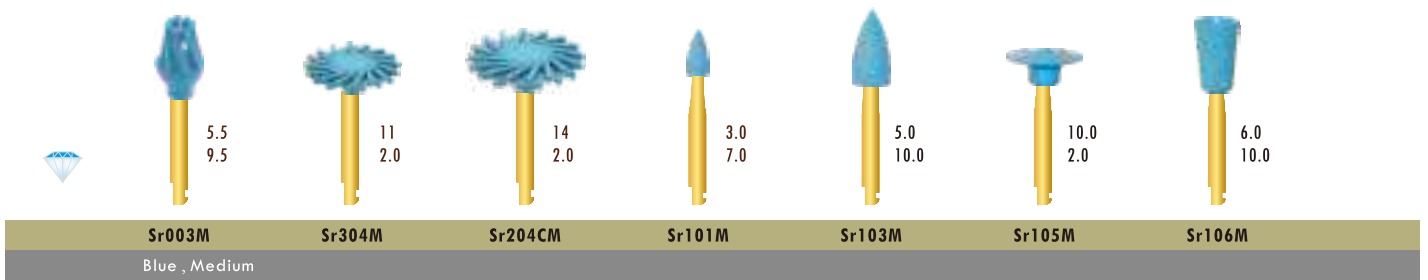Universally used to polish Zirconia, Lithium Disilicate, Ceramics, Composite etc.
广泛的用于抛光氧化锆、玻璃陶瓷、烤瓷、树脂等。



Cleaning Burs after debonding (脱胶后)洁牙磨头

Remove Calculus or Stains , No damage to the enamel.

3 pieces mixed



2.3
9.0



2.3
7.0



2.3
7.5



2.3
9.0



2.3
7.0



2.3
7.5

PC00C

PC02C

PC03C

PC00

PC02

PC03

Coarse

Medium

4



3.0
7.5



5.0
13.0



10.0
0.7



6.0
9.0

BR101

BR103

BR104

BR106

Silicon

矽粒子

25



SM103

6x10 mm



SM104

10x0.7 mm



SM106

6x9 mm

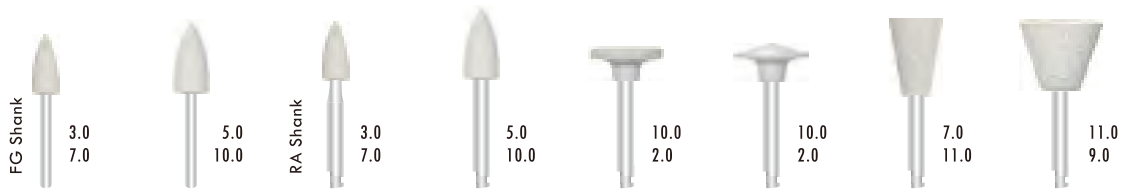


35RS

Silicon Polisher

For Ceramics

硅胶磨头
烤瓷抛光




White, Coarse

Cf101C Cf103C Cr101C Cr103C Cr104C Cr105C Cr106C Cr108C

Silicon Polisher

For Ceramics

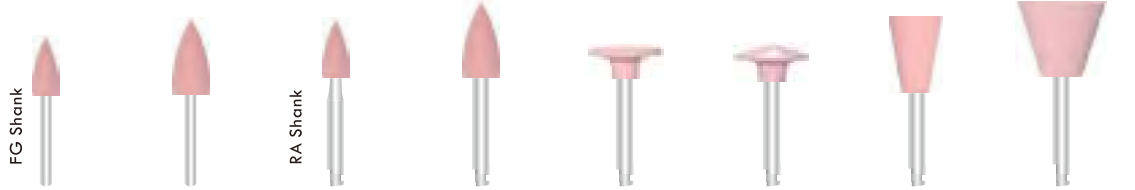硅胶磨头
烤瓷抛光




Pink, Medium

Cf101M Cf103M Cr101M Cr103M Cr104M Cr105M Cr106M Cr108M

Silicon Polisher

For Ceramics, Composite

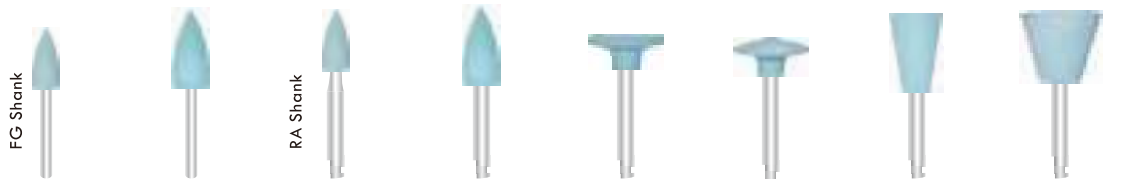硅胶磨头
烤瓷、树脂抛光



Blue, Fine

Cf101F Cf103F Cr101F Cr103F Cr104F Cr105F Cr106F Cr108F

Silicon Polisher

For Metal, Alloy

硅胶磨头
金属、合金抛光



Black, Coarse

Af101C Af103C Ar101C Ar103C Ar104C Ar105C Ar106C Ar108C

Silicon Polisher

For Metal, Alloy

硅胶磨头
金属、合金抛光



Brown, Medium

Af101M Af103M Ar101M Ar103M Ar104M Ar105M Ar106M Ar108M

Silicon Polisher

For Metal, Alloy, Composite

硅胶磨头
金属、合金、树脂抛光



Fine, Light Grey

Af101F Af103F Ar101F Ar103F Ar104F Ar105F Ar106F Ar108F

Silicon Polisher

For Composite

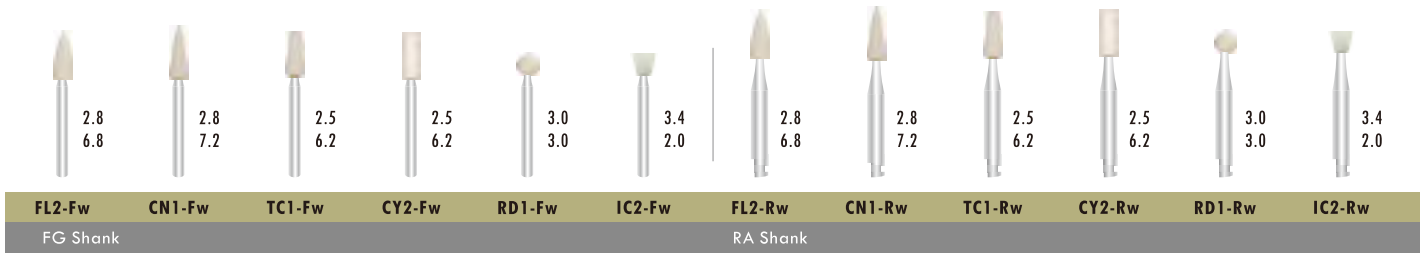
硅胶磨头
树脂抛光




Fine, Yellow

Rf101F Rf103F Rr101F Rr103F Rr104F Rr105F Rr106F Rr108F

White Stone 白车石 Finish, Polish Filling Composite 去树脂, 填充物及抛光。



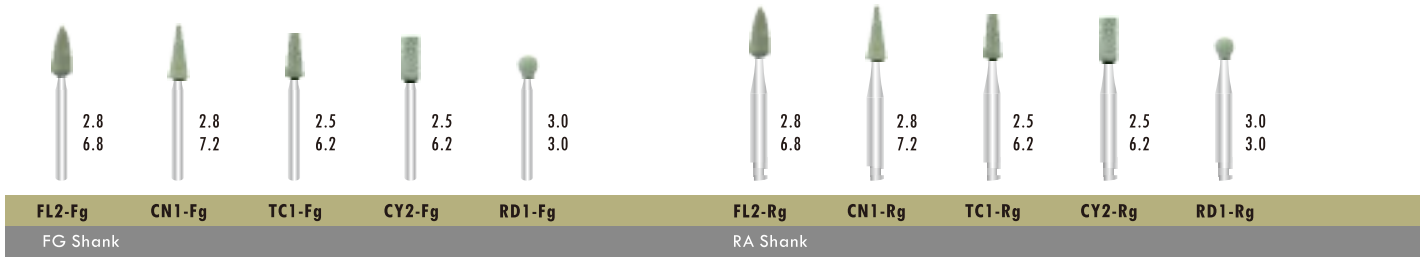Green Stone 绿车石 Grind Ceramics, Precious Metal 研磨烤瓷, 贵金属等。



Polishing Brushes

Nature Material, No hurt to teeth
Final polishing composite,
nature teeth or Ceramics
天然毛质、不伤牙
适用于树脂、真牙、烤瓷等



R1Ns R1Ry R1G R1T

Polishing Brushes

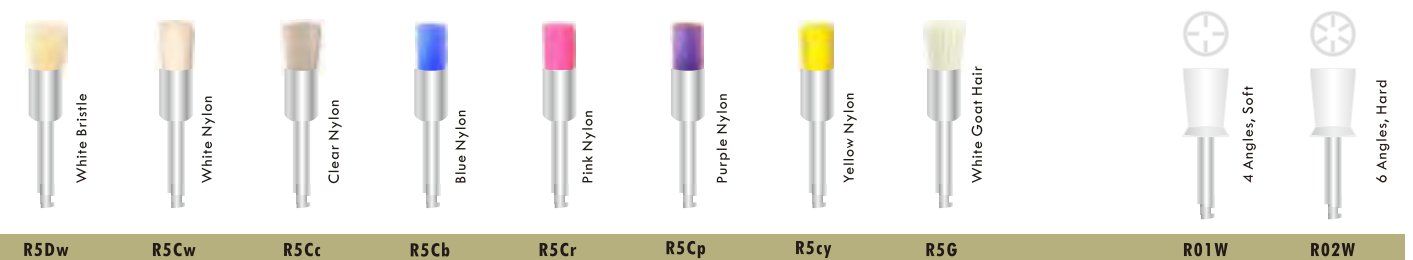
Abrasive Fiber
磨料丝 (含碳化硅
或氧化铝)
For Composite,
Ceramics



R1AI R4AI R5AI R7AI R1SiC R4SiC R5SiC R7SiC



R4Dw R4Cw R4Cc R4Cb R4Cr R4Cp R4Cy R4G R7Dw R7Cc R7G



R5Dw R5Cw R5Cc R5Cb R5Cr R5Cp R5cy R5G R01W R02W

Mandrels

夹针



31RS



31RSL
RA X Long shank



35RS



52RS



53RS



26RS
Adapter FC into RA



Disposable use
made of steel

51R



No. 201
for FG burs * 15



No. 211
for 10*FG + 5*RA



No. 221
for 5*FG + 10*RA



No. 241
for RA burs * 12



No. 202
for Root canal files * 18



No. 212
For HP burs * 18



No. 015A
for 18*FG



No. 015C
for 12*FG + 6*RA



No. 015E
for 18*RA



No. 015D
for 20*FG + 10*RA



No. 013E
for 16*FG + 9*RA



No. 022D
for 18*RA



No. 016B
for 8*HP



No. 023
132 * 107 * 65 mm
for FG, RA, HP, Files Mixed

Plastic Blocks, Sterilisable

塑料消毒盒

KIT 套装

Provide kinds of kits for clinic use , which offer the smartest operation ;
 Whatever for Ceramics, Zirconia , Lithium Disilicate, Composite, Metal or Alloy ;
 From Preparation to Final polishing ;
 High Efficiency to shining, No need Glaze after polishing .



BR1901

BR-49, BR-45, BR-46, BR-41, BR-40,
 BR-31, BR-28, BR-26, BR-25, BR-23 .

Diamond burs of Ball shape Assorted in FG shank
 圆球状金刚石车针混装(高速)



BR1902

BR-X45, BR-X46, BR-X41, BR-X40, BR-X31,
 BR-X40S, BR-X31S, BR-X41C, BR-X40C, BR-X31C


Diamond burs of Ball shape Assorted in FG LONG
 圆球状金刚石车针混装(高速加长)



BR1903

BR-R46, BR-R41, BR-R40, BR-R31, BR-R26,
 BR-R40F, BR-R31F, BR-R26F, BR-R40EF, BR-R26EF

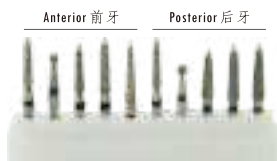Diamond burs of Ball shape Assorted in RA shank
 圆球状金刚石车针混装(低速)



4100

508T-014S, BR-41, 198-016S, 508T-016S, 508T-016F ;
 508T-014S, BR-S40, 198-021S, 508T-021S, 508T-021F

Safe End Diamond Preparation
 前后牙预备(安全端金刚石车针)



4101

TF-21, TF-10, TF-12, TC-10, WR-13
 SR-11, TR-11, EX-18, TR-12F, TC-11F

Crown & Bridge Preparation for Anterior
 前牙桥体/冠体预备



4102

EX-20S, SR-21, FO-34, TR-62C, TR-21
 TR-24, SO-21, EX-18, TR-21F, EX-21F

Crown & Bridge Preparation for Posterior
 后牙桥体/冠体预备



18041V

SI-45 x 2 ; SI-46 x 3,
 SI-47 x 3, SI-48 x 2 .

Diamond burs of Inverted Cone Assorted in FG shank
 常用倒锥金刚石车针混装(高速)



1805TA

TR-11 x 2, TR-12 x 2, TF-11 x 2 ;
 TF-12, TC-10, TC-27, TC-28 x 1

Regular Diamond burs of Taper Assorted in FG shank
 常用锥形金刚石车针混装 (高速)



1806

TC-11EF, 820-016EF, 820-031F,
 820-016F, TC-10F

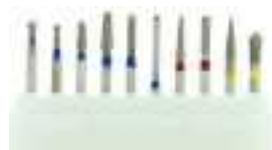
Interproximal Enamel Reduction (IPR) for Ortho
 邻面去釉套装



4103

BR-45, BR-40, BC-32, EX-31
 TF-20, SF-41, SF-21, SI-48


Cavity Preparation, Restoratives
 备洞



4104

SI-48, DI-42, TF-31, TF-22, SF-21
 EX-18, TF-21F, SF-21F, FO-22EF, EX-21EF

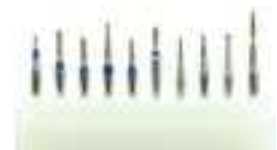Inlay & Onlay Preparation
 嵌体/高嵌体预备



4105

BR-49, SF-41, SR-41, TC-26, CR-41
 EX-41, FO-40EF / FG1/4, FG1/2, FG 1

Micro Invasive Surgery
 微创备牙



6000

BR-45, VP-20, VP-21, VP-22,
SR-41, TR-25, SO-20, TR-13F, CF101F

Porcelain Veneers' Preparation
瓷贴面备牙套装



6002

CD-50F, CD-51F, CD-52F, CD-53F, CD-54F
CD-55F, CD-56F, CD-57F, CD-58F, CD-59F

Children's Teeth Preparation
儿童牙预备



4106

BR-40EF, TC-11EF, TF-11EF, TF-12EF, TR-11EF,
EX-21EF, FO-40EF, FO-21EF, FO-22EF, FO-30EF

Ceramics / Composite Polishing
烤瓷牙/复合树脂抛光



FG008D

BR-49, BR-49, SI-49, SF-49, TF-49

Student or Beginner's Diamond FG 0.8mm
新手入门金刚石车针,高速



FG009D

001-009M, 001-009M, 010-009M,
109-009M, 170-009M

Student or Beginner's Diamond FG 0.9mm
新手入门金刚石车针,高速



RA008D

BR-R49, BR-R49, SI-R49, SF-R49, TF-R49

Student or Beginner's Diamond RA 0.8mm
新手入门金刚石车针,慢速



RA009D

001-R009M, 001-R009M, 010-R009M,
109-R009M, 170-R009M

Student or Beginner's Diamond RA 0.9mm
新手入门金刚石车针,慢速



TF7100

7004, 7404, 7214, 7714, 7903L

TC Trimming & Finishing burs, F6
钨钢精修车针, FG



1900 T

FG 1, FG 2, FG 3, FG 4, FG 5

Tungsten carbide burs of Ball shape Assorted in FG shank
圆球状钨钢车针混装(高速)



FG008T

FG 1, FG 1, FG 34, FG555, FG698

Tungsten Carbide of FG shank 0.8mm for Students
初学者用的钨钢车针 (头径0.8mm, 高速)



RA008T

RA 1, RA 1, RA 34, RA 555, RA 698

Tungsten Carbide of RA shank 0.8mm for Students
初学者用的钨钢车针 (头径0.8mm, 低速)



7109FW

FL2-Fw, CN1-Fw, RD1-Fw, RD2-Fw, TC1-Fw,

FG White Stone Burs, to Finish Composite
白车石, 高速



7109RW

FL2-Rw, CN1-Rw, RD1-Rw, RD2-Rw, TC1-Rw,

RA White Stone Burs, to Finish Composite
白车石, 低速



7109FG

FL2-Fg, CN1-Fg, RD1-Fg, RD2-Fg, TC1-Fg,

FG Green Stone burs, Finish Ceramics
绿车石, 高速



7109RG

FL2-Rg, CN1-Rg, RD1-Rg, RD2-Rg, TC1-Rg,

RA Green Stone burs, Finish Ceramics
绿车石, 低速



1906

Cr103C,Cr 104C, Cr106C;Cr103M, Cr104M, Cr106M, R4AL, R1AL,R5AL; R1Ry, R1G, R1T

Composite Finish and Polish RA
树脂成型和抛光,弯机柄



1907

Cr103C x 2,Cr 104C, Cr106C; , Cr103M, Cr104M; R4AL, R5AL; R4Cb, R5Cb, R4Cr, R5Cr

Composite Finish and Polish RA
树脂成型和抛光,弯机柄



QR2101

2 step : Sr003F : Pro-Polishing
Sr003UF: Polishing

Composite Polish, RA
树脂抛光,弯机柄



QR2102

2 step : Sr304F : Pro-Polishing
Sr304UF, Sr003UF: Polishing

Composite Polish, RA
树脂抛光,弯机柄



QR2103

3 step : Sr304CM: Smooth;
Sr003F : Pro-Polishing; Sr003UF: Polishing

Composite Polish, RA
树脂抛光,弯机柄



QR2104

3 step : Sr003M:Smooth;
Sr003F : Pro-Polishing; Sr003UF x 2: Polishing

Composite Polish, RA
树脂抛光,弯机柄



CP 6011

2 Step : Sr 204 M x 3: Pro-polish;
Sr 204 F x 3: Polish

Composite, Ceramic, Zirconia Polishers . RA
树脂、烤瓷、氧化锆抛光



CP 6001

3 Step: Sr 204 C x 2 : Smooth;
Sr 204 M x 2 : Pro-polish; Fine : Sr 204 F x 2: Polish

Composite, Ceramic, Zirconia Polishers . RA
树脂、烤瓷、氧化锆抛光



7115

Smooth : DCR101 , DCR102, DCR122, DCR105
Pro-polish: Sr101F, Sr103F, Sr105F, Sr106F
Polish: Sr101UF, Sr103UF, Sr105UF, Sr106UF

Diamond stones & polishers for Zirconia / Lithium Disilicate, 氧化锆,玻璃陶瓷研磨抛光



1904

X151, XL151, EC25, EC28, FG X6,
EX-28, BR-X31, TR-X12, TR-X15, TF-X13

Extraction Burs
拔阻生齿车针



1905

FG X2, FG X4, FG X6, FG X558, FG X701,
BR-X41, BR-X40, BR-X31, TF-X12, TR-X12

Oral Surgery burs FG Long
手术车针(高速加长)



4110V

Smooth: Rw01, Rw02 x 2
Pro-polish: Cr101M, Cr103M, Cr104M, Cr106M,
Polish: Cr101F, Cr103F, Cr104F, Cr106F

Composite Finish and Polish RA
树脂成型和抛光,弯机柄



4112

Smooth: Cr101C, Cr103C, Cr104C, Cr106C
Pro-polish : Cr101M, Cr103M, Cr104M, Cr106M
Polish: Cr101F, Cr103F, Cr104F, Cr106F

Porcelain / Natural Teeth Polishing RA
烤瓷牙,自然牙抛光,弯机柄



4114

Smooth: Ar101C, Ar103C, Ar104C, Ar106C
Pro-polish: Ar101M, Ar103M, Ar104M, Ar106M
Polish: Ar101F, Ar103F, Ar104F, Ar106F

Amalgam / Metal Polishing RA
银汞合金,金属抛光,弯机柄



4115

Smooth : Sr103M, Sr105M, Sr106M ; Pro-polish : Sr103F, Sr105F, Sr106F
Polish: Sr101EF, Sr103EF, Sr105EF, Sr106EF,
Final Polish: Sr101UF, Sr103UF

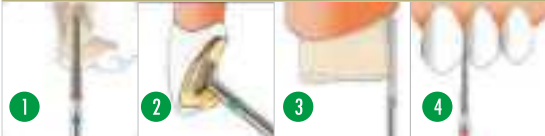
Diamond Polishers for Zirconia /Lithium Disilicate
氧化锆,玻璃陶瓷抛光,弯机柄

Crown & Bridge Preparation for Anterior
前牙桥体/冠体预备



4101A

TF-21, TF-10, TF-12 x 2, WR-13, TC-21,
TR-11, TC-10/EX-18, SR-11/ TR-12F, TC-11F



Crown & Bridge Preparation for Posterior
后牙桥体/冠体预备



4102A

EX-11, EX-20S, TR-22S, BR-31C, TR-62C, FO-34/
TR-24, TR-21/EX-18, SO-21/EX-21F, TR-21F



Cavity Preparation, Restoratives
备洞



4103A

BR-45, BR-40, BR-46, BR-31, BR-41/BC-32,
SI-47, SI-48/SF-41, SF-21/EX-31, TF-20

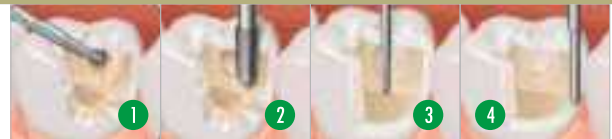


Inlay/Onlay Preparation
嵌体/高嵌体预备



4104A

D1-42, SI-48, TF-22, TF-31 / TF-21F, SF-21F /
TF-21, SF-21, EX-18, BC-32 / EX-21EF, FO-22EF




Porcelain Veneer's Preparation
瓷贴面备牙




6000A

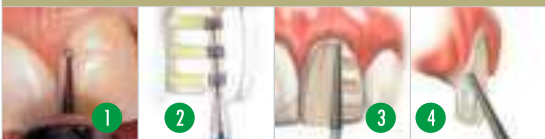SR-41, BR-45 / VP-20 x 2, VP-21 x 2, VP-22 x 2/
SO-20, TR-25, TR-13F, TC-28EF



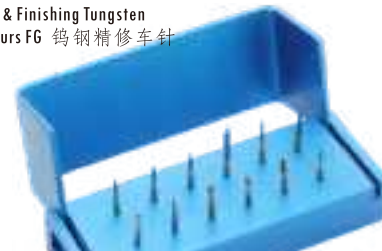Trimming & Finishing Tungsten
Carbide burs FG 钨钢精修车针



6003A

7004, 7005, 7006, 7404, 7406, 7114,
7214, 7714, 7901, 7902, 7903, 7903L

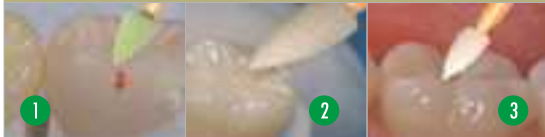


Zirconia, Lithium Disilicate Polishing kit, RA
氧化锆, 玻璃陶瓷抛光(低速)



7115A

Smooth : DCR101 , DCR102, DCR122, DCR105 ; Pro-polish: Sr101F, Sr103F, Sr105F, Sr106F
Polish: Sr101UF, Sr103UF, Sr105UF, Sr106UF




Zirconia, Lithium Disilicate High Gloss polishing kit, RA
氧化锆抛光(低速)



4115A

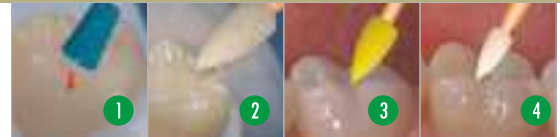Smooth : Sr103M, Sr105M, Sr106M ; Pro-polish : Sr103F, Sr105F, Sr106F
Polish: Sr101EF, Sr103EF, Sr105EF, Sr106EF, Final Polish: Sr101UF, Sr103UF



Composite or filling material Polisher RA
美学树脂、填充材料抛光



CP6001A

3 Step: Sr 204 C x 2 : Smooth;
Sr 204 M x 2 : Pro-polish ; Fine : Sr 204 F x 2 : Polish



Composite Polishing kit , RA
树脂抛光套装(低速)



4110A2

Smooth: FL2-Rwx 2 , CN1-Rw x 2 ; Pro-polish: Cr101M, Cr103M, Cr104M, Cr106M,
Polish: Cr101F, Cr103F, Cr104F, Cr106F

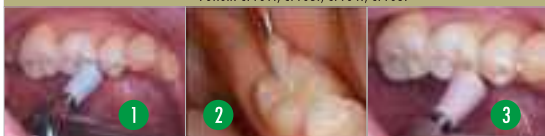


Ceramics / Porcelain Teeth Polish , RA
烤瓷牙抛光(低速)




4112A

Smooth: Cr101C, Cr103C, Cr104C, Cr106C
Pro-polish : Cr101M, Cr103M, Cr104M, Cr106M
Polish: Cr101F, Cr103F, Cr104F, Cr106F



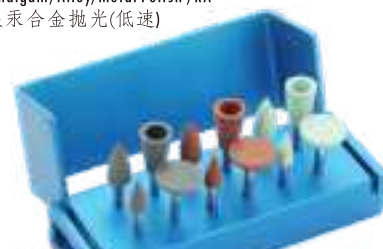Amalgam, Alloy, Metal Polish , RA
银汞合金抛光(低速)



4114A

Smooth : Ar101C, Ar103C, Ar104C, Ar106C
Pro-polish : Ar101M, Ar103M, Ar104M, Ar106M
Polish: Ar101F, Ar103F, Ar104F, Ar106F



Micro Invasive Surgery Burs
微创手术用车针



4105A

BR-49 x 2, SR-41 x 2, SF-41, EX-41, TC-26,
CR-41, FG1/4, FG1/2, FG 1, FO-40EF

Extra Fine Diamond burs for polishing
特细粒度金刚砂抛光车针,高速



4106A

CR-11EF, BR-40EF, TC-27EF, TC-11EF, TF-12EF, TF-11EF,
FO-40EF, TR-11EF, FO-22EF, FO-21EF, FO-30EF, EX-21EF

Composite, Ceramics Finishing kit, FG
树脂成型(高速)



0302A

FL2-Fw x 3; CN1-Fw x 3
FL2-Fg x 3; CN1-Fg x 3

Composite, Ceramics Finishing kit, RA
树脂成型(低速)



0306A

FL2-Rw x 3; CN1-Rw x 3
FL2-Rg x 3; CN1-Rg x 3

White Stone Burs kit, FG
白车石(高速)



302F1

FL2-Fw x 2, CN1-Fw x 2, RD1-Fw x 2
RD2-Fw x 2, TC1-Fw x 2, IC2-Fw x 2

White Stone Burs kit, RA
白车石(低速)



306R1

FL2-Rw x 2, FL3-Rw x 2, CN1-Rw x 2,
RD1-Rw x 2, RD2-Rw x 2, TC1-Rw x 2.

HP Kit to adjust Zirconia
氧化锆修整椅旁套装(HP)



Z31

ZF33, DY001EF to Smoothing
D203F for Pre-polishing; D203EF for Polishing

HP Kit to adjust Lithium Disilicate
玻璃陶瓷修整椅旁套装(HP)



L31

ZF33, DY001EF to Smoothing
D203F for Pre-polishing; D203UF for Polishing

Samples Collection
For Clinic Items

样品册

Flip-Cart Plastic Pages ;
Smart, Easily fill all burs ;
Economic and Recyclable ;
Nice-looking for showing to Customers ;
Convenient to take to anywhere .



No. 170B5

Empty Collection Wallet
样品本(不带车针)



4 Colors Available ,
100 holes per page,
5 pages per set,
Smart to FG or RA Burs



No. 1701D

Diamond Burs Collection Wallet
金刚石车针样品本



More than 300 specification ,
Diamond burs ,
FG, FG XL, FG 2XL,
FG Short, RA, RA Long



No. 1702T

Tungsten Carbide Collection Wallet
钨钢车针样品本



Nearly 200 specification ,
TC burs
FG, FG XL, FG Short, RA, RA Long
T & F Gold, Endo z etc.



No. 1703P

Polishing burs Collection Wallet
抛光器械样品本



Nearly 150 specifications ,
Silicon Rubber burs,
Diamond Abrasives & polishers,
Polishing brushes etc



Four Colors Available

Brown

Black

Blue

Yellow



Content 目录

Laboratory Part 技工类
Page 25-54 第25-54页

Fully sintered
Diamond bur HP
全砂烧结磨头

Page 28-29 第28-29页

CAD / CAM Milling Burs
铣刀

Page 32-34 第32-34页

Tungsten Carbide burs
钨钢车针

Page 35-38 第35-38页

DuCoBur [®] 陶瓷结合剂
金刚石磨头
豆酷巴 (全瓷、全锆磨头)

Diamond Stones & Polishers
Page 39-43 第39至43页

Silicon Polishers
硅胶抛光磨头

Page 44-45 第44-45页

Mounted stone
Sandpaper
车石、砂纸棒

Page 46 第46页

Polishing Brushes ,
Felt Bob & burMandrels
毛毡磨头、夹针、抛光刷

Page 47-49 第47-49页

Diamond burs 金刚石磨头

Page 26-27 第26-27页

Diamond Disc ,
Plaster Trimming wheel
金刚石切片,石膏模修整盘

Page 30-31 第30-31页



Laboratory

Kit 套装 Page 51-53 第51-53页

Samples Collection
样品册

Page 54 第54页

FlexSoft Partial Denture
隐形义齿胶

Page 55 第55页

**Diamond Bur HP
(Electroplated)
金钻石磨头 (电镀)**



Sharp and Durable,
Grind or polish Ceramics
锋利、耐用
磨削平稳、不偏摆
不留黑印、适合车瓷



Grit

- Black / Super Coarse
- Green / Coarse
- Blue / Medium
- Red / Fine

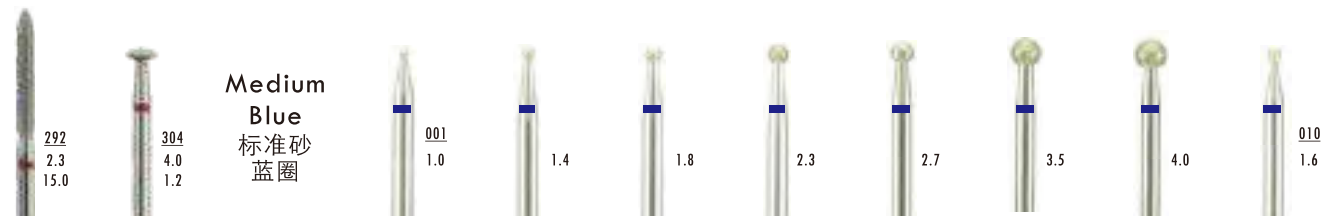
**Fine Red
细砂
红圈**



编号 No. FD30 FK33 FA23 FA33 FC18 FC25 FM27 FL33 FN23 FP33

This grit can grind Zirconia 可用于车氧化锆


**Medium Blue
标准砂
蓝圈**



编号 No.

Round / 圆球

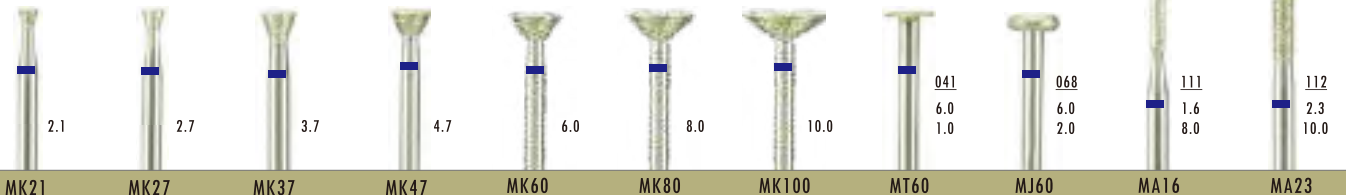
Inverted Cone



MK21 MK27 MK37 MK47 MK60 MK80 MK100 MT60 MJ60 MA16 MA23

Wheel / 平轮

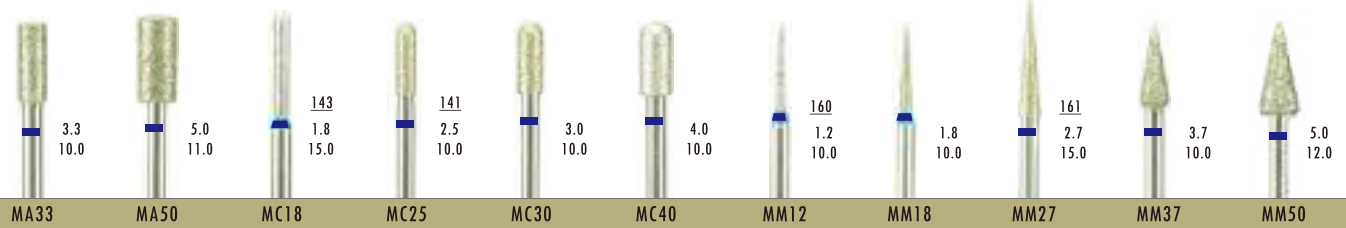
Cylinder / 平柱



MA33 MA50 MC18 MC25 MC30 MC40 MM12 MM18 MM27 MM37 MM50

Round Cylinder / 圆头柱

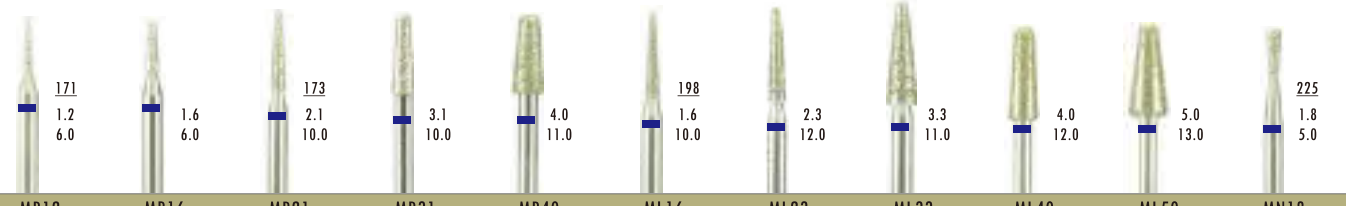Needle / 尖针





MB12 MB16 MB21 MB31 MB40 ML16 ML23 ML33 ML40 ML50 MN18

Taper / 平锥

Round Taper / 圆头锥

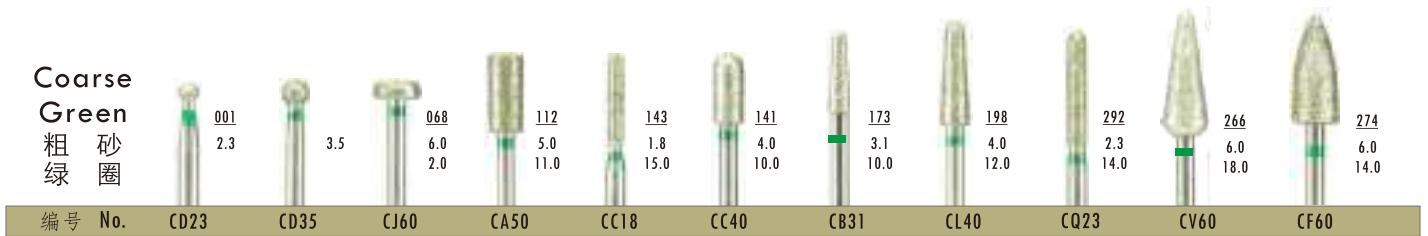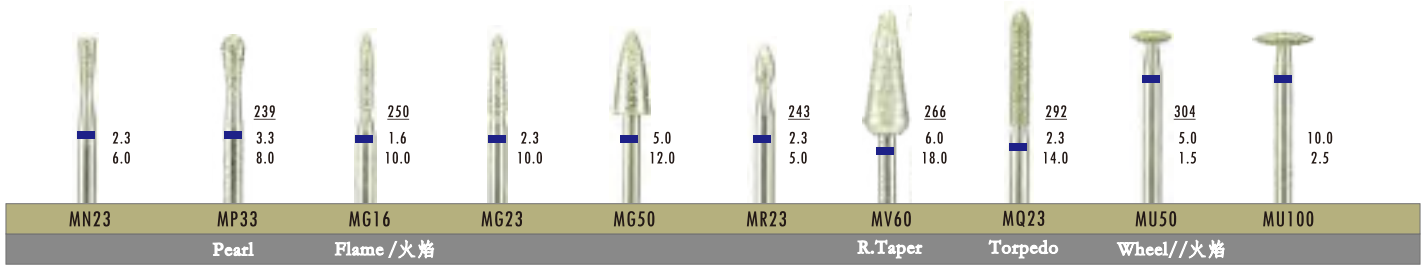Inverted Cone



Transparent Showing plates
展示台 120x120x8mm



No. 1216
7 x 7 holes
1.60mm, for FG burs



No. 1223
7 x 7 holes
2.35mm, for RA or HP burs.

Fully-Sintered Diamond Bur . HP

全砂烧结磨头



Extremely Durable,
Grind or polish smoothly
Especial for Metal

从里到外都又砂
平稳、耐用，适合车金



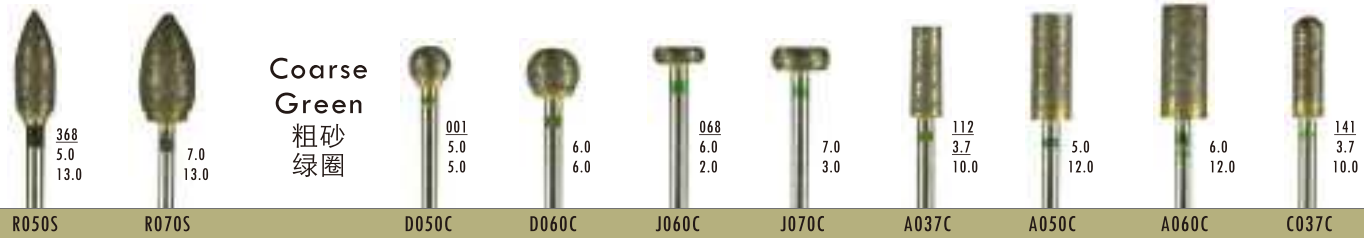
Grit

- Black/Super Coarse
- Green/Coarse
- Blue/Medium
- Red/Fine
- Yellow/Extra Fine


Super Coarse Black 特粗砂黑圈









Coarse Green 粗砂绿圈



Medium Blue 标准砂蓝圈

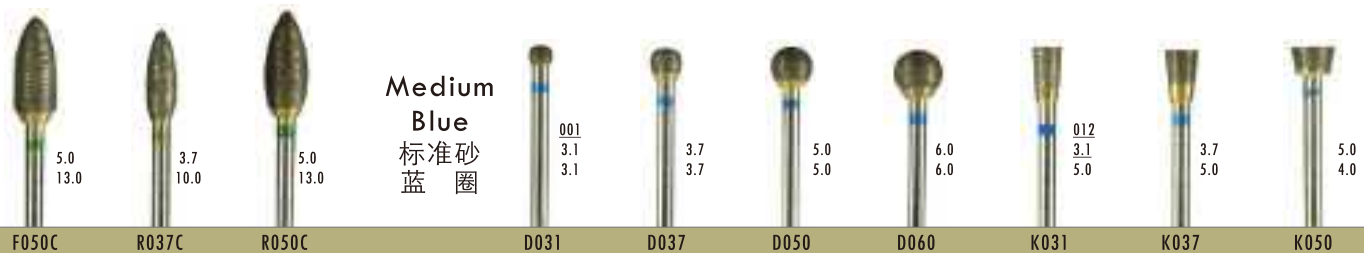

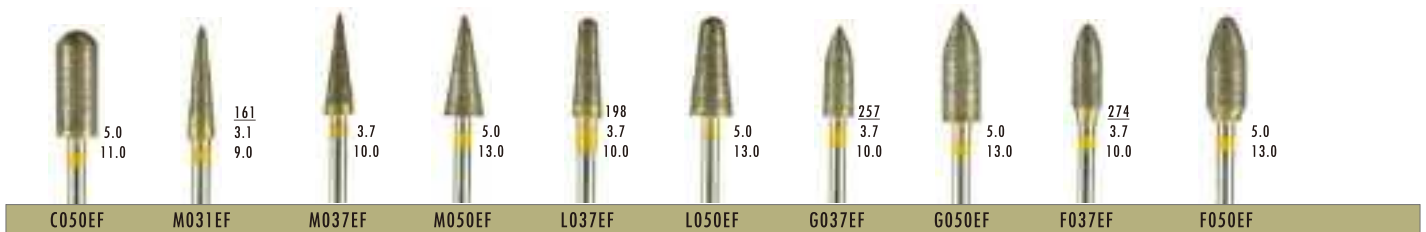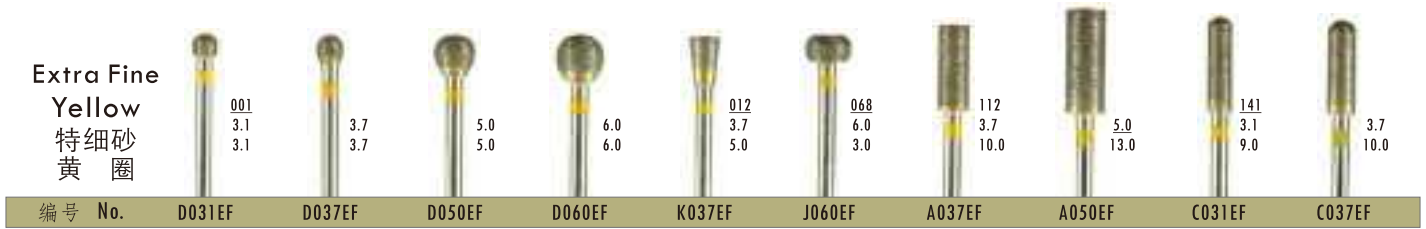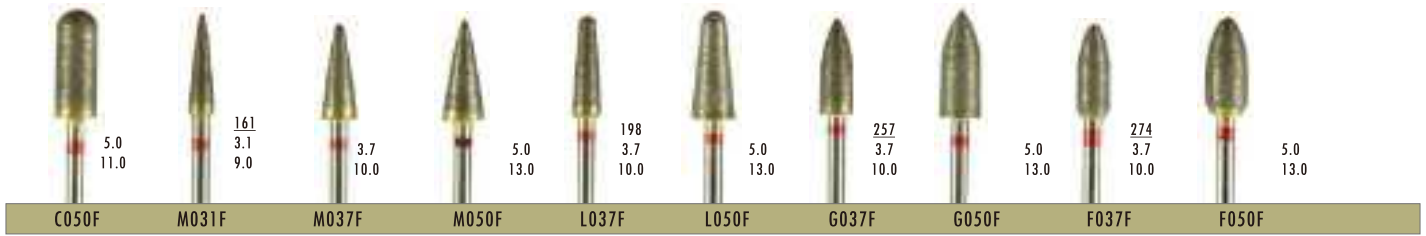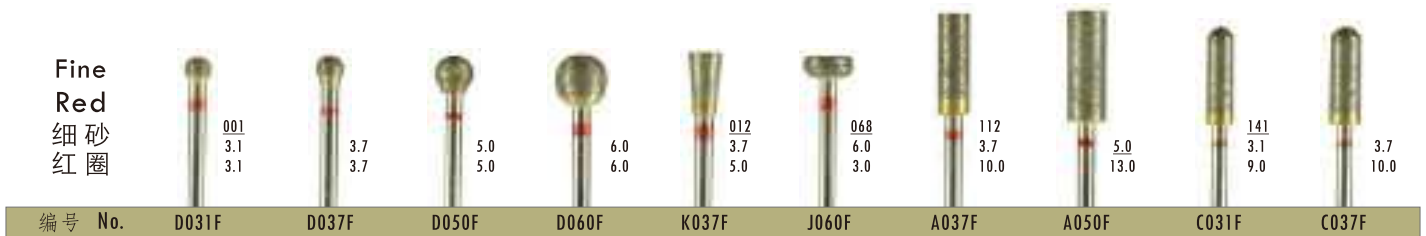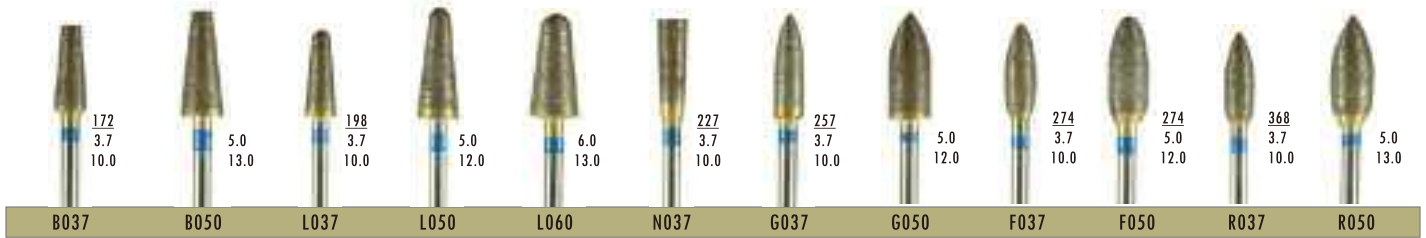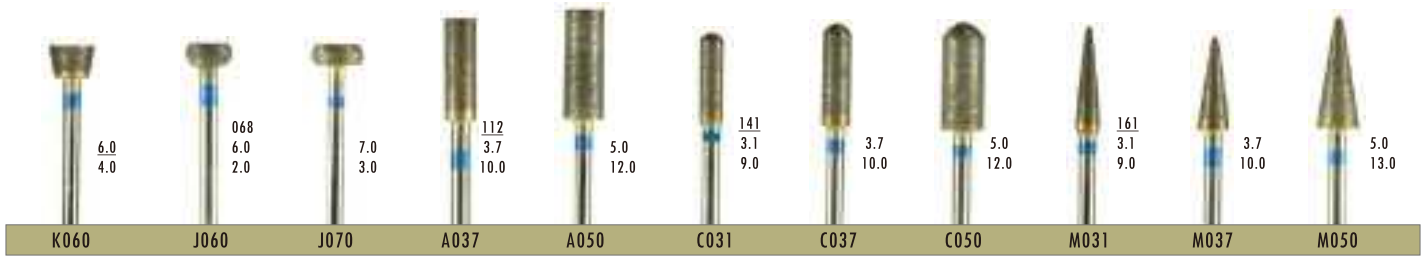


Diamond Disc

金刚石切片



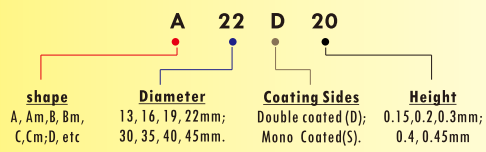How to order these Discs:



Extra Thin, 0.15mm
Double sided
Extra Fine
超薄, 双面砂



编号 No. A19D15 A22D15 Am13D15 Am16D15 Am19D15 Am22D15

Thin, 0.20 x 19mm
Double sided
Fine Grit
较薄, 双面砂



编号 No. A19D20 B19D20 C19D20 N19D20 Am19D20 Bm19D20 Em19D20 Hm19D20

Thin, 0.20 x 22mm
Double sided
Fine Grit
较薄, 双面砂



编号 No. Am22D20 Bm22D20 Cm22D20 Em22D20 Im22D20 Jm22D20 Hm22D20

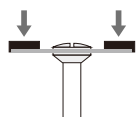


A22D20 B22D20 C22D20 E22D20 F22D20 H22D20 J22D20 N22D20

Hard, 0.3 - 0.4mm
Double sided
Medium Grit
厚硬, 双面砂



编号 No. Am19D30 Am22D30 Am22D40 A19D40 A22D40
19.0 x 0.3 mm 22.0 x 0.3 mm 22.0 x 0.4 mm 19.0 x 0.4 mm 22.0 x 0.4 mm




Single-sided
仅单面镀砂



编号 No. A22S20 Am22S20 Am22S30
22.0 x 0.2 mm, Fine 22.0 x 0.2 mm 22.0 x 0.3 mm, Medium

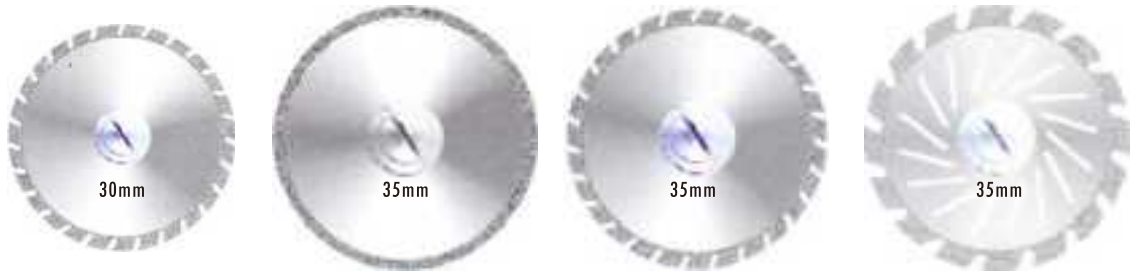
Fast Cutting
Plaster Disc ,
0.3mm,
Double sided
石膏模切片
平整、无偏摆




编号 No.	C30D30	A35D30	C35D30	J35D30
--------	--------	--------	--------	--------



A40D30	C40D30	J40D30
--------	--------	--------

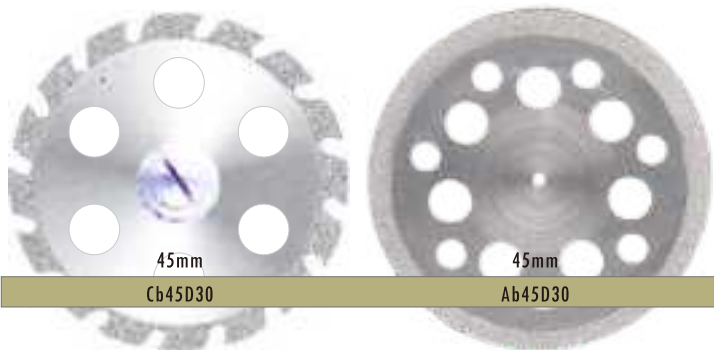

Cb45D30	Ab45D30
---------	---------

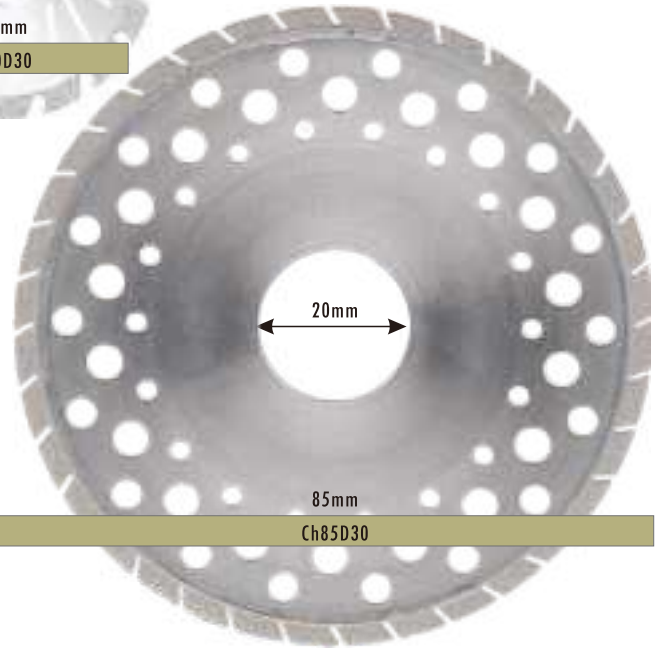

		85mm
--	--	------

Plaster Trimming wheel
石膏修整盘

Cuts Fast and last longer
Ensure immediate removal
Used dry or wet on all materials.

快速磨削,耐用;即装即用
干磨水磨都可以.



R10	R12	S10	C10
10' x 1/4' x 1' @ 250 x 6.3 x 25.4mm	12' x 1/4' x 1' @ 300 x 6.3 x 25.4mm	10' x 1/8' x 1' @ 250 x 3.0 x 25.4mm	10' x 1/4' x 1' @ 250 x 6.3 x 25.4mm
Resin Plastic base 塑料板, 轻便, 节能		Steel base 钢板	Silicon Carbide 碳化硅

Please notice Inner hole diameter ; 1.25' = 32 mm is also alternative .
重要提示: 请注意您设备转轴的直径, 并选择相应的内孔
我们可提供 25mm 和 32mm 的内孔供您选择.

Article No	Coating	D1	L1	L2	L3	D2
for Roland 00						
RL-00-06-ED	Dia.	0.6	9.0	—	40	3.0
RL-00-10-ED	Dia.	1.0	6.0	—	40	3.0
RL-00-25-ED	Dia.	2.5	15.0	—	40	3.0

Article No	Coating	D1	L1	L2	L3	D2
for Amann Girrbach 00						
AG-00-04-ED	Dia.	0.4	6.5	—	43	3.0
AG-00-10-ED	Dia.	1.0	14.0	—	43	3.0
AG-00-14-ED	Dia.	1.4	14.0	—	43	3.0
AG-00-18-ED	Dia.	1.8	14.0	—	43	3.0

Article No	Coating	D1	L1	L2	L3	D2
for VHF						
VH-00-06-ED	Dia.	0.6	6.0	—	35	3.0
VH-00-10-ED	Dia.	1.0	9.0	—	35	3.0
VH-00-25-ED	Dia.	2.5	16.0	—	35	3.0

Article No	Coating	D1	L1	L2	L3	D2
for Imes-Icore 350, 550, 750						
II-75-06-ED	Dia.	0.6	10.0	—	40	6.0
II-75-10-ED	Dia.	1.0	10.0	—	40	6.0
II-75-25-ED	Dia.	2.5	15.0	—	40	6.0

Article No	Coating	D1	L1	L2	L3	D2
for Imes-Icore 240, 250						
II-24-06-ED	Dia.	0.6	10.0	—	43	3.0
II-24-10-ED	Dia.	1.0	10.0	—	45	3.0
II-24-25-ED	Dia.	2.5	15.0	—	42	3.0

CAD / CAM Milling Burs for Lithium Disilicate 铣刀(加工玻璃陶瓷)

D1 : Cutting Edge Diameter (.mm) **L1** : Cutting Edge Length (.mm)
D2 : Shank Diameter (.mm) **L2** : Free Grinding to (.mm)
L3 : Total Length (.mm)

Article No	Coating	D1	L1	L2	L3	D2
for Sirona MCX5						
SR-X5-06-ED	Dia.	0.6	14.5	—	34.5	3.0
SR-X5-12-ED	Dia.	1.2	12.5	—	36.5	3.0
SR-X5-14-ED	Dia.	1.4	14.5	—	35.5	3.0
SR-X5-22-ED	Dia.	2.2	15.5	—	37.5	3.0

Article No	Coating	D1	L1	L2	L3	D2
for Sirona MCXL						
SR-XL-C125-ED	Dia.	—	13.7	—	38	3.5
Cylinder Pointed Bur 12S						
SR-XL-S125-ED	Dia.	—	13.7	—	38	3.5
Step Bur 12S						
SR-XL-S12-ED	Dia.	—	13.2	—	38	3.5
Step Bur 12						

Article No	Coating	D1	L1	L2	L3	D2
for Sirona MCXL						
SR-XL-CP12EF-ED	Dia.	—	19.0	—	40	3.5
Cylinder Pointed Bur 12EF						
SR-XL-C12EF-ED	Dia.	—	19.0	—	40	3.5
Cylinder Bur 12EF						

Article No	Coating	D1	L1	L2	L3	D2
for Sirona MCXL						
SR-XL-C20-ED	Dia.	—	21.5	—	46	3.5
Cylinder Pointed Bur 20						
SR-XL-S20-ED	Dia.	—	22.0	—	46.3	3.5
Step Bur 20						

Article No	Coating	D1	L1	L2	L3	D2
for Sirona Cerec 3						
SR-C3-CP-ED	Dia.	—	9.3	—	23.5	1.8
Cylinder Pointed Bur						
SR-C3-12S-ED	Dia.	—	13.2	—	25.2	1.8
Step Bur 12						
SR-C3-10S-ED	Dia.	—	11.5	—	23.3	1.8
Step Bur 10						

Article No	Coating	D1	L1	L2	L3	D2
for Dentium						
DT-00-06-ED	Dia.	0.6	19.8	—	50	6.0
DT-00-10-ED	Dia.	1.0	19.8	—	50	6.0
DT-00-20-ED	Dia.	2.0	19.8	—	50	6.0


Article No	Coating	D1	L1	L2	L3	D2
for Dentium						
DT-01-10-ED	Dia.	1.0	10.0	—	45	3.0
DT-01-20-ED	Dia.	2.0	12.0	—	45	3.0

Article No	Coating	D1	L1	L2	L3	D2
for Arum						
AR-00-10-ED	Dia.	1.0	10.5	—	50	6.0
AR-00-15-ED	Dia.	1.5	12.5	—	50	6.0
AR-00-20-ED	Dia.	2.0	12.5	—	50	6.0
AR-00-25-ED	Dia.	2.5	12.5	—	50	6.0

Article No	Coating	D1	L1	L2	L3	D2
for Arum						
AR-01-06-ED	Dia.	0.6	12.5	—	41	4.0
AR-01-10-ED	Dia.	1.0	16.0	—	44.5	4.0
AR-01-15-ED	Dia.	1.5	16.0	—	44.5	4.0
AR-01-25-ED	Dia.	2.5	16.0	—	44.5	4.0

CAD / CAM Milling Burs.

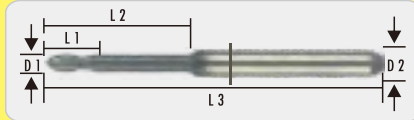
铣刀



D1 : Cutting Edge Diameter (.mm) **L1** : Cutting Edge Length (.mm)
D2 : Shank Diameter (.mm) **L2** : Free Grinding Length (.mm)
L3 : Total Length (.mm)

Coating

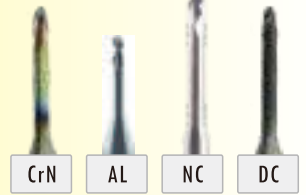CrN : for PMMA
 AL : for Metal
 NC : for Zirconia
 DC : for Zirconia,
 More Durable.



Article No	Coating	D1	L1	L2	L3	D2
for Roland 50						
RL-50-03-NC	NC	0.3	1.5	6.0	50.0	4.0
RL-50-06-CRN	CrN	0.6	1.5	6.0	50.0	4.0
RL-50-06-NC	NC					
RL-50-10A-CRN	CrN	1.0	3.5	16.0	50.0	4.0
RL-50-10A-NC	NC					
RL-50-20A-CRN	CrN	2.0	6.0	16.0	50.0	4.0
RL-50-20A-NC	NC					
RL-50-20A-DC	DC					
RL-50-10B-CRN	CrN	1.0	3.5	20.0	50.0	4.0
RL-50-10B-NC	NC					
RL-50-20B-CRN	CrN	2.0	6.0	25.0	52.0	4.0
RL-50-20B-NC	NC					
RL-50-20B-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for 3 M ESPE LAVA						
LV-00-10-CRN	CrN	1.0	3.5	16.0	38.0	3.0
LV-00-10-NC	NC					
LV-00-20-CRN	CrN	2.0	6.0	16.0	38.0	3.0
LV-00-20-NC	NC					
LV-00-20-DC	DC					
LV-01-10-CRN	CrN	1.0	2.0	25.0	50.0	3.0
LV-01-10-NC	NC					
LV-01-20-CRN	CrN	2.0	2.0	25.0	50.0	3.0
LV-01-20-NC	NC					
LV-01-20-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for RK						
RK-00-10-DC	DC	1.0	3.5	16.0	47.0	3.0
RK-00-20-DC	DC	2.0	3.5	16.0	47.0	3.0

Article No	Coating	D1	L1	L2	L3	D2
for Imes-Icore 240, 250						
II-24-06-CRN	CrN	0.6	1.5	10.0	48.0	3.0
II-24-06-NC	NC					
II-24-10-CRN	CrN	1.0	5.0	14.0	48.0	3.0
II-24-10-NC	NC					
II-24-25-CRN	CrN	2.5	8.0	20.0	48.0	3.0
II-24-25-NC	NC					
II-24-25-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Imes-Icore 350, 550, 750						
II-75-06-CRN	CrN	0.6	2.5	12.0	53.0	6.0
II-75-06-NC	NC					
II-75-10-CRN	CrN	1.0	6.0	14.0	53.0	6.0
II-75-10-NC	NC					
II-75-25-CRN	CrN	2.5	10.0	20.0	53.0	6.0
II-75-25-NC	NC					
II-75-25-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Amann Girrbach						
AG-00-06-CRN	CrN	0.6	1.5	13.0	47.0	3.0
AG-00-06-NC	NC					
AG-00-10-CRN	CrN	1.0	3.5	16.0	47.0	3.0
AG-00-10-NC	NC					
AG-00-25-CRN	CrN	2.5	5.0	16.0	47.0	3.0
AG-00-25-NC	NC					
AG-00-25-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Amann Girrbach - Metal mill						
AG-01-10-AL	Alnova	1.0	5.0	7.0	43.0	3.0
AG-01-20-AL	Alnova	2.0	5.0	7.0	43.0	3.0

Article No	Coating	D1	L1	L2	L3	D2
for VHF K1, K4						
VH-K4-06-CRN	CrN	0.6	2.0	7.0	35.0	3.0
VH-K4-06-NC	NC					
VH-K4-10-CRN	CrN	1.0	3.5	16.0	35.0	3.0
VH-K4-10-NC	NC					
VH-K4-20-CRN	CrN	2.0	5.0	16.0	35.0	3.0
VH-K4-20-NC	NC					

Article No	Coating	D1	L1	L2	L3	D2
for VHF K5						
VH-K5-03-NC	CrN	0.3	1.2	13.0	40.0	3.0
VH-K5-06-CRN	CrN	0.6	1.5	7.0	40.0	3.0
VH-K5-06-NC	NC					
VH-K5-10-CRN	CrN	1.0	3.5	14.0	40.0	3.0
VH-K5-10-NC	NC					
VH-K5-20-CRN	CrN	2.0	4.5	16.0	40.0	3.0
VH-K5-20-NC	NC					
VH-K5-20-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Dentium						
DT-00-10-NC	NC	1.0	3.5	16.0	50.0	6.0
DT-00-20-NC	NC	2.0	3.5	16.0	50.0	6.0
DT-00-20-DC	DC					
DT-01-05-NC	NC	0.5	2.0	8.0	43.0	3.0
DT-01-10-CrN	CrN	1.0	3.5	16.0	50.0	3.0
DT-01-10-NC	NC					
DT-01-20-CrN	CrN	2.0	3.5	18.0	50.0	3.0
DT-01-20-NC	NC					
DT-01-20-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Zirkon Zahn M3						
ZZ-03-05-CrN	CrN	0.5	1.5	6.0	57.0	3.0
ZZ-03-05-NC	NC					
ZZ-03-10-CrN	CrN	1.0	6.0	12.0	57.0	3.0
ZZ-03-10-NC	NC					
ZZ-03-20-CrN	CrN	2.0	10.0	18.0	57.0	3.0
ZZ-03-20-NC	NC					
ZZ-03-20-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Zirkon Zahn M5, M6, etc						
ZZ-06-05-NC	NC	0.5	1.5	3.0	50.0	6.0
ZZ-06-10-NC	NC	1.0	6.0	12.0	50.0	6.0
ZZ-06-20-NC	NC	2.0	10.0	18.0	50.0	6.0
ZZ-06-20-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Sirona MX 5						
SR-MX-05-NC	NC	0.5	1.5	5.0	42.0	3.0
SR-MX-10-NC	NC	1.0	4.5	18.0	43.0	3.0
SR-MX-25-NC	NC	2.5	7.0	25.0	43.0	3.0
SR-MX-25-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Cercon						
CC-00-05-NC	NC	0.5	2.0	10.0	55.0	3.5
CC-00-10-NC	NC	1.0	8.0	30.0	55.0	3.5
CC-00-20-NC	NC	2.0	10.0	30.0	55.0	3.5

Article No	Coating	D1	L1	L2	L3	D2
for Arum						
AR-00-10-CrN	CrN	1.0	2.0	16.0	55.0	6.0
AR-00-10-NC	NC					
AR-00-20-CrN	CrN	2.0	6.0	20.0	55.0	6.0
AR-00-20-NC	NC					
AR-00-20-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Wieland Mini						
WL-MI-07-CrN	CrN	0.7	1.5	7.0	35.0	3.0
WL-MI-07-NC	NC					
WL-MI-10-CrN	CrN	1.0	3.5	16.0	35.0	3.0
WL-MI-10-NC	NC					
WL-MI-25-CrN	CrN	2.5	3.5	20.0	35.0	3.0
WL-MI-25-NC	NC					

Article No	Coating	D1	L1	L2	L3	D2
for Wieland S1						
WL-S1-07-CrN	CrN	0.7	1.5	7.0	40.0	3.0
WL-S1-07-NC	NC					
WL-S1-10-CrN	CrN	1.0	3.5	16.0	40.0	3.0
WL-S1-10-NC	NC					
WL-S1-25-CrN	CrN	2.5	5.0	20.0	40.0	3.0
WL-S1-25-NC	NC					
WL-S1-25-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Kavo						
KA-00-10-DC	DC	1.0	2.0	16.0	53.0	6.0
KA-00-20-DC	DC	2.0	6.0	19.5	53.0	6.0

Article No	Coating	D1	L1	L2	L3	D2
for Yena						
YN-00-10-NC	NC	1.0	3.5	16.0	45.0	4.0
YN-00-20-NC	NC	2.0	3.5	16.0	45.0	4.0
YN-00-20-DC	DC					

Article No	Coating	D1	L1	L2	L3	D2
for Arum - Metal Mill						
AR-01-10-AL	Alnova	1.0	2.0	10.0	50.0	6.0
AR-01-15-AL	Alnova	1.5	3.0	10.0	50.0	6.0
AR-01-20-AL	Alnova	2.0	4.0	12.0	50.0	6.0
AR-01-30-AL	Alnova	3.0	6.0	12.0	50.0	6.0













Tungsten Carbide Burs . HP













钨钢车针



























Shape













- S1, ISO 001, Round
- S2, ISO 002, Inverted cone
- S21, S31, ISO 107, Cylinder
- S33, ISO 168, Taper
- S31R, ISO 137, Round, Cylinder













											
S1	001 0.5	001 0.6	001 0.7	001 0.8	001 0.9	001 1.0	001 1.2	001 1.4	001 1.6	001 1.8	
球型/Round											
编号 No.	1005	1006	1007	1008	1009	1010	1012	1014	1016	1018	
参考号 Ref.	HP1/4	HP1/2	HP1s	HP 1	HP 1b	HP 2	HP 3	HP 4	HP 5	HP 6	

											
001 2.1	001 2.3	001 2.5	001 2.7	001 3.1		010 0.8	010 1.0	010 1.2	010 1.4	010 1.6	
					倒锥/Inv. Cone						
1021	1023	1025	1027	1031		2008	2010	2012	2014	2016	
HP 7	HP 8	HP 9	HP 10	HP 12		HP 34	HP 35	HP 36	HP 37	HP 38	


											
010 1.8	010 2.1	010 2.3		107 0.8	107 0.9	107 1.0	107 1.2	107 1.4	107 1.6	107 1.8	
			圆柱/Cylinder								
2018	2021	2023		21008	21009	21010	21012	21014	21016	21018	
HP 39	HP 40	HP 41		HP 55	HP 56	HP 57	HP 58	HP 59	HP 60	HP 61	

											
107 2.1	107 2.3		107 0.8	107 0.9	107 1.0	107 1.2	107 1.4	107 1.6	107 1.8	107 2.1	
		圆柱/Cylinder									
21021	21023		31008	31009	31010	31012	31014	31016	31018	31021	
HP 62	HP 63		HP 555	HP 556	HP 557	HP 558	HP 559	HP 560	HP 561	HP 562	

											
107 2.3		168 0.8	168 0.9	168 1.0	168 1.2	168 1.6	168 2.1	168 2.3			171 0.9
	平锥/Taper									平锥加长/Taper L	
31023		33008	33009	33010	33012	33016	33021	33023			33L009
HP 563		HP 698	HP 699	HP 700	HP 701	HP 702	HP 703	HP 704			HP 699 L

											
171 1.0	171 1.2	171 1.6	171 2.1		137 1.0	137 1.2	137 1.4				
				圆头柱 Round Cyl.							
33L010	33L012	33L016	33L021		31R010	31R012	31R014				
HP 700 L	HP 701 L	HP 702 L	HP 703 L		HP 1557	HP 1558	HP 1559				

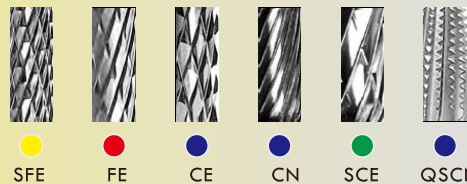Tungsten Carbide Cutter . HP 钨钢磨头



Super Fine Cross Cut (SFE) 特细双纹 黄圈 ISO 500.104.110.



编号 No.	C023SFE	F040SFE	F060SFE	G023SFE	L016SFE	L023SFE	L040SFE	L060SFE	M060SFE	P023SFE
参考号 Ref.	0230.023	7230.040	5430.060	1030.023	1230.016	0130.023	5730.040	5530.060	5930.060	7730.023

Fine Cross Cut (FE) 较细双纹 红圈 ISO 500.104.140.



编号 No.	Q023SFE	Q023SFE	C023FE	C060FE	F060FE	G040FE	G060FE	L016FE	L023FE
参考号 Ref.	0830.023	0430.023	0220.023	5820.023	5420.060	6120.040	5320.060	1220.016	0120.023

Medium Cross Cut (CE) 标准双纹 蓝圈 ISO 500.104.190.





编号 No.	L040FE	L060FE	M060FE	P023FE	P060FE	Q023FE	A023CE	C023CE	C060CE
参考号 Ref.	5720.040	5520.060	5920.060	7720.023	5120.060	0420.023	0310.023	0210.023	5810.060





编号 No.	D023CE	D030CE	D040CE	D050CE	F040CE	F060CE	G023CE	G040CE	G060CE	L023CE
参考号 Ref.	7110.023	7110.030	7110.040	7110.050	7210.040	5410.060	1010.023	6110.040	5310.060	0110.023

								
201 040 14.0	201 060 15.0	167 023 14.50	167 040 14.5	161 060 15.0	237 023 5.0	237 060 11.0	292 023 16.0	001 030 3.0
L040CE	L060CE	M023CE	M040CE	M060CE	P023CE	P060CE	Q023CE	D030CN
5710.040	5510.060	—	—	5910.060	7710.023	5110.060	0410.023	7170.030

Medium Plain Cut (CN)
标准单纹
蓝圈
ISO 500.104.175.

									
001 040 4.0	001 050 5.0	201 040 14.0	201 060 15.0	161 060 14.0	257 060 12.0	274 060 14.0	201 060 14.0	274 060 14.0	257 060 12.0
D040CN	D050CN	L040CN	L060CN	M060CN	G060CN	F060CN	C060SCE	F060SCE	G060SCE
7170.040	7170.050	5770.040	5570.060	5970.060	5370.060	5470.060	5850.060	5450.060	5350.060

Coarse Cross Cut (SCE)
粗糙双纹
绿圈
ISO 500.104.220

							
201 040 14.0	201 060 14.0	237 060 11.0	263 070 14.0	274 060 15.0	201 023 14.5	200 045 13.0	263 060 12.0
L040SCE	L060SCE	P060SCE	L070XCE	F060QSCN	L023QSCN	L045QSCN	LB060QSCN
5750.040	5550.060	5150.060	6055.070	7275.060	0175.023	5675.045	5275.060








Super Coarse Cross Cut (XCE)
特粗双纹
黑圈
ISO 500.104.223

Straight Edge with Spiral Cut (QSCN)
直刃螺旋纹
蓝圈
ISO 500.104.176

							
141 023 8.0	194 023 8.0	201 040 14.0	289 023 8.0	274 040 11.5	274 060 15.0	201 023 14.5	201 040 13.0
CS023GTi	LS023GTi	L040GTi	QS023GTi	F040QF	F060QF	L023QF	L040QF
1726.023	0726.023	6926.040	0826.023	7224.040	5424.060	0124.023	5724.040

Straight Edge with Cross Cut (GTi)
直刃双纹
黑圈
ISO 500.104.194

Spiral Cut Fine (QF)
较细螺旋纹
红圈
ISO 500.104.134

						
141 023 8.0	274 060 14.0	198 023 8.0	201 040 14.0	289 023 8.0	118 060 3.5	467 023 9.0
CS023QSF	F060QSF	LS023QSF	L040QSF	QS023QSF	ES801	AL023S4
1725.023	5425.060	0725.023	6925.040	0825.023	118.060	467.023

Spiral Cut Extra Fine (QSF)
特细螺旋纹
黄圈
ISO 500.104.137

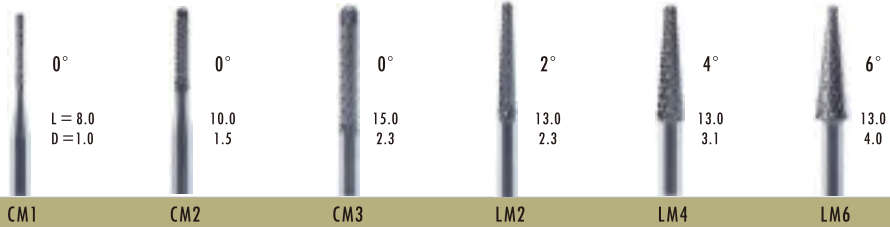
Tungsten Carbide Burs for Milling Machines

钨钢车针(精密研磨仪用)




Precisely match Various System
 Extra - long working life & high cutting capacity ;
 Used for trimming, polishing implant teeth shoulders ;
 Milling precious or Non-precious metals or Titanium .


Medium
 Cross Cut
 (CM, LM)
 标准双纹
 ISO 500.103.190



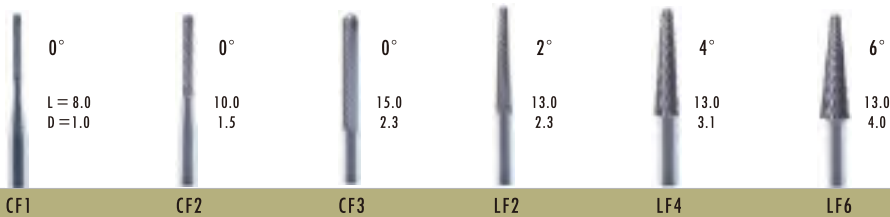
Fine
 Cross Cut
 (CF, LF)
 较细单纹
 ISO 500.103.135



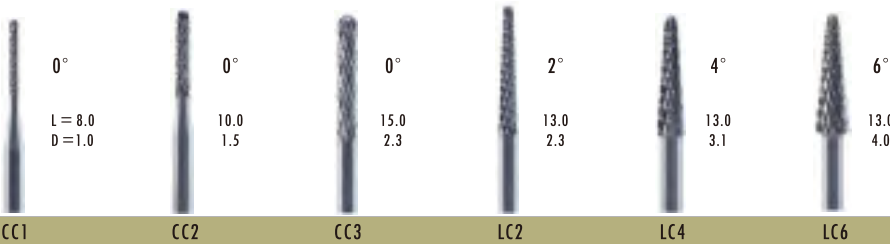Coarse
 Cross Cut
 (CC, LC)
 粗糙双纹
 ISO 500.103.220



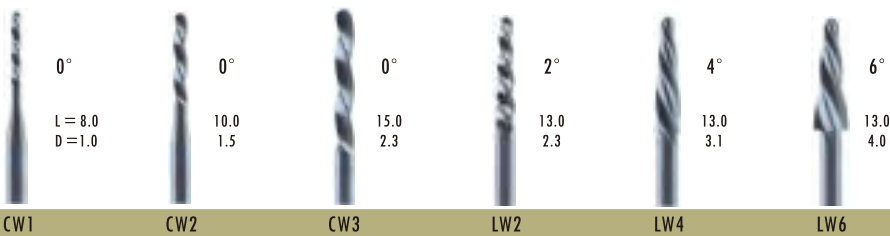Max Mill
 (CW, LW)
 蜡型针
 ISO 500.103.217



DuCoBur® 陶瓷结合剂金刚石磨头 (全瓷、全锆磨头)

Diamond Stones for Zirconia or AllCeramics (Ceramic Pregnant Diamond Grinding Burs)



Widely Used for Zirconia, Lithium Disilicate, Ceramics etc.
COMFORTABLE at the minimal pressure;
HEATLESS no need water cooling ;
MORE DURABLE under our suggested PRM

不产热、无需水冷 更耐磨、经济实惠 更锋利、手感好 不绷瓷、更轻松 型号齐、用途广、啥都磨得动

Medium Blue
蓝色、中等粒度



Item No.



DC009 DC010 DC010B DC011 DC014 DC015 DC016 DC017 DC003A

Size (mm) 4 x 4 4 x 7 4 x 11 3 x 7 2.5 x 4.5 3.5 x 10 2.5 x 7 3.5 x 11 6.5 x 2



DC001 DC002 DC002B DC012 DC012B DC013 DC004A

4 x 13 5 x 13 6.5 x 13 13 x 2 13 x 3.5 13 x 2 16 x 2/0.7



DC005A DC004 DC005 DC006 DC007

16 x 2 22 x 2/0.7 22 x 2 22 x 4.5 15 x 3.5

Medium Black
黑色、中等粒度



Item No.



DX001 DX002 DX002B DX012 DX012B DX013 DX004A

Size (mm) 4 x 13 5 x 13 6.5 x 13 13 x 2 13 x 3.5 13 x 2 16 x 2/0.7



DX005A DX004 DX005 DX006 DX007

16 x 2 22 x 2/0.7 22 x 2 22 x 4.5 15 x 3.5

Medium Yellow

黄色、中等粒度



Item No.	DY001	DY002	DY002B	DY012	DY012B	DY013	DY004A
Size (mm)	4 x 13	5 x 13	6.5 x 13	13 x 2	13 x 3.5	13 x 2	16 x 2/0.7



DY005A	DY004	DY005	DY006	DY007
16 x 2	22 x 2/0.7	22 x 2	22 x 4.5	15 x 3.5

Medium Green


绿色、中等粒度



Item No.	DG009	DG010	DG010B	DG011	DG014	DG015	DG016	DG017	DG003A
Size (mm)	4 x 4	4 x 7	4 x 11	3 x 7	2.5 x 4.5	3.5 x 10	2.5 x 7	3.5 x 11	6.5 x 2



DG001	DG002	DG002B	DG012	DG012B	DG013	DG004A
4 x 13	5 x 13	6.5 x 13	13 x 2	13 x 3.5	13 x 2	16 x 2/0.7

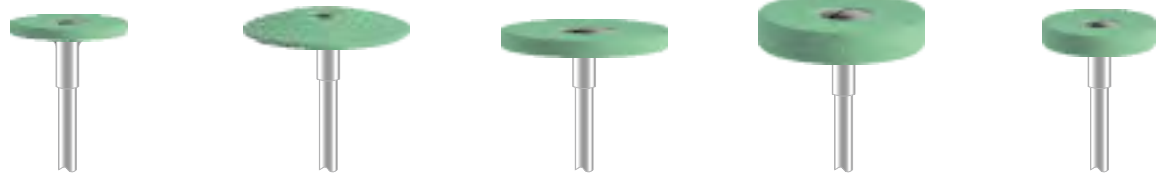


DG005A	DG004	DG005	DG006	DG007
16 x 2	22 x 2/0.7	22 x 2	22 x 4.5	15 x 3.5

Extra Fine
Yellow
Gold shank

Suitable to adjusting Lithium Disilicate, All Ceramics, Zirconia of Occlusal, Lingual, Edge, Intracoronal, Proximal etc. Dentists will feel very convenient if they can prepare a few pieces beside the chairs.

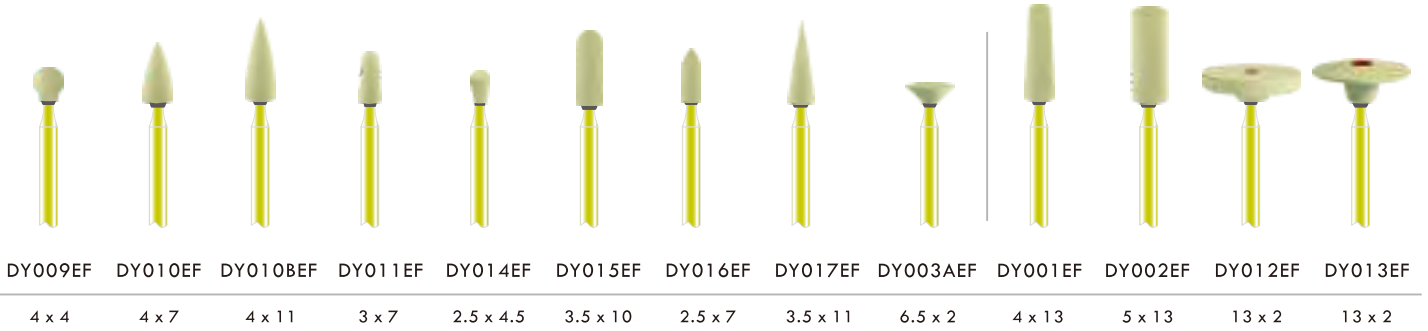专门针对诊所的全瓷、玻璃陶瓷打磨抛光、氧化锆牙调咬合、发育沟、窝沟、边缘、内冠、邻面、舌侧等。黄色、特细粒度。粒度细、不绷瓷。



DuCoBur® 豆酷巴
橡胶结合剂金刚石磨头
(全瓷、全锆抛光磨头)

Diamond Polishers for Zirconia or AllCeramics
(Rubber Pregnant Diamond Polishing Burs)



Widely Used for Zirconia, Lithium Disilicate, Ceramics etc.

COMFORTABLE at the minimal pressure;
HEATLESS no need water cooling ;
MORE DURABLE under our suggested PRM

不产热、无需水冷 更耐磨、经济实惠 手感好、很服帖 抛就亮、更轻松 型号齐、用途广、啥都抛得亮

Super Coarse
Black
黑色、特粗粒度



Item No.



D1S

D2AS

D2S

D3S

D3BS

Size (mm)

26 x 2

20 x 2

26 x 2

22 x 3

26 x 2

Coarse
Green
绿色、粗糙粒度



Item No.



D1C

D2AC

D2C

D3C

D3BC

Size (mm)



26 x 2

20 x 2

26 x 2


22 x 3

26 x 2

 Medium Blue
 蓝色、中等粒度
 1



Item No.	D1M	D2AM	D2M	D3M	D3BM
Size (mm)	26 x 2	20 x 2	26 x 2	22 x 3	26 x 2

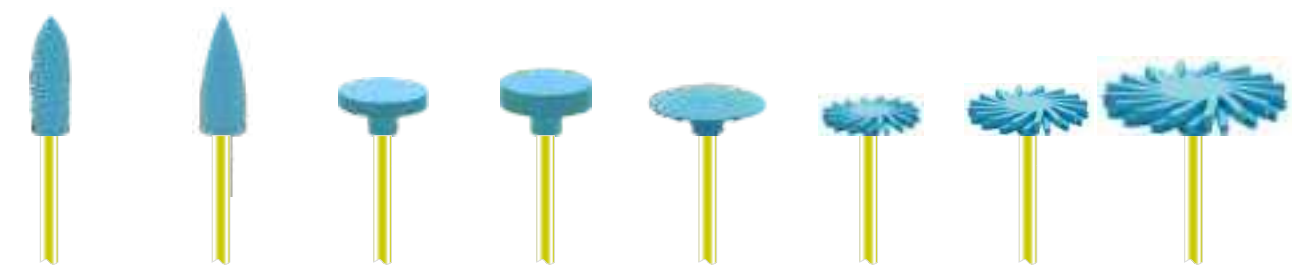


D201M	D203M	D205M	D205BM	D207M	D13M	D17M	D26M
5 x 16	6 x 16	12 x 2	12 x 4	12 x 2.5	14 x 2	17 x 2.2	26 x 2.2

 Fine Red
 红色、较细粒度
 1



Item No.	D1F	D2AF	D2F	D3F	D3BF
Size (mm)	26 x 2	20 x 2	26 x 2	22 x 3	26 x 2

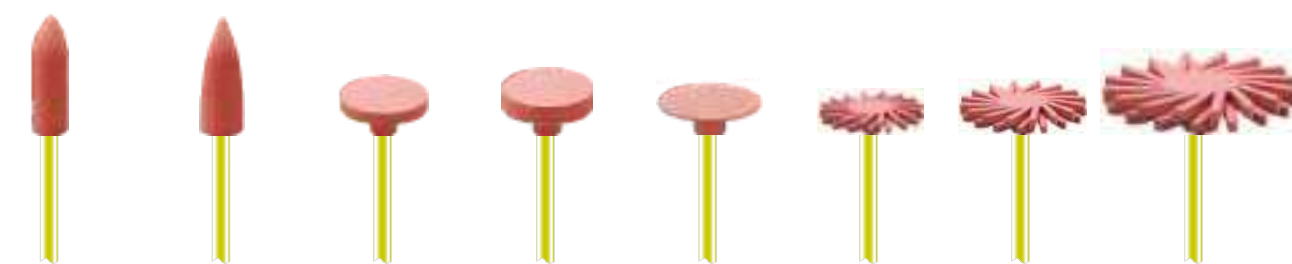

D201F	D203F	D205F	D205BF	D207F	D13F	D17F	D26F
5 x 16	6 x 16	12 x 2	12 x 4	12 x 2.5	14 x 2	17 x 2.2	26 x 2.2

Extra Fine
Yellow
黄色、特细粒度



Item No.	D1EF	D2AEF	D2EF	D3EF	D3BEF
Size (mm)	26 x 2	20 x 2	26 x 2	22 x 3	26 x 2



D201EF	D203EF	D205EF	D205BEF	D207EF	D13EF	D17EF	D26EF
5 x 16	6 x 16	12 x 2	12 x 4	12 x 2.5	14 x 2	17 x 2.2	26 x 2.2

Ultra Fine
Grey
灰白色、特细粒度



Item No.	D1UF	D2AUF	D2UF	D3UF	D3BUF
Size (mm)	26 x 2	20 x 2	26 x 2	22 x 3	26 x 2



D201UF	D203UF	D205UF	D205BUF	D207UF	D13UF	D17UF	D26UF
5 x 16	6 x 16	12 x 2	12 x 4	12 x 2.5	14 x 2	17 x 2.2	26 x 2.2

**Silicon Rubber Wheels ,
Points & Burs
For Polish Ceramics
胶轮、胶粒、硅胶磨头
(烤瓷牙抛光)**



White, Coarse , to adjust and remove Porcelain
白色, 粗糙, 用于烤瓷工件的修整。
Pink, Medium , to Smooth or polish Porcelain,
粉色, 中等, 用于烤瓷毛面平整和预抛光。
Blue, Fine , to high-lustre polishing Porcelain,
蓝色, 较细, 用于烤瓷工件的细抛光。

Coarse
White
白色、粗糙



100

Item No.

C 104 C

C 107 C

C 101 C

C203C

C201C

C205C

C207C

Size (mm)

22x 3

22 x 3

6 x 22

6 x 13

5 x 13

12 x 2

13 x 2

Medium
Pink
粉色、中等粒度



100

Item No.

C 104 M

C 107 M

C 101 M

C203M

C201M

C205M

C207M

Size (mm)

22x 3

22 x 3

6 x 22

6 x 13

5 x 13

12 x 2

13 x 2

Fine
Blue
蓝色、较细



100

Item No.

C 104 F

C 107 F

C 101 F

C203F

C201F

C205F

C207F

Size (mm)

22x 3

22 x 3

6 x 22

6 x 13

5 x 13

12 x 2

13 x 2

**Rubber Wheels , Points & Burs
For Polish Metal or Alloy
胶轮、胶粒、橡胶磨头
(金属、合金抛光)**



Black, Coarse , to adjust Cobalt-Chromium
Alloy and hard metal
黑色, 粗糙, 用于磨削钴铬合金和硬金属。
Brown, Medium , to remove metal and
semi-precious metal
棕色, 中等, 用于烤瓷毛面平整和预抛光。
Green, Fine , to fast smooth scratch or polish metal,
绿色, 较细, 用于快速去除划痕和抛光金属。

Coarse
Black
黑色、较粗



100

Item No.

A102 C

A103C

A104 C

A203C

A201C

A206C

A207C

Size (mm)

22x 3

22 x 3

6 x 22

4 x 12

5 x 13

12 x 2

13 x 2

Medium Brown
棕色、中等粒度
100



Item No.	A102 M	A103 M	A104 M	A203M	A201M	A206M	A207M
Size (mm)	22x 3	22 x 3	6 x 22	4 x 12	5 x 13	12 x 2	13 x 2

Fine Green
绿色、较细
100



Item No.	A102 F	A103 F	A104 F	A203F	A201F	A206F	A207F
Size (mm)	22x 3	22 x 3	6 x 22	4 x 12	5 x 13	12 x 2	13 x 2

Silicon Polishing Burs
For Composites, Resins
硅胶磨头
(树脂胶托抛光)



Green, Coarse, to reduce, finish Acrylic, Resin.
绿色, 粗糙, 用于修整和成型树脂胶托。
Gray, Medium, to Smooth Acrylic, Resin.
灰黑色, 中等, 用于平整和预抛光树脂胶托。
Yellow, Fine, to final polish Acrylic, Resin.
黄色, 较细, 用于树脂胶托的细抛光。

Coarse Green
绿色、粗糙
100



Item No.	Rh 101 C	Rh 103 C	Rh 106 C	R103C	R106C	R209C
Size (mm)	13 x 20	10x 24	10 x 20	10x 24	10 x 20	6 x 13

Medium Gray
灰黑、中等粒度
100



Item No.	Rh 101 M	Rh 103 M	Rh 106 M	R103M	R106M	R209M
Size (mm)	13 x 20	10x 24	10 x 20	10x 24	10 x 20	6 x 13


Fine Yellow
黄色、较细
100



Item No.	Rh 101 F	Rh 103 F	Rh 106 F	R103F	R106F	R209F	R202F
Size (mm)	13 x 20	10x 24	10 x 20	10x 24	10 x 20	6 x 13	4 x 12

Green Stones
Contour Ceramic


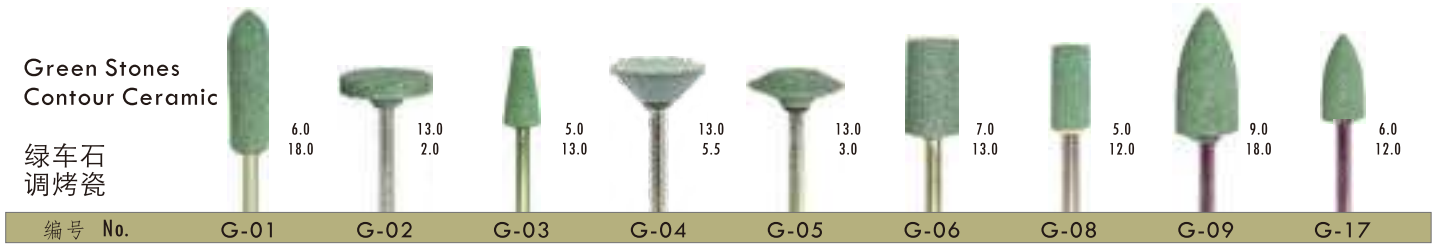
绿车石
调烤瓷




Pink Stones
for Metal
粉车石
调金属,合金



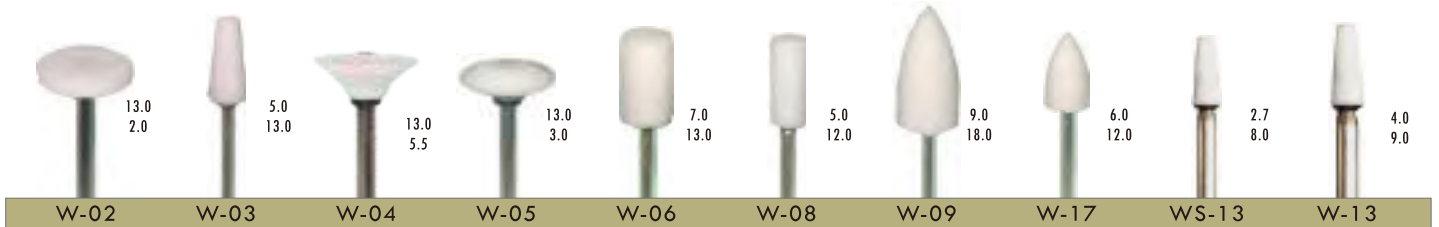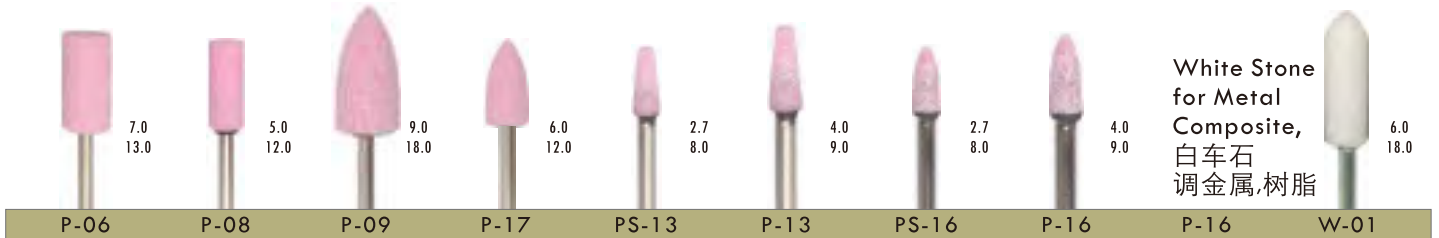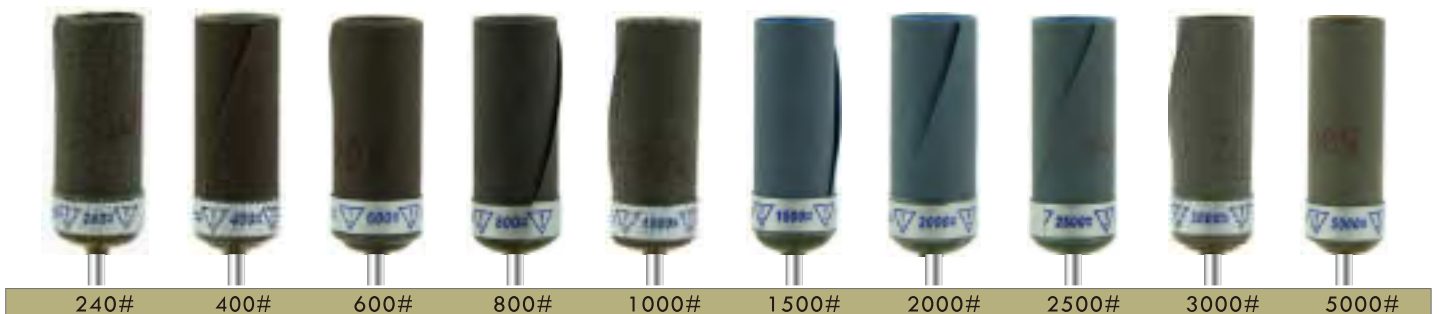White Stone
for Metal
Composite,
白车石
调金属,树脂








Mounted Sandpaper 砂纸棒



Polishing Brushes 抛光刷

				
百洁布 暗红,中粗	百洁布 红,硬	百洁布 棕,中粗	百洁布 紫,软	百洁布 蓝,普通
01MT25 Red.Soft	01Mp25 Red.coarse	01My25 Brown. Medium	01Ms25 Purple.fine	01M22 Blue

						
碳化硅 磨料丝	氧化铝 磨料丝	曲钢丝	曲铜丝	黑猪鬃	灰马毛	灰羊毛
01SiC22 Silicon Carbide	01Al22 Al. Oxide	01A22 Crimped steel wire	01B22 Crimped brass wire	01D22 Chungking Bristle	01E22 Grey Pony hair	01F22 Grey goat hair

					
白羊毛	灰羊毛 夹羊皮	黑猪鬃 19mm	灰马毛	灰羊毛	白羊毛
01G22 White goat hair	01Fr22 Goat hair / chamois	01D19 Chungking Bristle	01E19 Grey Pony hair	01F19 Grey goat hair	01G19 White goat hair

					
黄棉布	白棉布	橘红羊皮	鹅黄羊皮 双面绒	羊绒	白棉线
01P22 Yellow cloth	01Q22 White cloth	01R22 Chamois leather	01Ry22 Microfiber leather	01S22 Cashmere	01T22 Cotton yarn




						
黑猪鬃	灰马毛	白羊毛	曲钢丝	曲铜丝	黑猪鬃	灰羊毛
02D22 Chungking Bristle	02E22 Pony hair	02G22 White goat hair	04A12 Steel	04B12 Brass	04D12 Bristle	04F12 Goat Hair

						
碳化硅 磨料丝	氧化铝 磨料丝	曲钢丝	曲铜丝	黑猪鬃	灰马毛	灰羊毛
05SiC05 SiC	05Al05 Al. Oxide	05A05 Steel	05B05 Brass	05D05 Bristle	05E05 Pony Hair	05F05 Goat Hair

Polishing brush & Wheel
in Bristle, Cotton Cloth
抛光刷



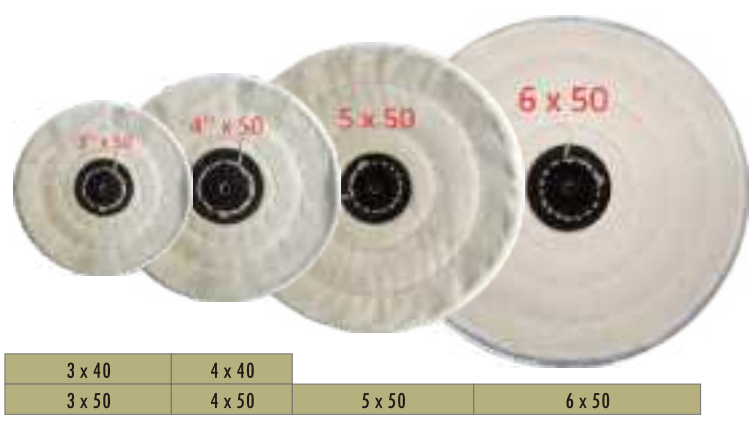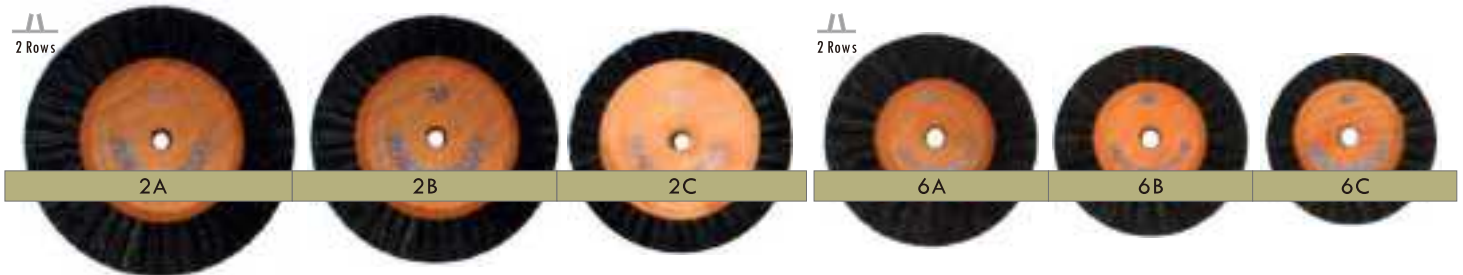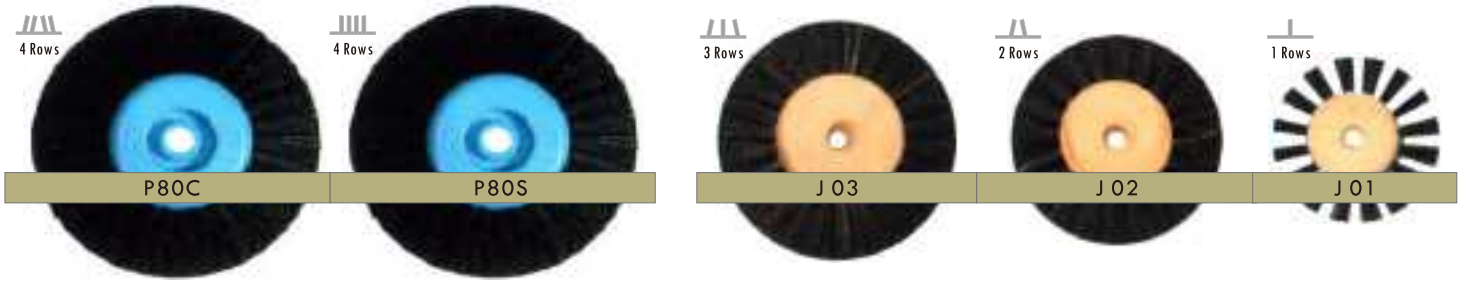
These brushes are made of bristle hair or White, yellow pure cotton cloth, assembled on polishing latch to finish or polish Plaster Moulds, Acrylic 台式抛光机配用的抛光刷/布轮,用于抛光塑料,树脂,模型等。



3 x 40	4 x 40		
3 x 50	4 x 50	5 x 50	6 x 50



4 x 50	5 x 50	6 x 50
--------	--------	--------

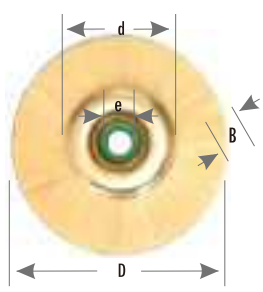


Item	2A	2B	2C	6A	6B	6C	J03	J02	J01	P80	03B48
D	79	73	68	64	60	54	70	58	51	80	48
d	48	48	48	35	35	35	35	30	25	40	25
e	6.5	6.5	6.5	6.5	6.5	6.5	8.5/5	8/5	7/5	6.5	6.35
B	16	16	16	16	16	16	24	16.5	15	36	9
Rows	2	2	2	2	2	2	3	2	1	4	1

Converging Straight

Felt Bobs & Burs 毛毡磨头




03N25	03Ns25	3P	2P	1P	1R
25 x 5 mm, 硬质	25 x 5 mm, 软质	15 x 25 mm	20 x 40 mm	22 x 50 mm	22 x 50 mm

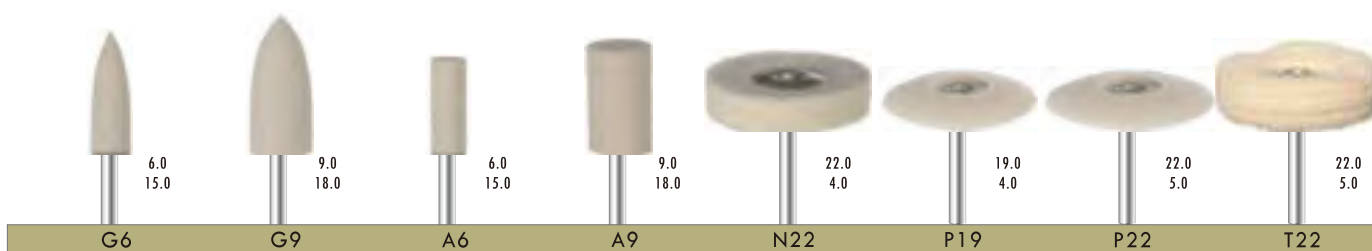


G6	G9	A6	A9	N22	P19	P22	T22
6.0 15.0	9.0 18.0	6.0 15.0	9.0 18.0	22.0 4.0	19.0 4.0	22.0 5.0	22.0 5.0

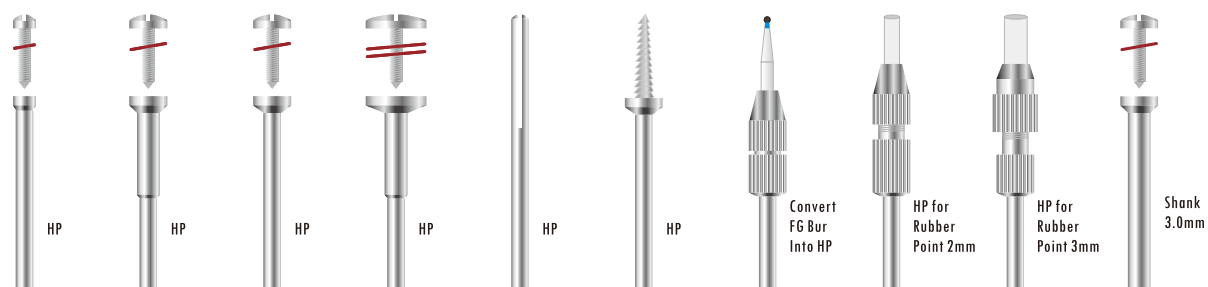Mandrels Stainless Reinforced 不锈钢夹针



编号 No.	32HS	51HS	52HS	81HS	25HS	53HS	55HS	D20HS	D30HS	52D
Head size	3.0mm	5.0mm	5.0mm	8.0mm	2.5mm	5.0mm	5.0mm	5.0mm	6.0mm	5.0mm

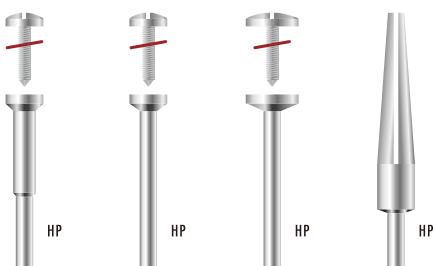Mandrels Steel 普通夹针



编号 No.	51H	52H	72H	54H
Head size	5.0mm	5.0mm	7.0mm	5.0mm

Hard Z Cutter 硬锆磨头

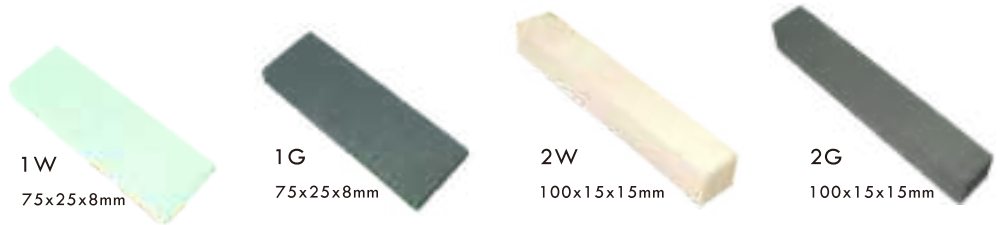
Made of Nano-grade Zirconia, used to grinding, finishing or polishing Ceramics and Acrylic, Resin, Zirconia etc. Also applied to Soft materials such as wood, Ceramic or stones etc

用于打磨,成型,抛光烤瓷,树脂,氧化锆(烧结前)。特殊的纳米级氧化锆胚体制而成。



Sharpening Stones 磨刀石

Adjust to minimum speed;
Touch the stones softly.
Using water will be better



Cleaning Brushes By Hand 清洁刷(手用)



Stainless Steel Disc 不锈钢小锯片

for Soft Plastic, Acrylic, Wood etc.



Stone Disc 刚玉车片



3Z09.1

Zirconia Luster Polishing Kit
氧化锆抛光套装



- + 2 pieces of Diamond Stones DG001 , DG012
- + 2 pieces of Diamond Polishers, Blue, D203M, D3M
- + 2 pieces of Diamond Polishers, Yellow, D203EF, D3EF

3Z00.2

Grinding and polishing Kit of Zirconia
氧化锆打磨抛光套装



- + 2 pieces of Tungsten carbide Burs HP 699, F040SFE
- + 1 set of Diamond Disc Em19D20 Mounted
- + 4 pieces of HP Diamond Burs: MK37, FM27, FD30, MD14
- + 4 pieces of Diamond Stones: DG017, DG002, DG004, DG005
- + 4 pieces of Diamond Polishers: D205M, D205F, D205EF, D201EF.

3Z00.3

Grinding and Polishing Kit for Zirconia
氧化锆打磨抛光套装



- + 1 piece HP Tungsten Carbide Burs HP 699
- + 5 pieces Diamond Burs: MG50, MD14, FD30, FM27, MK37;
- + 1 set of Diamond Disc Hm19D20 Mounted
- + 1 piece of Diamond Polisher D1M
- + 1 set of Silicon polisher C104C Mounted

3L04.2

Grinding and Polishing Kit for Lithium Disilicate
玻璃陶瓷打磨抛光套装



- + 1 set of Diamond Disc HM19D20, Mounted,
- + 3 pieces of Diamond Stones DG017, DG004, DY001EF
- + 3 pieces of HP Diamond Burs MD10, MB12, MR23
- + 3 pieces of Diamond Polisher D1M, D13F, D13UF
- + 1 piece of Silicon Polisher C205C

3L04.3

Grinding and Polishing Kit for Lithium Disilicate
玻璃陶瓷打磨抛光套装



- + 1 set of Diamond Disc HM19D20, Mounted,
- + 2 pieces of Diamond Stones DG004, DY001EF,
- + 3 pieces of HP Diamond Burs MR23, ML16, MM27,
- + 1 piece of Diamond Polisher D1M
- + 1 piece of Silicon Polisher C205C

3MC02.3

Adjusting and Polishing Kit for Ceramics
烤瓷工件打磨抛光套装



- + 1 piece of HP Diamond Burs MR23
- + 1 piece of HP Silicon Polisher C203C
- + 1 piece of Diamond Polisher D26f
- + 2 pieces of Polishing Brushes 01G22, N22

3MR09.1

Adjust & Polish Invisible Braces
隐形牙套调整抛光套装



- + 1 piece of Tungsten Carbide Burs AL02354
- + 1 set of Diamond Disc C19D20, Mounted
- + 2 pieces of Polishing Brushes 01Ms25, 01S22

3A00.1

Popular HP Kit for grinding or polish Metal or Alloy
常用打磨抛光金属、合金的套装



* Universally grind, remove, polish and finish of Alloy, metal, Jewellery, stones etc. It widely applied to Laboratory, DIY or Hobbies.
+ 5 pieces of HP sintered diamond burs ,
+ 4 pieces of HP Tungsten Carbide Cutters
+ 10 pieces of HP Tungsten Carbide burs

+ 4 pieces of White stones
+ 2 pieces of Pink stones
+ 10 pieces of HP Silicon Polishers
+ 3 pieces of Rubber wheel mounted
+ 1 pieces of Woolen Felt wheel mounted ,
+ 5 pieces of HP Polishing brushes

3MC00.1

Popular HP Kit for grinding or polish Ceramics / Porcelain
常用打磨抛光陶瓷的套装



* Universally grind, remove, polish and finish of Ceramics or Porcelains etc. It widely applied to Laboratory, DIY or Hobbies.
+ 10 pieces of HP Diamond burs ,
+ 4 pieces of Diamond Discs mounted
+ 6 pieces of HP Green Stones

+ 3 pieces of Silicon Rubber wheel mounted
+ 8 pieces of Silicon Rubber Polishers
+ 4 pieces of Polishing brushes

3LR00.1

Popular HP Kit for grinding or polish Acrylic, Resin
常用打磨抛光树脂胶托的套装



* Universally grind, remove, polish and finish of Acrylic, Resin etc. It widely applied to Laboratory, DIY or Hobbies.
+ 2 pieces of HP Tungsten Carbide Cutters
+ 32 pieces of HP Silicon rubber Polishers

+ 10 pieces of Polishing brushes
+ 3 pieces of Woolen Felt polishing brushes

3Z00.1

Popular HP Kit for grinding or polish Zirconia, Lithium Disilicate
常用打磨抛光氧化锆、玻璃陶瓷的套装



* Universally grind, polish and finish of Zirconia and Lithium Disilicate etc. It widely applied to Laboratory, DIY or Hobbies.
+ 4 pieces of HP Diamond Burs
+ 2 pieces of Diamond Discs Mounted

+ 6 pieces of Diamond Stones
+ 4 pieces of Diamond polishing wheels
+ 2 pieces of Diamond polishing burs
+ 5 pieces of Diamond Twist polishers
+ 1 pieces of Polishing brush .

ZS100

Hot HP Mounted stones mixed
车石套装,常用。




* 100 Pcs, Hot Mounted stones in Green, Pink, White , applied to Alloy, metal, ceramics, Jewellery, stones, hard wood etc.
+ 30 pieces of HP Green Stones
+ 30 pieces of HP Pink Stones

+ 30 pieces of White Stones
+ 10 pieces of Brown Stones

ZD 101

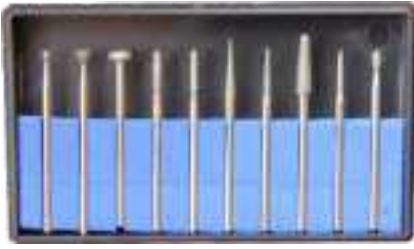Popular HP diamond burs KIT
常用金刚石磨头套装



MD23, MK47, MJ60, MA23, MC18, MM27, MB21, ML40, MG16, MR23

10 Hot Spec. of HP Diamond, universally grind, remove, polish and finish of Alloy, metal, Jewellery, stones etc. It widely applied to Laboratory, DIY or Hobbies.

用途: 选取10支最常用的金刚石磨头, 形状全, 大小适中, 可以用来处理陶瓷工件。

ZD 101F

Popular HP diamond burs KIT
常用金刚石磨头套装



FD23, FK47, FJ60, FA23, FC18, FM27, FB21, FL40, FG16, FR23

10 Hot Spec. of HP Diamond of Fine grit, universally grind, remove, polish and finish of Zirconia, Lithium Disilicate, Porcelain etc. It widely applied to Laboratory, DIY or Hobbies.

用途: 选取10支最常用的金刚石磨头, 形状全, 大小适中, 可以用来处理氧化锆、玻璃陶瓷, 烤瓷等工件。

ZD 102

Small & Tiny Cut, HP diamond burs KIT
金刚石磨头, 小号套装



MD10, MD14, MD23, MK16, MA16, MA23, MC25, MM12, MM18, MM27, MB12, MB16, MB21, ML16, ML23, MN18, MG16, MR23, MQ23, MU60.

20 Spec. small or sharp burs. To create tiny hole, slim curves or tiny, deep grind smartly. Used to Alloy, metal, ceramics, Jewellery, stones, hard wood etc. It widely applied to Laboratory, DIY or Hobbies.

用途: 20支金刚石磨头, 用来开小洞, 孔, 细边, 细槽, 或者小的磨削。很轻巧, 灵活。适用于各种金属, 合金, 陶瓷, 珠宝, 宝石, 木头等。

ZD 103

Fast Trimming, HP diamond burs KIT
金刚石磨头套装, 快速磨削



CQ23, CL40, SD35, SA50, SC40, SP70, SG50, SV70, SV60, SF60

20 Spec. HP Diamond for fast trimming, Grinding of Plastic, Jewellery, stones, Hard wood, etc. Save time when removing and work efficiently.

用途: 10支常用的金刚石磨头, 头部相对比较大, 可以快速磨削塑料, 木头, 陶瓷和石头等等, 适用于牙科, 珠宝加工和DIY爱好者。

ZD 104

Hot Fully-Sintered diamond KITs
全砂烧结磨头套装, 常用



C 050S, F050S, L05S, B050C, A037, G050, C037F, F050F, M050EF, G050EF

10 spec. Hot Fully-sintered Diamond HP burs, universal grind, remove, polish and finish of Alloy, metal, Jewellery, stones etc. It widely applied to Lab, DIY or Hobbies.

用途: 10个常用全砂烧结磨头, 5种不同粗细。广泛适用于研磨, 抛光金属, 合金, 珠宝, 宝石。经久耐用, 方便实惠。

ZD 105

Fast Trimming Fully-Sintered diamond KIT
全砂烧结磨头套装, 大磨削



K080S, K100S, J100S, C070S, M070S, B070S, L070S, G060S, F070S, R070S

If workpieces are extremely hard and waste too much time when grinding, this fast trimming kit may help to work efficiently. Applied to Alloy, metal, Jewellery, stones, hard wood, etc.

用途: 10个大头的全砂烧结磨头, 用于粗磨, 去块, 效率高, 经济实惠。

ZT 101

Popular HP Tungsten carbide burs KIT
常用低速钨钢针



HP 1/4, HP 1/2, HP 1, HP 2, HP 3, HP 4, HP 5, HP8, HP 34, HP 36, HP 38, HP 41, HP 556, HP 557, HP 560, HP 563, HP 699, HP 701, HP 702, HP 704

20 Spec. HP Tungsten Carbide burs are used to create tiny hole, slim curves or tiny, deep grind very smartly. It is applied to Alloy, metal, ceramics, Jewellery, stones, hard wood etc.

用途: 20个型号钨钢车针, 用于雕牙花, 钻小孔, 洞, 开细槽, 细线等, 金属, 合金, 宝石, 坚木都可以用。

ZT 102

Hot Tungsten carbide cutters Kit
钨钢磨头套装, 常用



D050SP, L040SC, L060SC, F060SC, B023SC, L023SC, L040FC, C023FC, L040EFC, F060CC

10 spec. Tungsten carbide Cutter, universal grind, remove, polish and finish of Acrylic, Resin, stones, soft metal or wood.

用途: 10个常用钨钢磨头, 形状全, 大小适中。主要用来研磨, 抛光塑料, 树脂, 软金属, 石头等, 方便实惠。

ZF 101

Acrylic Polishing Kit
树脂基托抛光套装



F060SC, L040SC, U12 x 2, U13 x 2, U16 x 2

Simple and straight way to adjust Acrylic without the problem associated with pumice.

用途: 2个型号的钨钢磨头快速磨削, 三个不同粒度的硅胶磨头依次进行抛光。

ZA 101

Non-precise metal, Alloy Adjust and polish
非贵金属, 铬钴调整与抛光



W-01, W-09, W-08, W-16, W-03, W-06
A201C, A203C, A206C, A201F, A 203F

To adjust and polish Cobalt-chromium, nickel chromium Alloy. White stones for adjusting or grinding, Rubber burs for polishing.

用途: 试用于铬钴、镍铬合金等非贵金属的研磨和抛光
6支白车石用来研磨和调整, 5支硅胶磨头用来抛光。

ZP 100

Sandpaper Polishing burs Kit
砂纸棒套装



240#, 400#, 600#, 800#, 1000#, 1500#, 2000#, 2500#, 3000#, 5000#

10 different Grit Sandpaper burs are used for polishing Alloy, metal, Jewellery stones, hard woods etc. Working step by step, mirror surface in comming soon, water cooling or dry are acceptable.

用途: 10个不同粒度的砂纸棒, 从240目到5000目不等, 请根据需要, 选取合适的粒度, 耐心细致的使用。

Z 2001

Popular Felt Polishing bur Kit
毛毡抛光磨头套装



12 pieces of High quality Felt polishing burs are widely used to mould, Alloy, metal, Jewellery, stones, Hard woods etc.

12 pieces of High quality Felt polishing burs are widely used to mould, Alloy, metal, Jewellery, stones, Hard woods etc.

用途: 12支常用的毛毡磨头, 广泛适用于各种模具, 金属, 陶瓷, 珠宝, 宝石和木头等。

Samples Collection

样品册

- Flip-Cart Plastic Pages ;
- Smart, Easily fill all burs ;
- Economic and Recyclable ;
- Nice-looking for showing to Customers;
- Convenient to take to anywhere .



No. 171B5

Empty Collection Wallet
样品本(不带车针)



4 Colors Available ,
100 holes per page,
5 pages per set,
Smart to HP Burs



No. 1704D

Sintered Diamond Burs Collection
全砂烧结磨头样品本



115 specifications ,
Fully Sintered Diamond,
HP shanks



No. 1705D

HP Diamond Burs, Mounted Stones
金刚石磨头, 车石样品本



More than 100 specifications ,
Electroplating Diamond burs
HP Shanks



No. 1707T

Tungsten carbide Burs Collection
高品质钨钢磨头样品本



More than 100 specifications ,
High Quality TC burs,
HP Shanks




No. 1708P

Polishing Burs, Mandrels Collection
抛光磨头, 夹针样品本



More than 100 specifications ,
Silicon, Rubber, Diamond burs
HP Shanks



隐形义齿胶


FlexSoft Partial Dentures





The alternative to conventional Partial Dentures ;
Flexible . Soft. Easy-to-finish.Easy-to-polish
弹性好, 易成型, 易抛亮, 美观大方。



编号 No.	规格 Spec.
D2	Crystal ; 18 gs
C1	Ligth Pink ; 9 gs
C2	Ligth Pink ; 18 gs
B1	Standard Pink ; 9 gs
B2	Standard Pink ; 18 gs
A1	Deep Pink ; 9 gs
A2	Deep Pink ; 18 gs
For bulk pak of 1 kilo.	
C3	Ligth Pink ; 1000 gs
B3	Standard Pink ; 1000 gs
A3	Deep Pink ; 1000 gs

 Crystal (D)
(透明)

 Light pink (C)
(浅粉红)

 Standard pink (B)
(粉红)

 Deep pink (A)
(深粉红)



Instructions for Use

Suitable for: Partial and complete dentures, occlusal splints and nightguards.

Not Suitable for: Small unilateral restorations.

Warnings: Use eye protection and heat-resistant gloves when injecting.

Use proper ventilation, vacuum system and masks when grinding prosthodontic resins.

Precautions: Cartridge must be injected within 1 hour of opening the foil pouch.

I. PREPARING THE MASTER MODEL FOR DUPLICATION

All undercuts that will affect the path of insertion should be blocked-out with wax prior to duplication. Do not relieve necklines of teeth or areas desirable for retention. For free-end saddles, it is not necessary to remove the distal undercut of the last natural tooth, unless the area is too deep. For free-end saddle cases where anterior teeth are also being replaced, undercuts in the anterior should also be blocked out.

II. SET-UP

Duplicate blocked-out master model. Articulate the duplicate model and the master model (you will need the master model articulated at finishing steps to check occlusion). Set up the case on the duplicate model using only wax (no base plate). Allow 1 mm of thickness under the teeth to flow to the labial side and into the retention diatorics in the teeth. Use a small round carbide bur to create mechanical retention in the teeth, from mesial to distal and connected to the center diatoric hole.

III. WAXING

Important: Wax the case with anatomical features to reduce grinding during the finishing stage.

Recommended thickness of wax-up:

1. Palatal thickness 1.2 - 1.5 mm
2. Buccal & Labial Flange 1.2 - 2.0 mm
3. Clasps 1.0 - 1.5 mm
4. Lingual Flange 1.2 - 2.0 mm
5. Lingual Major Connector (Mandibular) 1.7 - 2.0 mm

IV. INVESTING AND SPRUING

Invest the case in the center of the flask, aiming the posterior portion toward the sprue opening. After the stone (or plaster) has set, use two sprues from the flask opening to the most posterior portion of the wax on both sides (see diagram). Sprues should be approximately 7 mm. width. Do not create sharp angles with the sprues, as the heated resin must flow in a straight or slightly curved path. Apply a thin layer of petroleum jelly on the entire stone surface. Then invest top half of flask.

V. REMOVING WAX

After the stone has set, place flask in boiling water for 8 minutes. Open the flask and wash out the wax, making sure all wax residues are eliminated. Allow the flask to cool. Apply thin layers of separator on the entire surface of the stone and allow it to dry completely. Do not apply separator on teeth. When dry, close flask.

VI. INJECTING

Wear safety goggles and heat-protective gloves for this procedure. Spray inside of the cartridge sleeve with mold release spray and place it into furnace. Turn on the furnace and allow it to reach approximately 550°F (290°C), preheating cartridge sleeve for approximately 15 minutes.

Remove cartridge from foil-lined bag (do not remove cartridge from bag until this step). Spray the cartridge and bronze disc with mold release spray. Insert the cartridge into the cartridge sleeve with the cap end first. Insert the bronze disc into the cartridge sleeve against the end of the cartridge. Use the thick bronze disc for medium tubes and the thin bronze disc for the large tubes. Heat cartridge for 12 minutes.

Position the flask with the sprue end up, directly under the shaft of the press. After 12 minutes, remove the cartridge sleeve from the furnace. Keep the cartridge sleeve in a horizontal position until resting on top of the flask, and then turn the cartridge sleeve into a vertical position. Immediately use the press to inject the heated cartridge into the flask using a steady fast motion. Continue applying pressure until tube bursts and continue with maximum downward pressure for 3 minutes.

Remove the cartridge sleeve from the flask with a slight twisting motion. After removing the cartridge sleeve, let the flask cool for at least 30 minutes. Do not put flask into a cold water bath. Use a plaster nipper to de-flask. Do not use a hammer. Remove the compressed cartridge from the cartridge sleeve with a knock-out punch. Important: Cartridge sleeve must be pre-heated for 10 minutes BEFORE the next injection.

VII. FINISHING AND POLISHING

Remove appliance from the model. Cut off the sprues with a cut-off wheel. Trim the edges of appliance with a coarse grinding wheel. Avoid overheating. Use rubber wheels to smooth scratches. Place the appliance on the articulated master model to check the fit, retention, thickness and occlusion. Clean the tissue surface with a small bristle brush.

To polish, use a coarse dental abrasive, such as pumice, with a wet rag wheel. Apply final shine with a hi-shine compound on low speed with a final buff wheel. Be careful not to overheat appliance.

铝桶


Empty Aluminum Cartridge

编号 No.	规格 D x H, mm
25085	25 x 85 mm
25100	25 x 100 mm
25CAP	24.2 x 16 mm, Cap
22065	22 x 65 mm
22085	22 x 85 mm
22100	22 x 100 mm
22120	22 x 120 mm
22CAP	21.2 x 16 mm, Cap





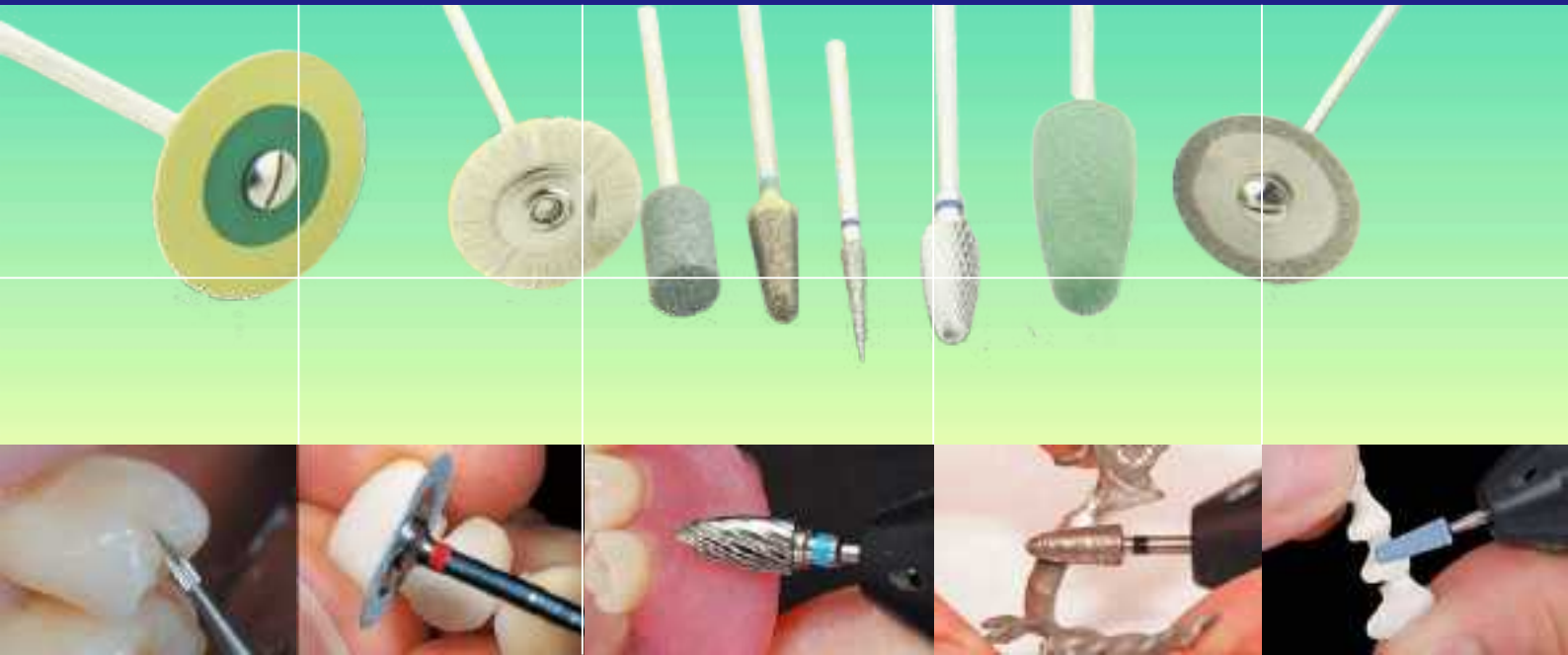BurDental



B&D Technologies Co., Limited

WebSite: www.burdental.com

Tel/WhatsApp: 0086-1802 8579 361

Address: 4th Floor, Blk 2, Zhao Yuan Business Plaza, Pan Yu Distr., GuangZhou, Guangdong Province, China 510430

Email: info@burdental.com

Skype: burdental